

600

SERIES

CLD



Chemiluminescent Analyzer

USER'S MANUAL

The Model 600 CLD Series Instruments starting with Serial Number UO6081 have several new Hardware and Software features. For a complete explanation, see section 13.5 starting on page 70



Safety Alert
Caution or Warning



Temperature Hazard
Caution or Warning



Electrical Shock Hazard
Caution or Warning

Safety Information in this Manual


Note, caution and warning symbols appear on the instrument and throughout this manual to draw your attention to important operational and safety information.

A “**NOTE**” marks a short message to alert you to an important detail.


A “**CAUTION**” safety alert appears with information that is important for protecting your equipment and performance.

A “**WARNING**” safety alert appears with information that is important for protecting you, others and equipment from damage. Pay very close attention to all warnings that apply to your application.




The  symbol (an exclamation point in a triangle) precedes a general CAUTION or WARNING statement.



The  symbol (wavy vertical lines with an under score in a triangle) precedes an elevated temperature hazard CAUTION or WARNING statement.



The  symbol (a lightning bolt in a triangle) precedes an electric shock hazard CAUTION or WARNING statement.

Some or all of the above symbols may appear in this manual or on the equipment. This manual should be consulted whenever one of these symbols is encountered on the equipment.

ALWAYS REMOVE POWER BEFORE CONNECTING OR DISCONNECTING SIGNAL CABLES OR WHEN SERVICING THE EQUIPMENT.

The 600 series CLD instruments meet or exceed the following directives and standards.

Application of Council Directive(s):

Electrical Safety:

Low Voltage Directive 73/23/EEC

Electromagnetic Compatibility:

EMC Directive 89/336/EEC

Standard(s) to which Conformity is Declared:

Electrical Safety:

*Standard for Electrical Equipment for Measurement, Control, and Laboratory Use
[EN 61010-1:2001 (2nd Edition)]*

Electromagnetic Compatibility:

EN 61326:1997 Electrical equipment for measurement, control and laboratory use - EMC requirements (Amendment A1: 1998 to EN 61326:1997; Amendment A2:2001 to EN 61326:1997)



POSSIBLE EXPLOSION HAZARD

Do not apply power to the analyzer or attempt to energize the ozone supply or converter until **ALL** leak checks have been performed and until the analyzer environment has been determined to be non-hazardous.

This analyzer is designed for use in a **NON-HAZARDOUS** environment.

This analyzer is designed for use with a **HAZARDOUS** sample.



Tampering or use of substitute components may cause a safety hazard. Use only factory authorized replacement parts.



Do not operate without the cover secured. Servicing requires access to live electrical components which can cause death or serious injury. Refer servicing to qualified service personnel. For safety and proper performance, this instrument must be connected to a properly grounded three-wire receptacle.

TABLE OF CONTENTS

| <u>Section</u> | <u>Title</u> | <u>Page</u> |
|-----------------------|--|--------------------|
| 1. | INTRODUCTION ----- | 1 |
| 1.1. | Overview----- | 1 |
| 1.2. | Unpacking Instructions----- | 1 |
| 1.3. | Reporting Damage----- | 1 |
| 1.4. | Contact Information----- | 1 |
| 1.5. | Warranty Certificate----- | 2 |
| 2. | FEATURES AND SPECIFICATIONS ----- | 3 |
| 3.1. | Descriptions----- | 3 |
| 2.2. | Features----- | 3 |
| 2.3. | Specifications----- | 4 |
| 3. | INSTALLATION ----- | 5 |
| 3.1. | General----- | 5 |
| 3.2. | Site and Mounting ----- | 5 |
| 3.3. | Electrical----- | 5-6 |
| 3.4. | OutputConnections ----- | 6 |
| 3.5. | Required Gases----- | 7 |
| 3.6. | Gas Handling Equipment----- | 7 |
| 3.7. | Gas Connections ----- | 5 |
| 3.8. | Sampling Requirements ----- | 7 |
| 3.8.1. | Filtration ----- | 7 |
| 3.8.2. | Condensation ----- | 7 |
| 3.8.3. | Presence of Corrosive Gases----- | 7 |
| 3.8.4. | Gas Temperature ----- | 7 |
| 3.8.5. | Flow Rate ----- | 7-8 |
| 3.8.6. | Sample Gas Outlet ----- | 8 |
| 4. | BASIC OPERATION ----- | 9 |
| 4.1. | Display ----- | 9 |
| 4.2. | Menu Tree----- | 10 |
| 4.3. | Keyboard ----- | 11 |
| 4.4. | Operation with Cursor Keys ----- | 11 |
| 4.5. | Operation with Function Keys----- | 11 |
| 4.6. | Read and Change Parameters ----- | 11 |
| 5. | OPERATING STRUCTURE ----- | 12 |
| 5.1. | Main Menu ----- | 12-13 |
| 6. | MENU STRUCTURE ----- | 13 |

| | | |
|-----------|---|-----------|
| 6.1. | User Functions (Level 1) ----- | 13 |
| 6.2. | Advanced User Functions (Level 2) ----- | 13 |
| 6.3. | Maintenance User Functions (Level 3) ----- | 13 |
| 6.4. | System User Functions (Level 4)----- | 13 |
| 7. | MAIN MENU ----- | 15 |
| 7.1. | F1 Measurements ----- | 15 |
| 7.1.1. | F1 NO or NOx Measurement ----- | 16 |
| 7.1.2. | F2 NO + NO2 Measurement ----- | 15-16 |
| 7.1.3. | F3 Diagnostics----- | 16 |
| 7.1.4. | Arrow Keys to Cycle Ranges and Auto Range ----- | 17 |
| 7.2. | F2 Purge Analyzer ----- | 18 |
| 7.3. | F3 Diagnostics ----- | 18 |
| 7.4. | F4 Calibrations ----- | 19 |
| 7.4.1. | F1 Automatic Calibration ----- | 19-20 |
| 7.4.2. | F2 Manual Calibration ----- | 20-21 |
| 7.4.3. | F3 Display Deviations----- | 21-22 |
| 7.4.4. | F4 Check Calibration Values----- | 22 |
| 7.4.5. | F5 Reset Calibration Values ----- | 22 |
| 7.4.6. | F6 Range Select ----- | 22 |
| 7.5. | F5 Setup ----- | 23 |
| 7.5.1. | F1 Span Gas Concentrations ----- | 24 |
| 7.5.2. | F2 Calibration Settings ----- | 24 |
| 7.5.2.1. | F1 Times ----- | 24 |
| 7.5.2.2. | F2 Measuring Deviations----- | 25 |
| 7.5.2.3. | F3 Deviations----- | 25 |
| 7.5.2.4. | F4 Calibration via Valves----- | 25 |
| 7.5.2.5. | F5 Calibration via Probe----- | 25 |
| 7.5.3. | F3 Range Limits ----- | 25 |
| 7.5.3.1. | F1 Change Upper Range Limits----- | 25 |
| 7.5.3.2. | F2 Change Auto Range Limits----- | 26 |
| 7.5.4. | F4 Change Alarms----- | 26 |
| 7.5.5. | F5 Password ----- | 27 |
| 7.5.5.1. | F1 Enter Password ----- | 27 |
| 7.5.5.2. | F2 Change Password ----- | 27 |
| 7.5.5.3. | F3 Reset Passwords ----- | 28 |
| 7.5.6. | F6 Linearization ----- | 28 |
| 7.5.7. | F7 System Settings----- | 29 |
| 7.5.7.1. | F1 Clock Setup ----- | 29-30 |
| 7.5.7.2. | F2 TCP/IP Address ----- | 31 |
| 7.5.7.3. | F3 Output Signal Assignments ----- | 31 |
| 7.5.5.4. | F4 Output Ranges ----- | 31 |
| 7.5.7.5. | F5 Protocol On/Off ----- | 32 |

| | | |
|------------|---|-----------|
| 7.5.7.6. | F6 Set Language | 32 |
| 7.5.7.7. | F7 Automatic Setup | 32 |
| 7.5.8. | F8 Measure Settings | 33 |
| 7.5.8.1. | F1 Set Purge Time Before Measure | 33 |
| 7.5.8.2. | F2 Set Converter Efficiency | 33 |
| 7.5.8.3. | F3 Set Time Constant | 34 |
| 7.5.8.4. | F4 Set Purge Time Before Calibration | 34 |
| 7.5.10. | F10 Display Software Version | 34-35 |
| 7.6. | F7 Remote/Manual Operation | 35 |
| 7.7. | F8 Standby | 35 |
| 8. | ANALYZER COMPONENTS | 36 |
| 8.1. | Rear Panel Connections | 36 |
| 8.1.1. | Main Analog Connections | 37 |
| 8.1.2. | Auxiliary Analog Connections | 37 |
| 8.1.3. | Digital (RS-232) Connections | 38 |
| 8.1.4. | Digital (TCP/IP) Connections | 38 |
| 8.2. | Internal Components | 39 |
| 8.3. | Main Electronics Board (Potentiometers) | 40 |
| 8.4. | Main Electronics Board (Connections) | 41 |
| 8.5. | Reaction Chamber | 42 |
| 8.6. | Relay Board (Connections) | 43 |
| 9. | OPERATION | 44 |
| 9.1. | Preparation | 44 |
| 9.2. | Operation | 44-45 |
| 9.3. | Shut Down Procedure | 45 |
| 10. | FUNCTIONAL DESCRIPTION | 46 |
| 10.1. | Operating Principle | 46 |
| 10.2. | Reaction Chamber Assembly | 46 |
| 10.3. | Flow System | 46 |
| 10.4. | Electronics | 46 |
| 10.5. | Main Electronics Board | 47 |
| 10.6. | Relay Board | 47 |
| 11. | REACTION CHAMBER | 48 |
| 11.1. | Reaction Chamber Disassembly | 48 |
| 11.2. | Reaction Chamber Assembly | 48 |
| 12. | TROUBLESHOOTING | 49 |
| 12.1. | Ozone Air/O ₂ Supply | 49 |
| 12.2. | Sample Supply | 49 |
| 12.3. | NO/NO _x Converter | 49 |
| 13. | APPENDIX 1 | 50 |
| 13.1. | AK Protocol | 51-62 |

| | | |
|-------|-------------------------------|--------------|
| 13.2. | I/O Lists ----- | 63 |
| 13.3. | Flow Diagrams ----- | 64-67 |
| 13.4. | Electrical Block Diagram----- | 68 |
| 13.5 | Serial Number UO6081----- | 70 |

1. Introduction

1.1. Overview

Congratulations and thank you! You have just purchased one of the most reliable gas analyzers in the world. Before using the analyzer, please familiarize yourself with its operation by reading this manual. If you have any questions, please do not hesitate to call California Analytical Instruments for assistance. We want you to be a member of our thousands of satisfied customers.

1.2. Unpacking Instructions

Open the shipping container and carefully remove the analyzer from the packing materials. Inspect the instrument for any sign of damage. Remove the Top Cover retaining screws. Visually check for loose parts or connectors that are not properly seated. Verify all circuit boards and circuit board connections are secure. If all internal components look normal, re-install the cover.

1.3. Reporting Damage

Should there be any apparent damage to either the inside or outside of the instrument due to shipping or handling, immediately notify the shipper. The shipping container or packing materials should be retained for inspection by the shipper.

1.4. Contact Information

California Analytical Instruments, Inc.
1312 West Grove Avenue
Orange, CA 92865
713 974-5560
Fax 713 921-2531
Website: www.gasanalyzers.com

1.5. Warranty Certificate

Subject to the exceptions and upon the conditions stated below, California Analytical Instruments (CAI) warrants that the products sold under this sales order shall be free from defects in workmanship and materials for one year after delivery of the product to the original Buyer by CAI and if any such product should prove to be defective within such one year period, CAI agrees, at its option, either (i) to correct by repair or, at CAI's election, by replacement with equivalent product any such defective product, provided that investigation and factory inspection discloses that such defect developed under normal and proper uses, or (ii) to refund the purchase price. The exceptions and conditions mentioned above are as follows:

- a) components or accessories manufactured by CAI which by their nature are not intended to and will not function for one year are warranted only to give reasonable service for a reasonable time; which constitutes reasonable time and reasonable services shall be determined solely by CAI. A complete list of such components and accessories is maintained at the factory;
- b) CAI makes no warranty with respect to components or accessories not manufactured by it; in the event of defect in any such component or accessory CAI will give reasonable assistance to Buyer in obtaining from the respective manufacturer whatever adjustment is authorized by the manufacturer's warranty;
- c) any product claimed to be defective must be returned to the factory transportation charges prepaid and CAI will return the repaired or replaced product freight collect;
- d) if the product claimed to be defective requires on-site repair, such warranty labor will be provided at no charge; however, transportation and living expenses will be charged to Buyer;
- e) if the product is a consumable or the like, it is warranted only to conform to the quantity and content and for the period (but not in excess of one year) stated on the label at the time of delivery or 90 days;
- f) CAI may from time to time provide a special printed warranty with respect to a certain product, and where applicable, such warranty shall be deemed incorporated herein by reference;
- g) CAI shall be released from all obligations under all warranties, either expressed or implied, if any product covered hereby is repaired or modified by persons other than its own authorized service personnel unless such repair by others is made with the written consent of CAI.

IT IS EXPRESSLY AGREED THAT THE ABOVE WARRANTY SHALL BE IN LIEU OF ALL WARRANTIES OF FITNESS AND OF THE WARRANTY OF MERCHANTABILITY AND THAT CAI SHALL HAVE NO LIABILITY FOR SPECIAL OR CONSEQUENTIAL DAMAGES OF ANY KIND OR FROM ANY CAUSE WHATSOEVER ARISING OUT OF THE MANUFACTURE USE, SALE, HANDLING, REPAIR, MAINTENANCE OR REPLACEMENT OF ANY OF THE PRODUCTS SOLD UNDER THIS SALES ORDER. SOME STATES DO NOT ALLOW THE EXCLUSION OR LIMITATION OF INCIDENTAL OR CONSEQUENTIAL DAMAGES, SO THAT THE ABOVE LIMITATIONS OR EXCLUSIONS MAY NOT APPLY. THIS WARRANTY GIVES YOU SPECIFIC LEGAL RIGHTS, AND YOU MAY ALSO HAVE OTHER RIGHTS, WHICH VARY FROM STATE TO STATE.

Representations and warranties made by any person, including dealers and representatives of CAI, which are inconsistent, or in conflict with the terms of this warranty, shall not be binding upon CAI unless reduced to writing and approved by an expressly authorized officer of CAI.

2. Features

2.1. Description

The CAI Model 600 CLD Analyzer is a highly sensitive chemiluminescent (CLD) gas analyzer for measuring oxides of nitrogen gas concentrations in industrial and vehicle emission applications.

2.2. Features-General

The Model 600 CLD analyzer has a 3 by 5 inch liquid crystal display and a 20 key data/operation input keyboard. The 16 bit microprocessor control board consists of the MSR-Card with 16 digital inputs, 16 digital outputs, 16 analog inputs and 4 analog outputs. The analyzer can be manually operated from the keypad or remotely via TCP/IP or RS-232C communications. After turning on the analyzer, it needs at least 30 seconds for initialization. During this time, the screen is illuminated. The analyzer is available with an optional internal heated sample pump and internal zero and span solenoids.

- + **IMPORTANT TIP:** When the analyzer is powered up, it defaults to access level 1 (User). To operate ALL parameters, check the access level. See Section 5.5.5.

The contents of this operators manual include:

- **Specifications**
- **Installation Requirements, Mechanical & Electrical**
- **Operation & Calibration Instructions**
- **Reaction Chamber Description with Procedures for Disassembly of its Component Parts**
- **Function Explanation of the Electronic Circuitry**
- **Electrical Block Diagram**

NOTE: A detailed service and factory setup instruction manual for purchase will be available in April 2003 and will include full schematics.

Model 600 CLD Specifications

| | | |
|--|---|--|
| DETECTOR | Chemiluminescence (CLD) Photodiode (thermally stabilized with Peltier Cooler) | |
| NO/NO_x RANGES | 0-1 to 3,000 ppm NO or NO _x (Four user programmable ranges) (Higher Ranges Available upon Request) | |
| RESPONSE TIME | T90 < 2 Seconds to 60 Seconds Adjustable | |
| RESOLUTION | 10 ppb NO/NO _x (Displays 5 significant digits) | |
| REPEATABILITY | Better than 0.5% of Full Scale | |
| LINEARITY | Better than 0.5% of Full Scale | |
| NOISE | Less than 1% of Full Scale | |
| ZERO & SPAN DRIFT | Less than 1% of Full Scale per 24 Hours | |
| ZERO & SPAN ADJUSTMENT | Via front panel, TCP/IP or RS-232 | |
| NH₃, HCN & SO₂ EFFECT | Not detectable with 100 ppm | |
| CO₂ EFFECT | Less than 0.5% with 10% CO ₂ | |
| FLOW CONTROL | Electronic Proportional Pressure Controller | |
| SAMPLE FLOW RATE | .5 to 3.0 LPM (See footnote below) | |
| CONVERTER | Vitreous Carbon Material @ 205°C > 98% efficiency | |
| OZONATOR | Ultraviolet Lamp | |
| AIR OR O₂ REQUIREMENTS | Less than 0.01 ppm NO _x at 350 cc/Min. @ 25 psig (Dew Point < -35°C) | |
| NO/NO_x Control | Manual/Remote/Auto Cycle (Remote NO _x mode by dry contact closure) | |
| OUTPUTS AVAILABLE | TCP/IP, RS232, Four Scalable Analog 0-10 V / 4-20 mA Maximum | |
| DISCRETE ALARMS (Local & Remote Adjustable) | General Fault/ TTL Logic (Ground True) Calibration Failure/ TTL Logic (Ground True) High Concentration (2 each)/ TTL Logic (Ground True) | |
| DIGITAL DIAGNOSTICS | Control Voltages Temperatures | Pressures Flow Parameters |
| KEYPAD DISPLAYS | Factory Settings TCP/IP Address Passwords (4) | Scalable Analog Output Voltages Full Scale Range Select Auto Cal Times |
| SPECIAL FEATURES | Calculated NO ₂ derived from NO _x converter efficiency Auto Ranging Auto Calibration (adjustable through internal clock) Less than 3 cc Gold Plated Reaction Chamber | |
| DISPLAY | 3" x 5" Back lit LCD | |
| SAMPLE TEMPERATURE | Up to 50°C Noncondensing | |
| AMBIENT TEMPERATURE | 5 to 40°C | |
| AMBIENT HUMIDITY | Less than 90% RH Noncondensing | |
| WARM-UP TIME | 1 Hour (Typical) | |
| FITTINGS | 1/4 Inch Tube | |
| POWER REQUIREMENTS | 115V 60Hz (Option: 230V 50 Hz) , ± 10%, 500 W | |
| DIMENSIONS | 5¼ H × 19 W × 23 D (Inches) | |

Note: .5 to 1.5 l/min flow rate options available upon request only

SPECIFICATIONS ARE SUBJECT TO CHANGE WITHOUT NOTICE

3. Installation

3.1. General

The instrument is designed for industrial applications. These installation instructions are for a typical site. Any questions regarding specific installation situations should be directed to Technical Service of California Analytical Instruments, Inc.

3.2. Site and Mounting

NOTE: The following precautions must be carefully observed:

1. Select a site free from direct sunlight, radiation from a high temperature surface, or abrupt temperature variations.
2. This analyzer is not suitable for installation outdoors.
3. Select a site where the air is clean. Avoid exposing the instrument to corrosive or combustible gases.
4. The instrument must not be subject to severe vibration. If severe vibration is present, use isolation mounts.
5. The instrument is designed for rack-mounting. Optional rack mount slides are available.
6. Do not install near equipment emitting electromagnetic interference (EMI).


NOTE: A rear supporting brace or equivalent is required if the optional rack mount slides were not purchased.



The power on/off switch is accessible from the rear of the instrument only. DO NOT mount such that the power on/off switch is inaccessible.

3.3. Electrical

All wiring is connected at the rear of the instrument. The AC power is connected to the power/fuse/switch as shown below:

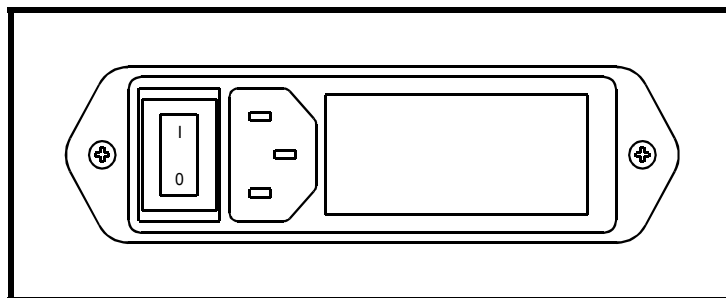

AC Power Switch, Connector, and Fuse.

NOTE: A defective ground may affect the operation of the instrument. The output voltages are connected per Table 8.1.1. Shielded wiring is recommended for output signals.

3.4. Analog Output Connections (Appendix)

See Appendix for connector pinouts located on the analyzer rear panel. Remote range identification and range selection are obtained via the rear panel connections. When a range is selected, the corresponding control line is pulled low to zero VDC. Ranges not selected will remain at approximately 5 VDC. When remote range control is selected on the front panel switch, a contact closure is provided at the rear panel connector. Remote range selection is made by connection of the control line for the desired range to the analyzers zero VDC line provided in the connector. Five VDC is also provided. Remote NOx On is selected by connection to the common line. This contact closure turns on the NOx function by flowing the sample first through the NO/NOx converter.

3.5. Gases

1. Air or O₂ (Ozone Air, < 1 ppm C) in pressurized cylinder.
2. Nitrogen or (zero air) in pressurized cylinder.
3. Standard span gas(es) near full scale concentration with a nitrogen balance, in a pressurized, certified cylinder.

3.6. Gas Handling Equipment

1. Pressure regulators for zero gas (Air or N₂), ozone supply (air or O₂) and span gas cylinders.
2. Corrosive resistant gas tubing.

NOTE

High levels of Ammonia (greater than 10 PPM NH₃) may reduce the NO₂ to NO Converter's conversion efficiency to a level that is unacceptable. It is therefore recommended that the customer purchase a commercially available NH₃ scrubber and install it in the path of the sample gas prior to its introduction into the analyzer.

3.7. Gas Connections

The tubing from the sampling system to the gas analyzer should be made from corrosive-resistant material such as Teflon or stainless steel. Even when the gases being sampled are corrosive themselves, rubber or soft vinyl tubing should not be used since readings may be inaccurate due to gas absorption into the piping material. To obtain fast response, the tube should be as short as possible. Optimum tube internal diameter is 0.16 inch (4 mm). Couplings to the instrument are ¼ Inch tube.

NOTE

Be sure tubing and joints are clean.
Dust entering the instrument may cause it to malfunction.

3.8. Sampling Requirements

3.8.1. Filtration

Dust must be eliminated completely. Use filters as necessary. The final filter must be capable of removing particles larger than 4 microns.

3.8.2. Condensation

Dew point of the sample gases must be lower than the instrument temperature to prevent accidental condensation within the instrument. Bypass the sample through a dehumidifier to reduce the dew point to about 2 to 4°C or less. If the sample contains an acid mist, use an acid mist filter, cooler or similar device to remove all traces of the mist.

3.8.3. Presence of Corrosive Gases

Useful service life of the instrument will be shortened if high concentrations of corrosive gases such as Cl₂, SO₂, F₂, HCl, etc., are present in the sampled gas.

3.8.4. Gas Temperature

When measuring high temperature gases, take care that the maximum rating of the instrument 104 °F (50 °C) is not exceeded.

3.8.5 Pressure and Flow Rates

The air or oxygen supply entering the instrument is controlled by an electronically controlled proportional flow (EPC) controller. The regulator is factory adjusted for optimum analyzer performance. The ozone supply (Air or O₂) air cylinder pressure should be set at approximately 25 PSIG. The sample entering the instrument is controlled by a factory set precision electronically controlled proportional flow (EPC) controller. The EPC is factory set for optimum analyzer performance as indicated by the sample pressure. If the analyzer does not contain the optional internal sample pump, the sample gas entering the instrument should be at a pressure between 10 and 25 PSIG with a flow capacity at a minimum of 3 liters/min. If the analyzer contains the optional sample pump, do not apply a pressurized sample. The optional pump is capable of drawing a sample through a ¼ inch heated sample line of approximately 75 feet. The calibration/span gas cylinder pressures should be set at 25 PSIG for delivery into the optional zero and span inlets located on the rear panel.

NOTE: If the analyzer contains an optional internal sample pump, the introduction of a pressurized sample gas in excess of 1.5 PSIG will damage the pump.

3.8.6. Sample Gas Bypass Outlet (Vent)


A sample gas bypass outlet connector is located on the rear panel (¼ Inch Tube). Pressure at this outlet should be kept at atmospheric level. **ANY** backpressure will cause an error in reading. The vent outlet is located on the rear panel and may contain high levels of ozone which should be vented away from the instrument.

4. Basic Operation

The operation of the digital microprocessor conforms to the guidelines of the AK committee, originally developed in the German automotive industry. Via the serial port of the MSR-Card, the analyzer can be remote-controlled by a master computer. The serial communication fully corresponds to the specifications of the AK protocol. TCP/IP communication is also available.

4.1. Display

The analyzer's LCD display includes 16 lines with 30 characters each. The display also has background lighting that can be switched on and off via the Display key on the keyboard. The following example shows the measurement screen which is formatted into 4 information areas.

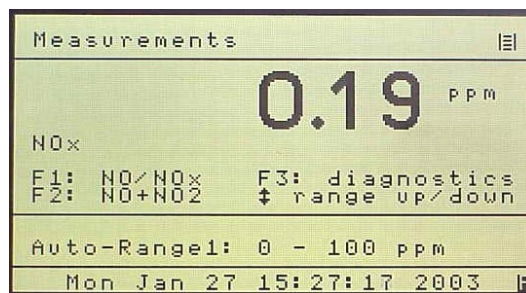

Measurement Screen

THE TOP INFORMATION AREA CONTAINS:

The AK Protocol Information. This capability is for advanced uses and may be toggled on and off in the setup screen, F5. Next to the symbol for the active operating mode, the device status is indicated. The status field is also displayed on all other screens.

- SARE Autorange enabled
- SMGA Measuring gas is flowing
- SMAN Device is in manual operation status

The level of Password Entry is shown on the right with 1 to 4 horizontal lines.

THE LARGE INFORMATION AREA CONTAINS:

The Concentration of the gas sample and mode of operation.

THE THIRD INFORMATION AREA CONTAINS:

The help information for the parameter selected, ranges, etc.


THE LOWER INFORMATION AREA CONTAINS:

The time & date and any error condition.

The symbol in the bottom right corner indicates the keyboard mode. In the example shown, the keyboard is in the function key mode. For input fields, the mode is usually switched to numerical input. Then, an N appears in the lower right of the screen. This symbol is displayed on all screens.


4.2. Menu Tree

Main Menu (from “Main” Key) (p12)

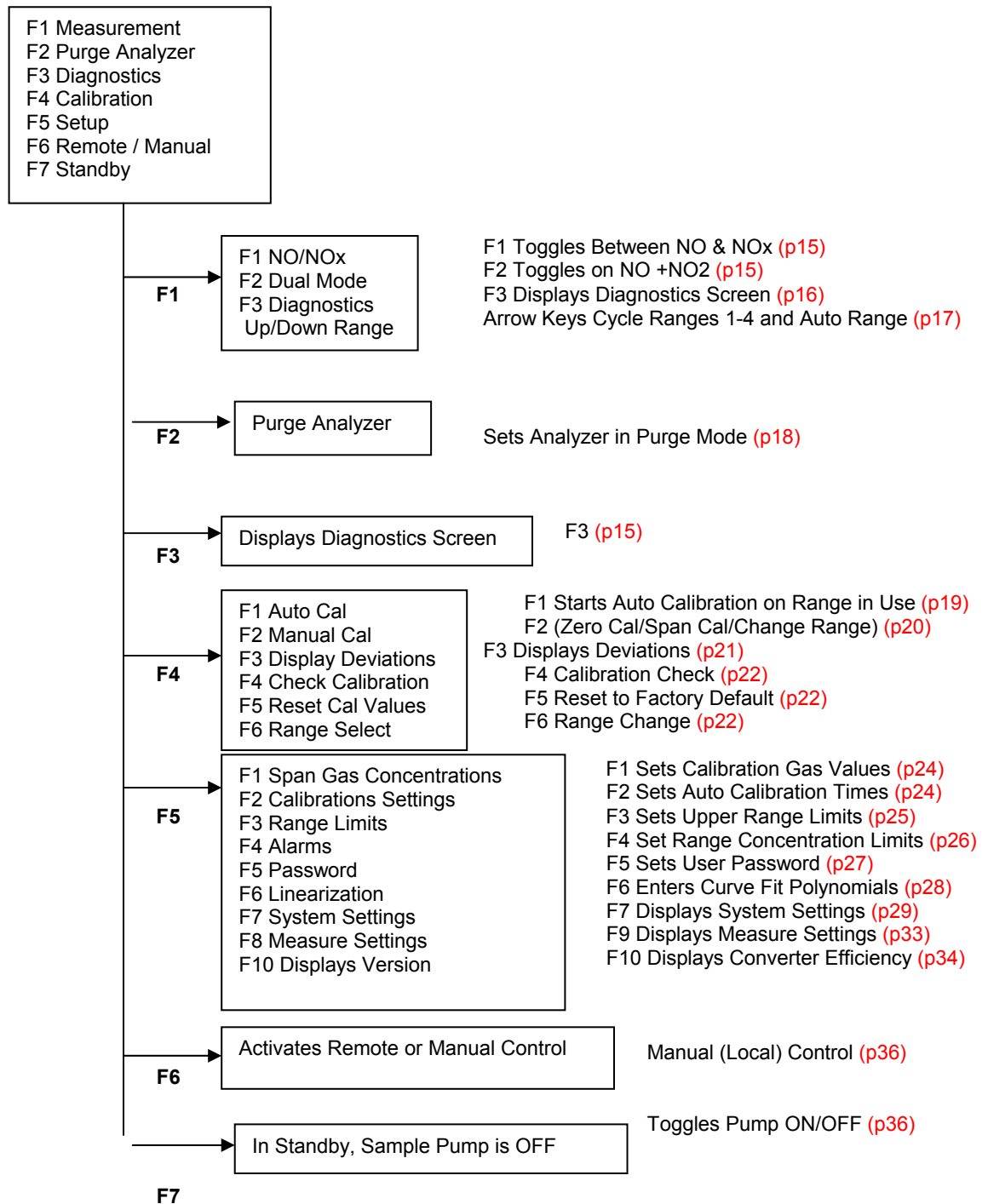

4.2. Keyboard

The keyboard looks as follows:

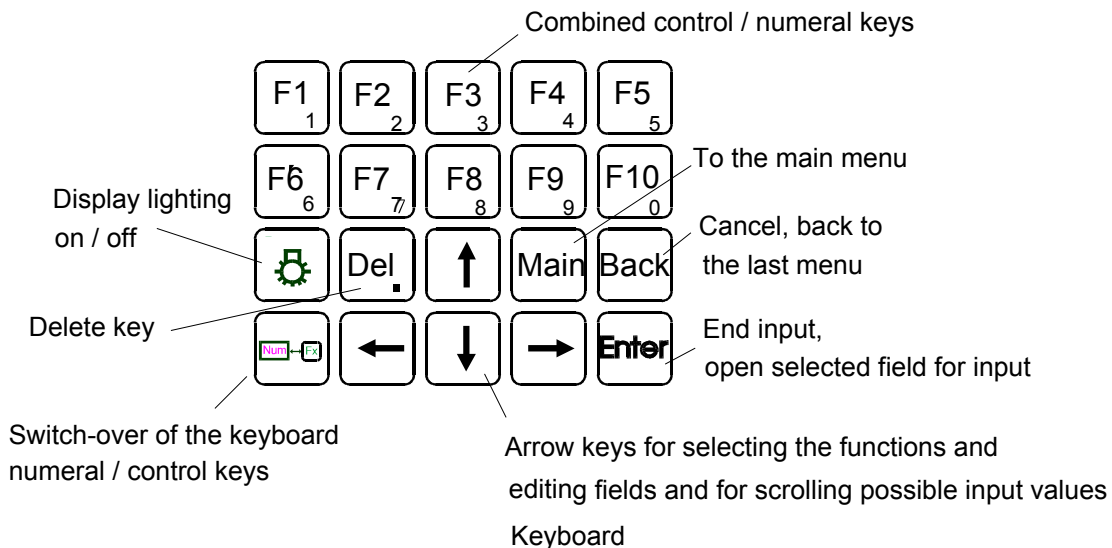

4.3. Operation with the Cursor Keys and the Enter Key

When operating the unit with the cursor keys, you select the various functions with the up/down cursor keys and start them with the Enter key. This method is particularly suitable for less proficient users since the system displays a short on-line help for nearly every function selected. The actual cursor position is shown as a black horizontal bar.

TIP: If you are not yet familiar with the screens and their fields, just press any cursor key after a screen appears. This moves the cursor from field to field and displays the corresponding online help.

4.4. Operation with the Function Keys

When using the function keys (F1 through F10), functions are directly accessed by pressing their corresponding function keys. This method is suitable for the advanced user since it is faster than the operation with the cursor keys.

4.5. Read/Change Parameters

To read and/or change parameters, you must switch to the parameter input mode by pressing the Enter key after calling the corresponding parameter screen. The input cursor (horizontal bar under the first character) then appears in the active edit field (black background). The cursor can be positioned with the right and left cursor keys, and the value displayed (number or letter) can be changed with the up and down cursor keys or entered directly. Every input has to be concluded by pressing the Enter key again, which causes the cursor to disappear.

5. Operating Structure

The analyzer's operation can be divided into 4 operating levels. The current level is always displayed as a stack of 1 to 4 horizontal bars in the top right corner of the screen. In the access level menu, you can choose between the following operating levels:

| | | |
|----|---------------|---------------------|
| F1 | User | (operating level 1) |
| F2 | Advanced user | (operating level 2) |
| F3 | Maintenance | (operating level 3) |
| F4 | System user | (operating level 4) |

A password can be assigned to each operating level. Only the system user, who normally has the highest operating priority, can assign the password. At the factory, the default passwords for the CAI analyzers are set as follows:


| | |
|----------------|-----|
| User: | 111 |
| Advanced user: | 222 |
| Maintenance: | 333 |
| System: | 444 |

The default setting can be changed only by the system user. This manual is written to include all information for the advanced system user.

TIP: Because of the user settings, some of the parameters shown in this manual may not appear on your analyzer. Check the access level.

5.1. The Main Menu

Upon power up, the CAI logo is first displayed and then the main menu appears as below:

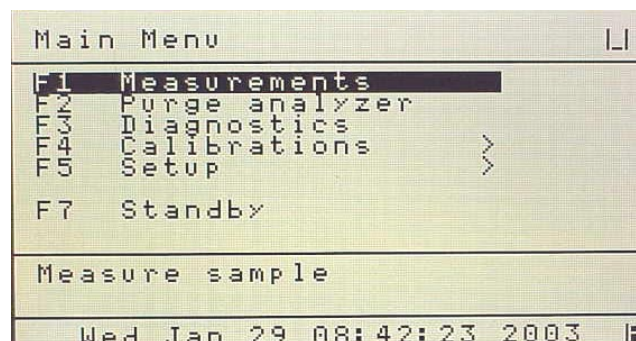

NOTE:
Access
Level
Indication

Main Menu on Power Up Screen


NOTE: F6 is not available because, on initial start up, the analyzer reverts to ONLY Level 1 access. See Section 7.5.5 for Password information.

All functions can be selected with the cursor keys and activated by pressing the Enter

key, or directly with the function keys F1 through F7. A ">" to the right of a function means that one or more sub-menus are available. If this sign is missing, the function starts immediately after the activation.

NOTE: Access level is 4.

NOTE:
Access
Level
Indication

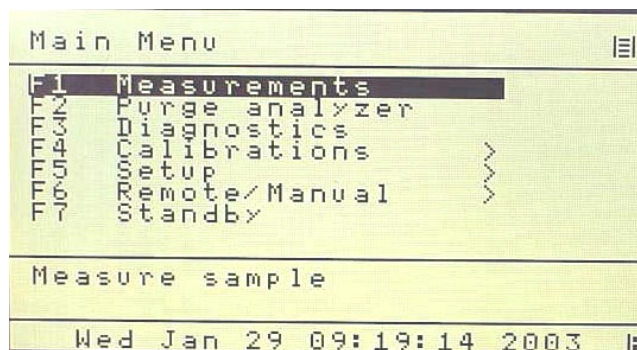

Main User Menu (Level 4)

6. Menu Structure

There are 4 operating levels based on the level of your password. This section shows the access rights of the single levels.

6.1. User Functions (Level 1)

Main Menu

F1 : Measurements
F2 : Purge Analyzer
F3 : Diagnostics
F4 : Calibrations
F5 : Setup
F7 : Standby

F5 :Setup

F5 : Password
F10:Version

F5 :Password

F1 :Enter password

6.2. Advanced User Functions (Level 2)

Main Menu

F1 : Measurements
F2 : Purge Analyzer
F3 : Diagnostics
F4 : Calibrations
F5 : Setup
F7 : Standby

F5 : Setup

F3 : Range Limits
F5 : Password
F10:Version

F5 :Password

F1 :Enter password

6.3. Maintenance Functions (Level 3)

Main Menu

F1 : Measurements
F2 : Purge Analyzer
F3 : Diagnostics
F4 : Calibrations
F5 : Setup
F7 : Standby

F5 : Setup

F1 : Span Gas Conc.
F3 : Range limits
F5 : Password
F7 : System Settings
F8 : Measure Settings
F10:Version

F5 :Password

F1 :Enter password
F2 :Reset password

F7:System Settings


F1 : Real Time Clock
F5 : Status Line on/off
F7 : Auto Startup

6.4. System User Functions (Level 4)

All functions described in this manual may be accessed from Level 4.


7. Main Menu Function Descriptions

7.1. F1 Measurements

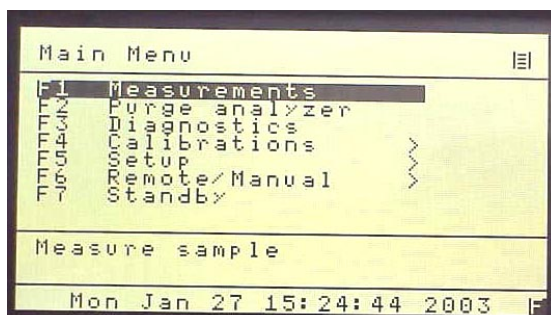

7.1.1. F1 NO or NOx Measurement

The measurements screen is activated by pressing F1 on the Main Menu screen. The NO/NOx content is displayed in ppm. Pressing F1 switches between measuring the sample gas for NOx or NO only. When the converter is off, only NO is measured. When the converter is on, NOx is measured.

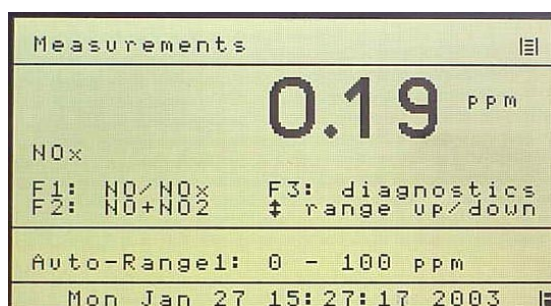

Measurements Screen

7.1.2. F2 NO + NOx Measurement

The F2 function activates the "hold and sample" feature which allows the analyzer to automatically switch between NO and NOx measurement. The time duration for the sample read is set up in the Setup Menu. The analyzer will read and display the NO (converter is bypassed) value. At the predetermined time, it will switch to the NOx mode (through converter) and read and display the NOx value, while the last 15 second NO average is displayed. The top value will be "real time" values and will change between NO and NOx. The difference between the two average values is shown as NO2. All three values are sent to the analog and digital outputs.

| Measurements | | | |
|--------------------------|---------|-----|--|
| NOx | 0 . 1 7 | ppm | |
| NO | 0 . 1 7 | ppm | |
| NO2 | 0 . 0 0 | ppm | |
| F1: NO/NOx | | | |
| Auto-Range1: 0 - 100 ppm | | | |
| Mon Jan 27 15:27:47 2003 | | | |

NO/NOx/NO2 Screen

7.1.3. F3 Diagnostics

F3 activates the diagnostic screen where pressures, flow rates, temperatures and EPC control voltages are displayed in real time. The units are psig, degrees C, ml/min. and voltage. Use the arrow key to switch between diagnostic screens.

| Diagnostics I | | | |
|----------------------------|--|--|--|
| Temperatures: Cell : 66.27 | | | |
| [°C] Diode : -5.21 | | | |
| Conv. : 205.00 | | | |
| page 1/2 | | | |
| next screen > | | | |
| Mon Jan 27 15:31:47 2003 | | | |


First Diagnostics Screen

| Diagnostics II | | | |
|--------------------------|--|--|--|
| EPC Coil: Sample: 5.76 | | | |
| [V] Ozone : 9.51 | | | |
| Pressures: Sample: 3.85 | | | |
| [PSIG] Air : -0.05 | | | |
| Gas Flow: Sample: 60.24 | | | |
| mL/Min Air : 0.00 | | | |
| page 2/2 | | | |
| < previous screen | | | |
| Mon Jan 27 15:32:03 2003 | | | |


Second Diagnostics Screen

7.1.4. Range Select

With the arrow keys, the ranges 1 to 4 can be selected and locked in which will disable the auto range capability. Continue pressing the arrow keys will recycle the analyzer back to auto range. The range and/or auto range is displayed on the measurement screen. If the limits are exceeded while not in the auto range mode, a warning "Over Range" appears on the screen.

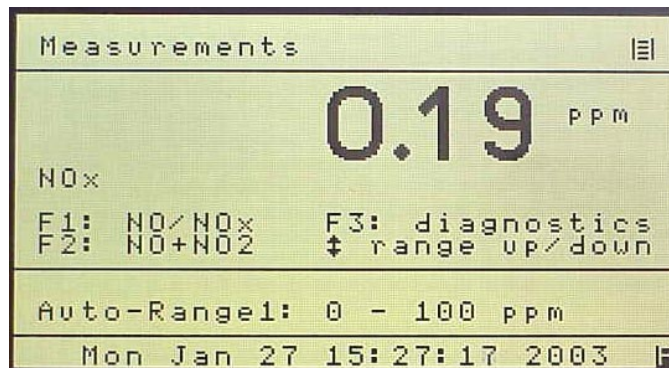


Set to Auto-Range

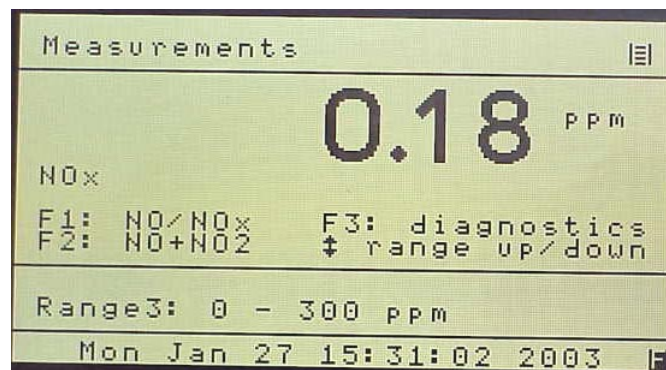


Analyzer Set to Range 3

7.2. F2 Purge Analyzer

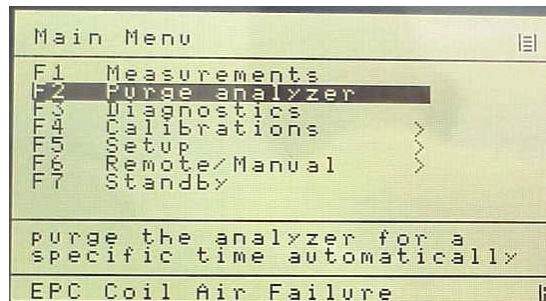

Main Menu (User Level 4)

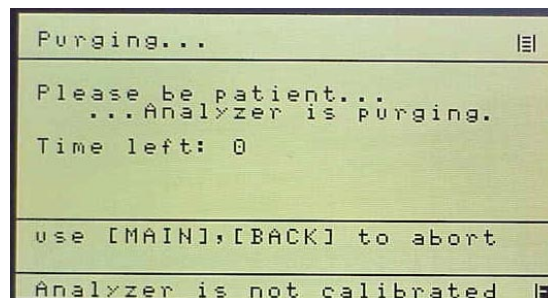


Purge Screen

F2 from the Main Menu activates the Purge (analyzer) function if equipped.

7.3. F3 Diagnostics


F3 from the Main Menu activates the Diagnostics function. As described in Section 7.1.3, F3 brings up the two diagnostics screens. The Diagnostics screens may be brought up from **EITHER** the Main Menu or the Measurements screen.

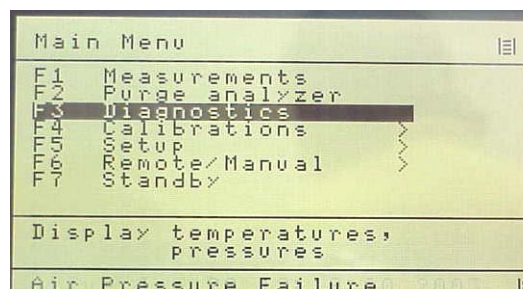


Main User Menu (Level 4)

7.4. F4 Calibrations

F4 from the Main Menu activates the Calibrations screen. Calibrations may be automatic or manual. Deviations can also be displayed. Calibration values can be reset to default values and the range to be calibrated can be changed.

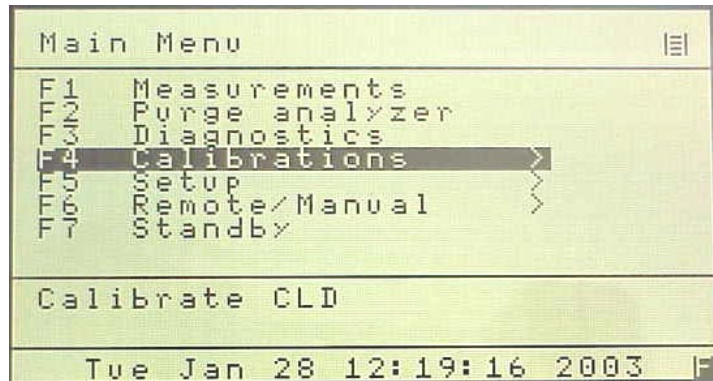

Main Menu (User Level 4)

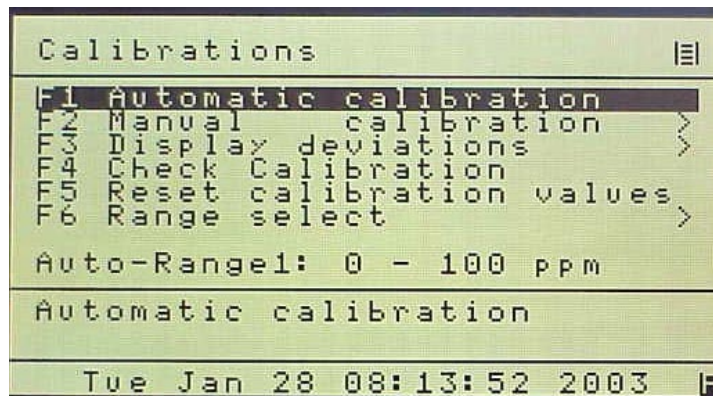


Calibration Screen

7.4.1. F1 Automatic Calibration

From the Calibrations screen, F1 starts automatic calibration. If auto range is selected, the actual range in use will be calibrated. Auto calibration works as follows: First zero gas is purged a certain time, called purge-time. Then the measurement begins. The measured value must be a minimum-time, called measuring-time and within an upper and a lower limit to be saved as new offset value. The maximum length of measuring time is 9 seconds. If the measured value was constant during calibration time, it is checked to determine if this value deviates from the preceding value. If the deviations are too large, a warning "Deviation error!" appears and the user can choose if the new value is saved or not. At last, the zero gas is flown a further time, verifying time, so it can be checked if the signal is still constant. All of these times can be changed. After zero gas calibration, the same happens with span gas. During auto calibration "Calibration in progress" is displayed. It also shows,

which gas is flowing and which time runs. When auto calibration has finished it is displayed. If the span value of the selected range is 0 (see section 5.6.1), then it will not be calibrated. If one range is calibrated and the span value for the lower ranges is zero, calibration parameters will be copied to this range. To calibrate all ranges with the same span gas, you must enter the gas concentration in the Span Gas Calibration screen for ALL RANGES. You must also calibrate each range. Offsets and scalors are NOT copied to other ranges.

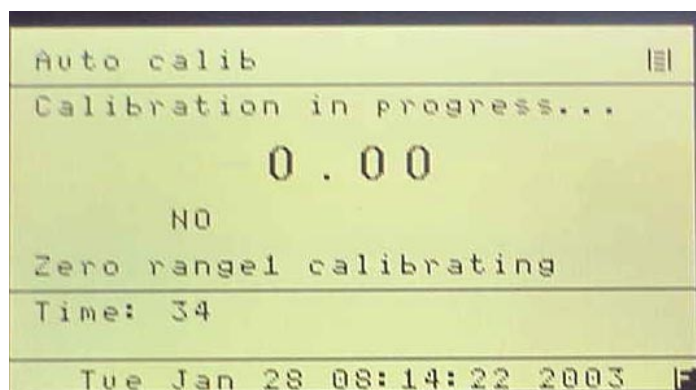


Auto Calibration Screen

7.4.2. F2 Manual Calibration



From the Calibration screen, F2 starts manual calibration. If auto range is selected, calibration is not possible, and the appropriate range can be selected. In the manual calibrations menu, three options are possible:

- F1 Flow zero gas
- F2 Flow span gas
- F3 Range select

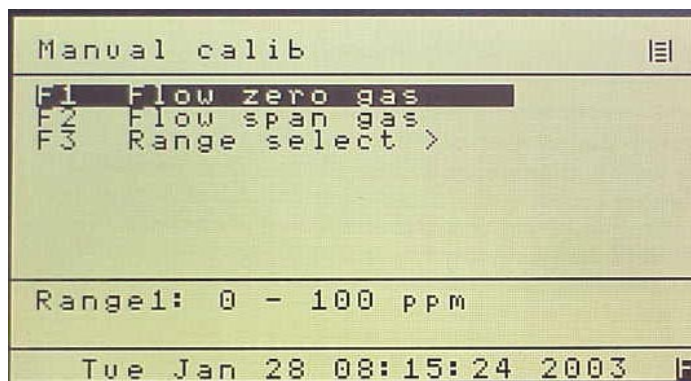

Manual Calibration Screen

When zero or span gas is flown, the measured value can be saved by pressing F1. If the screen is left by pressing the buttons "Main" or "Back", the measured value is not saved. Solenoids are closed by pressing F2. From the manual calibration menu, the range to calibrate can be chosen by pressing F3.

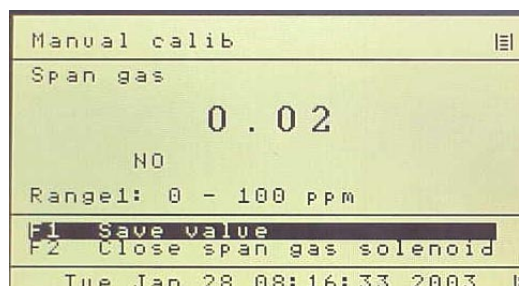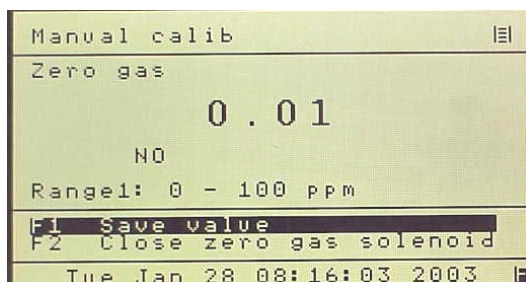


Manual Zero and Span Calibration Screens

5.4.3. F3 Display Deviations


After every calibration, the deviations are calculated for zero and for span gas.

- F1 shows zero gas deviations
- F2 shows span gas deviations
- F3 Deviations of zero gas during verifying
- F4 Deviations of span gas during verifying

F1 and F2 deviations are displayed in percent.

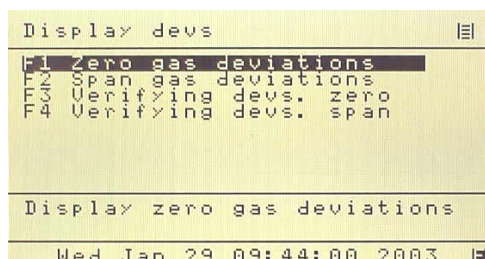

Deviation Screen

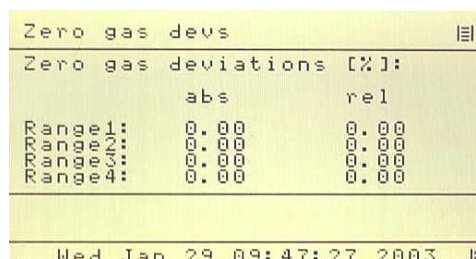

Zero Gas Deviations

During calibration there is a verification for zero and span gas. With option F3 and F4 you can view the deviations during the verification time. Absolute deviation is the absolute average difference from the saved value in ppm. Relative deviation is the absolute average difference in percent, related to the range limit.

7.4.3.1 Absolute Zero Gas Deviation

Absolute zero gas deviation is zero gas content calculated by the factory polynom related to the range limit of the calibrated range.

7.4.3.2. Relative Zero Gas Deviation

Relative zero gas deviation is the actual deviation minus the deviation of the previous calibration related to the range limit of the calibrated range.

7.4.3.3. Absolute Span Gas Deviation

Absolute span gas deviation is span gas bottle value minus span gas value calculated by the factory-polynom related to the range limit of the calibrated range.

7.4.3.4. Relative Span Gas Deviation

Relative span gas deviation is the actual deviation minus the deviation of the previous calibration related to the range limit of the calibrated range.

7.4.4. F4 Check Calibration

There is a default calibration. Pressing F4, activates an automatic zero and span check for verification.

7.4.4. F5 Reset Calibration Values


There is a default calibration. Pressing F5, a new screen appears and asks if the user is sure to reset calibration values to the default calibration values. F1 confirms and the calibration values are reset to default calibration values. F2 leaves this menu without resetting to default values. This function will overwrite all calibrations with factory values. Also the linearization polynom will be overwritten with the factory values.

7.4.5. F6 Range Select


This allows a range change to be activated from the calibration menu.

7.5. F5 Setup


From the Main Menu, F5 brings up the setup menu. Span gas concentrations, calibration settings, range limits, alarms, password, linearization, system and measure settings can be changed. The Setup menu begins as shown below. A description of each parameter is shown in the information box. NOTE: Use the down arrow key to obtain the additional setup parameters.

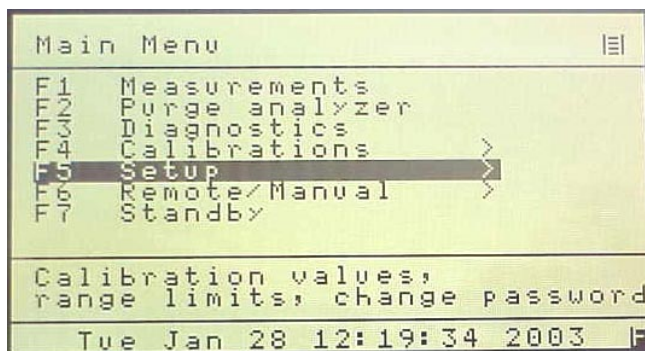

Main Menu (User Level 4)

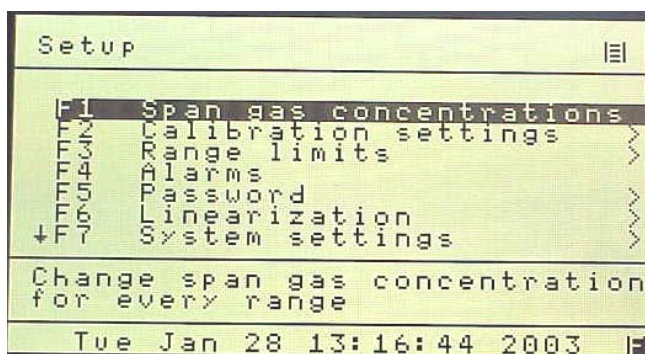


Setup Menu 1

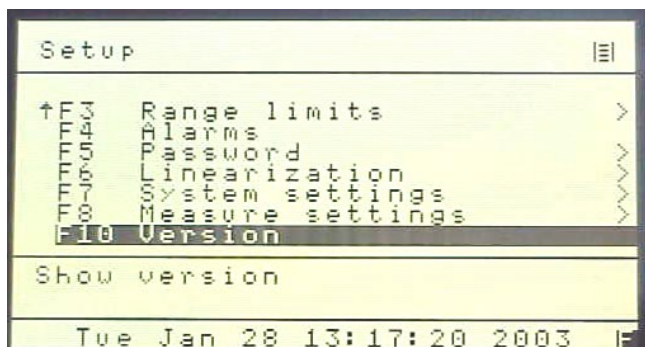

Setup Menu 2

7.5.1. F1 Span Gas Concentration

For calibration, it is necessary to input the concentration of the span gas in ppm. For every range, the span gas concentration can be changed. After pressing F1 in the setup menu, a screen appears in which changes can be made. Select with the cursor buttons the range to change. The selected field turns black. To change parameters, switch to parameter input mode by pressing the Enter key. The input cursor (horizontal bar under the first character) then appears in the active edit field (black background). The cursor can be positioned with the right and left cursor keys, and the value displayed (number or letter) can be changed with the up and down cursor keys or entered directly. Every input has to be concluded by pressing the Enter key again. Then the input cursor disappears and a new range can be selected. The changes are saved by leaving the screen by pressing "Main" or "Back". At the right side of the screen, the range limits of the 4 ranges are displayed. They cannot be changed in this screen.

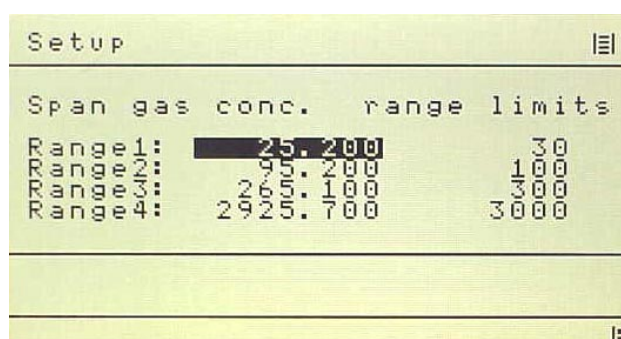

| Setup | | ≡ |
|----------------|--------------|------|
| Span gas conc. | range limits | |
| Range1: | 25.200 | 30 |
| Range2: | 95.200 | 100 |
| Range3: | 265.100 | 300 |
| Range4: | 2925.700 | 3000 |

Change Span Gas Settings

7.5.2. F2 Calibration Settings

In the calibration settings menu, times, deviations and methods can be changed.

| | | |
|--|----------------------|---|
| Setup | | ≡ |
| F1 | Times | |
| F2 | Measuring deviations | |
| F3 | Deviations | |
| F4 | Calib. via valves | |
| F5 | Calib. via probe | |
| Change purge-, measuring-, verifying-, calibration time | | |
| Wed Jan 29 11:08:11 2003 | | |

Change Auto Calibration Settings

7.5.2.1 F1 Times

There are four times (in seconds) for auto calibration that can be changed. Purge, measuring, calibration and verifying time. Changes are made and saved as above.

7.5.2.2 F2 Measuring Deviations

During auto calibration, the measured value is only saved if it is within a certain time within an upper and a lower limit. These two limits format a working window. In the setup menu the deviation is in percent.

7.5.2.3 F3 Deviations

Here you can change absolute and relative deviation in percent. After auto calibration, it is checked to assure the deviations are within this limit. If the deviations are not in this limit, a warning "Deviation error!" appears.

7.5.2.4 F4 Calibrations via Valves


Calibrations can be made by using the solenoids for zero and span gas or by using the pump. Calibration via valves means that the zero gas is flown by the zero gas solenoid and the span gas is flown by the span gas solenoid.

7.5.2.5 F5 Calibration via Probe

Calibration via probe means that the zero and the sample gas is flown by the pump, the solenoids for zero and span gas are not used.

7.5.3. F3 Range Limits


There are 4 different ranges. The user can define the upper range limits in ppm.

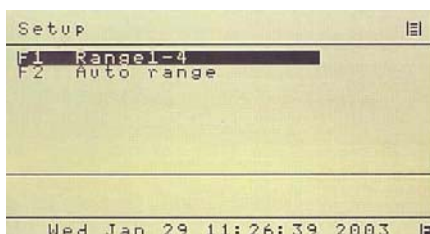

Change Range Limits

7.5.3.1 F1 Range 1-4 (Change Upper Range Limits)

In this menu the upper range limits can be changed. The new settings are saved by pressing MAIN or BACK. The auto range limits are automatically adapted. This means that if the upper range limit of range 1 for example has reached 90% of the upper range limit in the auto range mode, it is switched automatically to the second range.

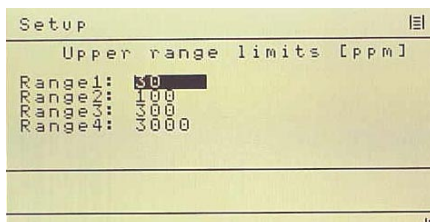

Change Upper Range Limits

7.5.3.2 F2 Change Auto Range Limits

Although the auto range limits are adapted automatically, it is possible to define them manually. Up means the value when the next higher range is selected in auto range mode, down the value when the next lower range is selected.

| Setup | | | |
|---------|------|-----|-------|
| | down | up | [ppm] |
| Range1: | | 27 | |
| Range2: | 27 | 90 | |
| Range3: | 90 | 270 | |
| Range4: | 270 | | |

Change Auto Range Limits

7.5.4. F4 Alarms

Error reports are always displayed in the lowest line of the screen. There are two pressures, three temperatures, one concentration and two voltages with alarm limits that can be defined. The user can define the range limits and, if exceeded, will display an error-message.

| Alarms I | | | | |
|--------------------------|--------|--------|--------|----------|
| | | Min | Max | |
| -Temps: | Cell : | 65.00 | 67.00 | |
| [°C] | Diode: | -5.50 | -4.50 | |
| | Conv.: | 204.00 | 207.00 | |
| F1 - next page | | | | page 1/2 |
| Tue Jan 28 14:10:36 2003 | | | | |


Set Temperature Alarms

| Alarms II | | | | |
|----------------------|--------|-------|-------|----------|
| | | Min | Max | |
| -Conc. alarms: | | 500 | 750 | |
| -Press: | Sample | 3.82 | 3.88 | |
| [PSIG] | Air | 14.00 | 16.00 | |
| -EPC: | Sample | 1.00 | 5.00 | |
| [V] | Air | 1.00 | 7.00 | |
| F1 - previous page | | | | page 2/2 |
| Air Pressure Failure | | | | |

Set Concentration, Pressure and Voltage Alarms

7.5.5. F5 Password


After turning on the analyzer, you are in access level 1. To change the access level or to change the passwords, press F5 (Setup) in the main menu and Press F5 (Password) again. The following screen appears:

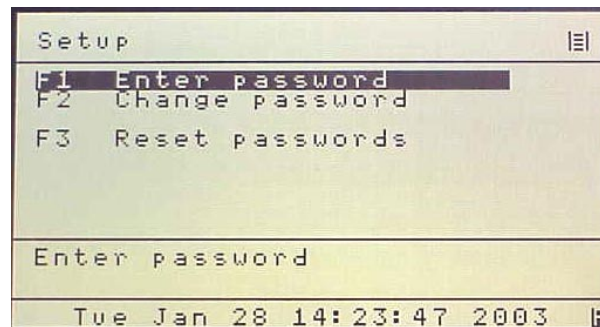

Enter / Change Password

7.5.5.1 F1 Enter Password

To change access level, press F1. The following screen appears:

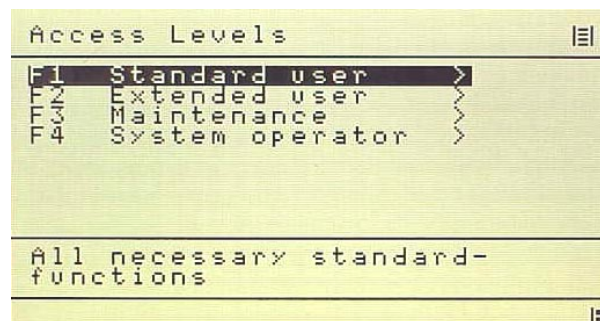

Access Level Screen

F1 to F4 selects an access level. Move the cursor to the access level to be modified. You must enter the correct password for the access level desired. The passwords for the various operation levels consist of three numbers that must be entered on the numeric keypad. If the code word is incorrect, you are asked to re-enter the codeword.

+ **IMPORTANT TIP:** When a new analyzer is powered up, it defaults to access level 1 (User). To operate ALL parameters and gain complete access, select F4. Press the Enter key twice and enter 444.

7.5.5.2 F2 Change Password

The passwords can only be changed, if you are in access level 4. After F2, enter your new 3 digit passwords.


+ **IMPORTANT TIP:** You MUST remember and record this new password. If this is lost, you will need to consult the factory for the default password !!

7.5.5.3 F3 Reset Passwords


The passwords can only be changed, if you are in access level 4. Reset passwords will revert back to the factory defaults.

7.5.6. F6 Linearization


Pressing F6 on the Setup screen brings up the Linearization screen. The analyzer can be linearized by a polynomial with 5 coefficients. By pressing F1, these 5 coefficients can be changed for each range. By pressing F2, the raw value can be displayed. This is the value before linearization and offset span correction. There are two values on the screen: The value at the top is the linearized, offset-span-corrected value, and the other value is the raw-value.

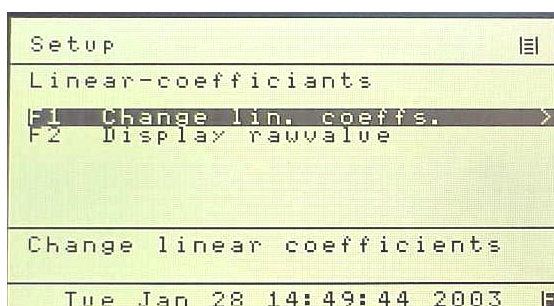


Linearization Screen Linearization

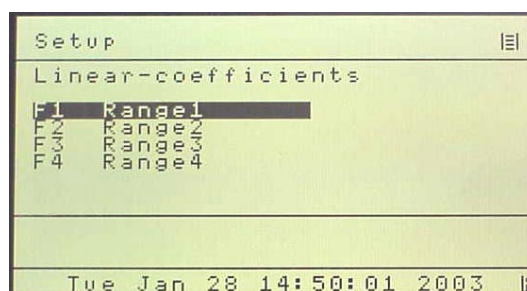

Coefficients Range Select

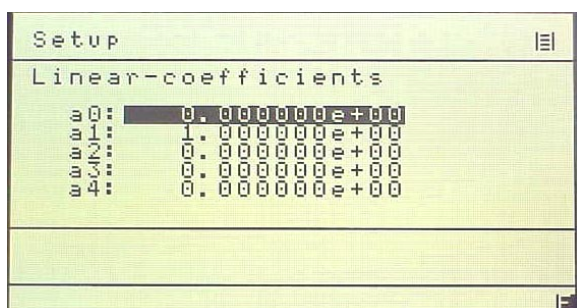


Change Linearization Coefficients
of Selected Range

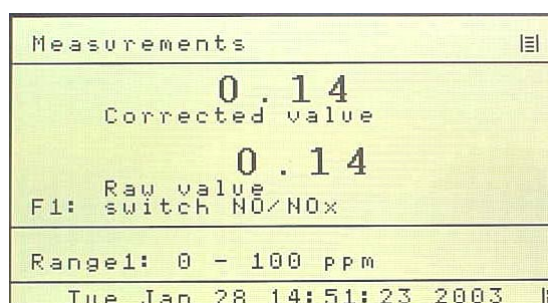

Example of Linearized and
Raw Data with F2

7.5.7. F7 System Settings


This screen allows all the system settings to be displayed and modified.

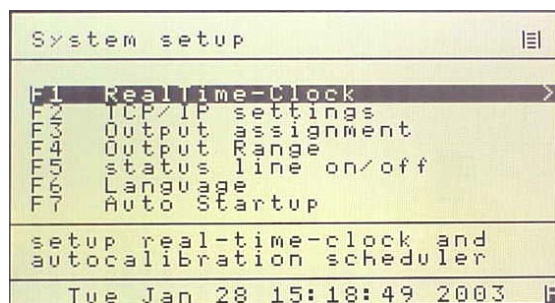

System Setup Screen


7.5.7.1 F1 Real Time Clock

This brings up the clock time set screen, auto cal and auto cal enable screens.

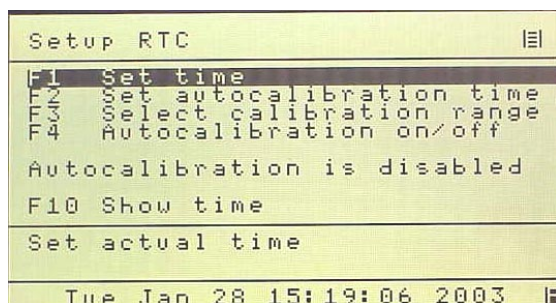

Clock and Timing Setup Screen

F1 brings up the clock set screen

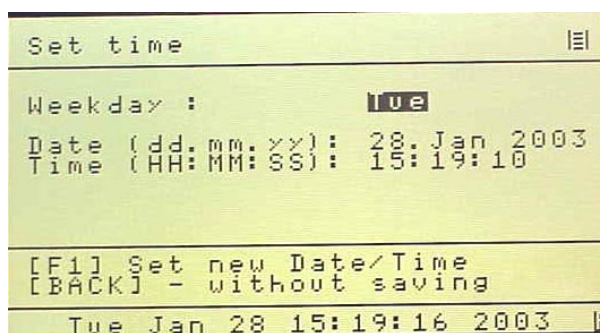

Set Clock Screen

The current time may be set by using the cursor to highlight the entry and using the numeric keys to change the values.

F2 brings up the auto cal time set. As above, the date and times can be set by using the cursor to highlight the entry and using the numeric keys to change the values. F3 Sets autocalibration ranges.

```
schedule |
Starttime : 16:00 on Tue
Date      : 28 Jan 2003
Every     : 1 hour(s).
F1 - change to weekly
F2 - change to daily
F3 - change to hourly
press MAIN or BACK to exit
F1, F2, F3 to save changes
Tue Jan 28 15:19:33 2003 |
```

Set Auto Cal Timing

```
Range selection |
Enter Range [0..4]: 4
Converter Mode   :
If 0 all ranges will be
used for autocalibration
Tue Jan 28 15:20:00 2003 |
```

Set Auto Cal Ranges

```
Setup RTC |
F1 Set time
F2 Set autocalibration time
F3 Select calibration range
F4 Autocalibration on/off
Autocalibration is disabled
F10 Show time
enable/disable
autocalibration
Air Pressure Failure |
```

F4 Toggles Auto Cal ON of OFF.

7.5.7.2 F2 Displays TCP/IP Address

```
TCP/IP setup
IP-address: 192.000.000.229
Netmask   : 255.255.255.000
Port      : 7700
WinIfPort : 2000

HWaddress : 00.E0.4B.01.9D.D4

enter IP-Address
take effect after reboot

Tue Jan 28 15:48:23 2003
```

TCP/IP Address

7.5.7.3 F3 Displays Output Signal Assignments

(Used to Adjust Analog Output Channels)

```
assignment
Output  Signal      Error
       on AOut
1      Realtime     Off
2      off
3      off
4      off

Air Pressure Failure
```

Output Assignments

7.5.7.4 F4 Displays Output Ranges

(Used to Adjust Scale of Analog Output Channels)

```
output ranges
Output  Limit      Mode
       lower  upper
1       0       0     V
2       0       0     V
3       0       0     mA
4       0       0     mA


enter 0 in both fields to
use the default range.

Air Pressure Failure
```


Output Ranges

7.5.7.5 F5 Turns Status Line On or Off

The status line displays the AK Protocol action on the top line of the display.

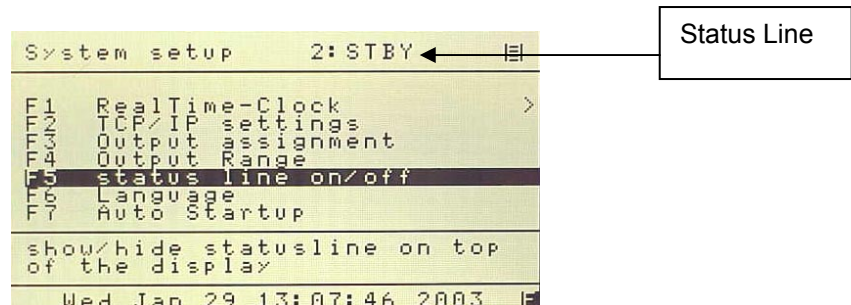


7.5.7.6 F6 Language

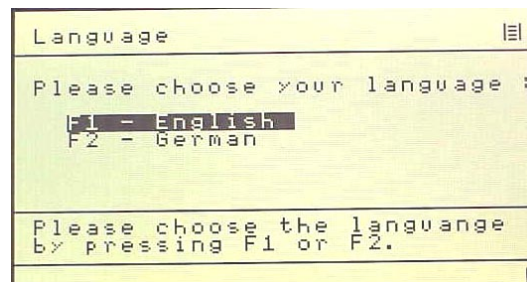

Select Language

7.5.7.7 F7 Automatic Setup


This screen brings up the automatic startup parameters. If activated, the analyzer will automatically start up the autocalibration cycle upon power on. The function is toggled on and off with the Enter key. The cycle timing, number of cals, range and NO/Nox mode may be set. After calibration, the analyzer enters the sample mode and outputs a digital signal. This is very useful in unattended applications.

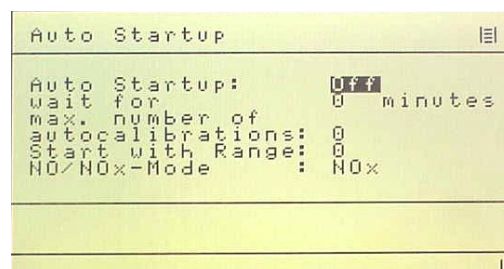


Automatic Startup Parameters

7.5.8. F8 Measure Settings

This screen allows several of the system settings to be displayed and modified.

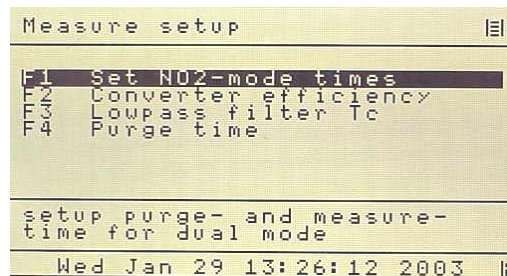


Menu Settings Screen

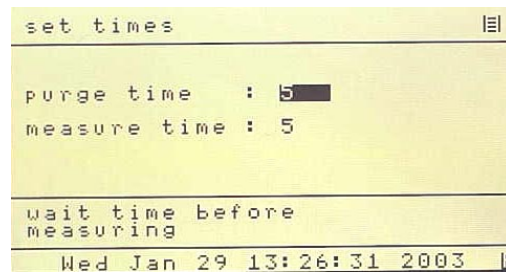7.5.8.1 F1 Set NO2 Mode Times

Set NO2 Purge and Measure Time

7.5.8.2 F2 Converter Efficiency


F2 on the Menu Settings screen allows the NO to NO2 converter efficiency to be set to the actual measured converter efficiency. A value of 100% equals 1.00. F2 will prompt the operator through the NOx efficiency test using a NOx generator.

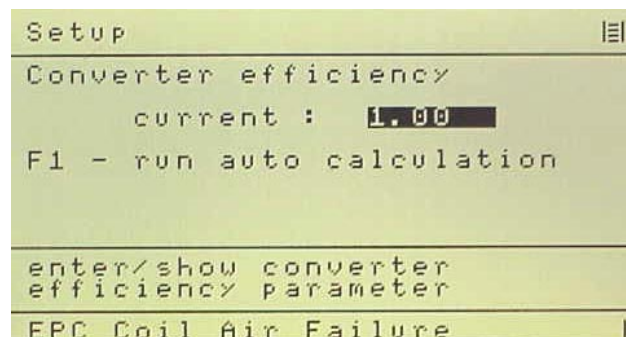

Set Converter Efficiency

7.5.8.3 F3 Low Pass Filter Time Constant


F3 on the Menu Settings screen allows the software time constant to be set between 1 and 60 seconds. This is very useful in eliminating noise when measuring low level concentrations.

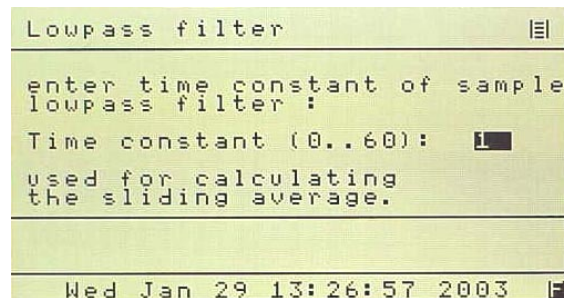

Set Time Constant

7.5.8.4 F4 Purge Time


F4 on the Menu Settings screen sets the purge time before continuing with a zero or span calibration.

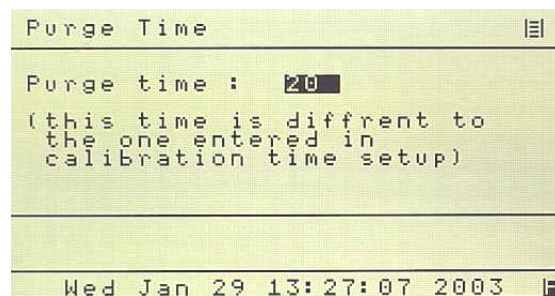


Set Purge Time

7.5.10. F10 Displays the Current Analyzer and Software Versions

This displays the analyzer's information, including the factory recommended air and sample pressure settings.

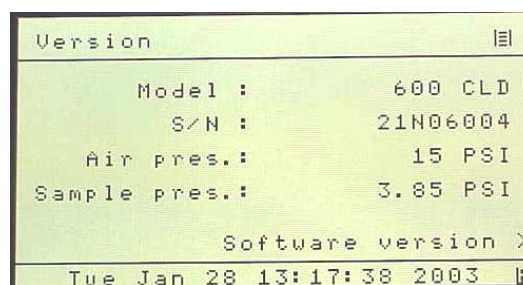


Analyzer Information Version

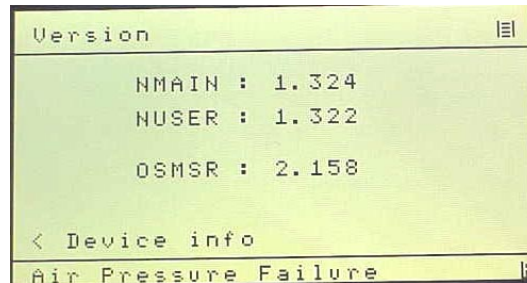

Software Version


7.6. F7 Remote / Manual Control

The analyzer can be remote-controlled by either a master computer or via contact closures. The TCP/IP and serial communication fully corresponds to the specifications of the AK protocol. To change remote/manual control, press F6 in the main menu. This toggles between remote and manual control.

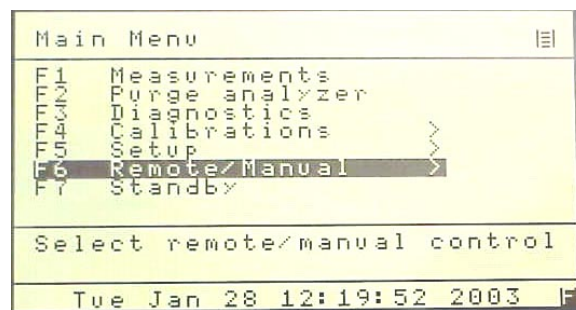

Main Menu (User Level 4)

7.7. F8 Standby

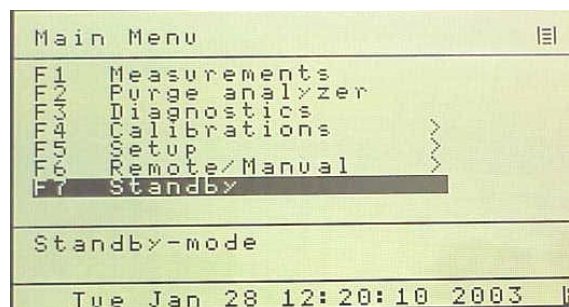

Main Menu (User Level 4)

In Standby mode, pump is turned off and the solenoids are closed. The CAI logo is displayed.

8. Analyzer Components

8.1. Rear Panel

The following details the rear panel connections:

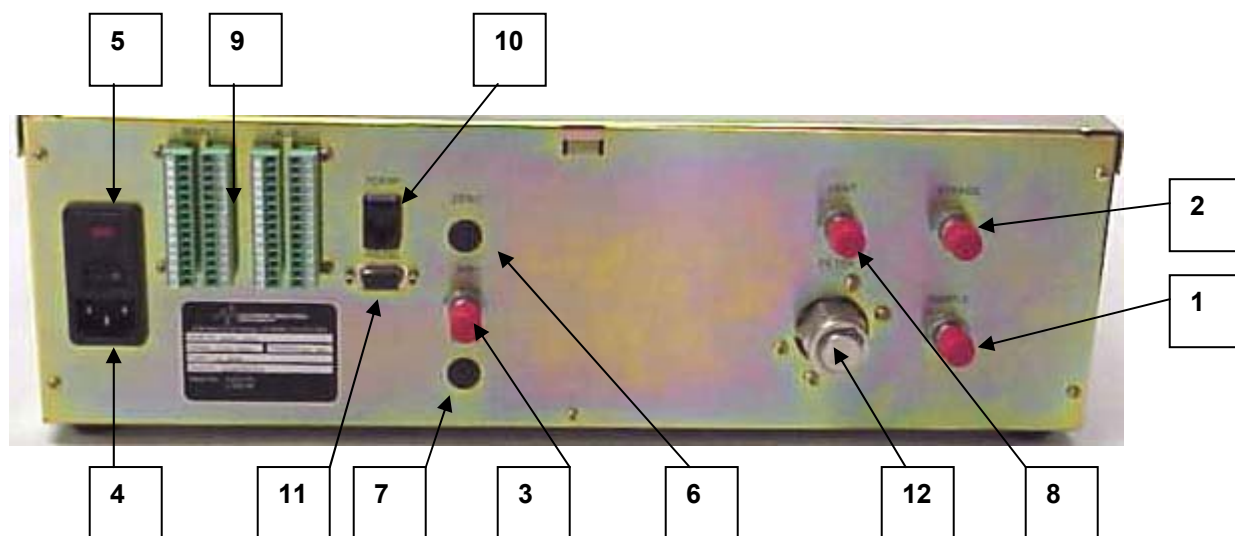


Figure 32: Rear Panel

1. **Sample Gas Inlet:** Feeds sample gas to the analyzer. ¼ Inch Tube.
2. **Sample Gas Bypass Outlet (Vent):** Exhaust for sample. ¼ Inch Tube.
3. **Ozone Air Inlet:** For feeding hydrocarbon free air or oxygen to the ozone generator.
4. **Power Entry Module:** Power connection, power switch, fuse compartment (2 Amp).
5. **Rear Panel Power ON/OFF Switch:** Turns ON/OFF line power to instrument.
6. **Zero Gas Inlet:** For feeding hydrocarbon free zero air to the analyzer.
7. **Span Gas Inlet:** For feeding calibration gas to the analyzer.
8. **Vent:** Exhaust from reaction chamber, ¼ inch tube fitting.
9. **Output Connectors:** Analog Outputs and Remote Functions.
10. **TCP/IP Connection:** Connect Network Connector.
11. **Serial Connector:** Connect Serial Connector
12. **Filter:** Analyzer Filter Housing

8.1.1. Main Connector (Standard 28 Pin Connector)

| <u>Pin</u> | <u>Signal</u> | <u>Function</u> | <u>Pin</u> | <u>Signal</u> | <u>Function</u> |
|------------|----------------|------------------|------------|----------------|-----------------|
| 1 | Analog Output | Ground (Analog) | 15 | Digital Input | Control Range 3 |
| 2 | Analog Output | Realtime | 16 | Digital Input | Control Range 4 |
| 3 | Analog Output | NO | 17 | Digital Input | Auto Cal |
| 4 | Analog Output | NOx | 18 | Digital Input | Calibrate |
| 5 | Analog Output | NO2 | 19 | Digital Input | Zero |
| 6 | Digital Output | Ground (Digital) | 20 | Digital Input | Span |
| 7 | Digital Output | Sense AutoRange | 21 | Digital Input | Pump |
| 8 | Digital Output | Sense Range 1 | 22 | Digital Input | Zero Gas Flow |
| 9 | Digital Output | Sense Range 2 | 23 | Digital Output | Span Gas Flow |
| 10 | Digital Output | Sense Range 3 | 24 | Digital Output | Sample Gas Flow |
| 11 | Digital Output | Sense Range 4 | 25 | Digital Output | Local/Remote |
| 12 | Digital Input | Set Auto Range | 26 | Digital Output | Read Cal Mode |
| 13 | Digital Input | Control Range 1 | 27 | Digital Output | Reserved |
| 13 | Digital Input | Control Range 2 | 28 | Digital Output | Reserved |

8.1.2. Auxiliary Connector (Standard 28 Pin Connector)

| <u>Pin</u> | <u>Signal</u> | <u>Function</u> | <u>Pin</u> | <u>Signal</u> | <u>Function</u> |
|------------|----------------|-------------------|------------|----------------|-------------------|
| 1 | Analog Input | Ground | 15 | Digital Output | Ground (Alarm) |
| 2 | Analog Input | External Analog 1 | 16 | Digital Output | Calibrate Alarm 1 |
| 3 | Analog Input | External Analog 2 | 17 | Digital Output | Reserved |
| 4 | Analog Input | Spare Analog | 18 | Digital Output | Reserved |
| 5 | Analog Input | Spare Analog | 19 | Digital Output | Reserved |
| 6 | Digital Output | Ground (Alarm) | 20 | Digital Output | Read Wet Mode |
| 7 | Digital Output | General Alarm | 21 | Digital Output | Read Overflow |
| 8 | Digital Output | Ch 1 Conc Alarm | 22 | Digital Output | Read NO Mode |
| 9 | Digital Output | Ch 2 Conc Alarm | 23 | Digital Input | Set Wet Mode |
| 10 | Digital Output | Reserved | 24 | Digital Input | Set Overflow Mode |
| 11 | Digital Output | Reserved | 25 | Digital Input | Set NO Mode |
| 12 | Digital Input | Reserved | 26 | DI/DO | Spare |
| 13 | Digital Input | Reserved | 27 | DI/DO | Spare |
| 13 | Digital Input | Reserved | 28 | DI/DO | Spare |

NOTE: Analog outputs 0-10 VDC Maximum and Digital outputs are 0-5 VDC Maximum. Analog inputs are 0-10 VDC Maximum.

8.1.3. Digital Outputs – RS-232 (Standard 9 Pin DIN Connector)


| <u>Pin</u> | <u>Function</u> |
|------------|-------------------------|
| 1 | DCD Carrier Detect |
| 2 | RxD Receive Data |
| 3 | TxD Transmit Data |
| 4 | DTR Data Terminal Ready |
| 5 | Ground |
| 6 | DSR Data Set Ready |
| 7 | RTS Ready to Send |
| 8 | CTS Clear to Send |
| 9 | RI Ring Indicator |

8.1.4. Digital Outputs – TCP/IP (8 Pin RJ-47 Connector)

| <u>Pin</u> | <u>Function</u> |
|------------|-----------------|
| 1 | TDX+ |
| 2 | TDX- |
| 3 | RXD+ |
| 4 | Open |
| 5 | Open |
| 6 | RXD- |
| 7 | LNLED |
| 8 | LNLED |

+ **IMPORTANT TIP:** For direct connect to a PC a crossover cable is required. Connection to a hub requires a straight cable.


8.2. Internal Component Locations

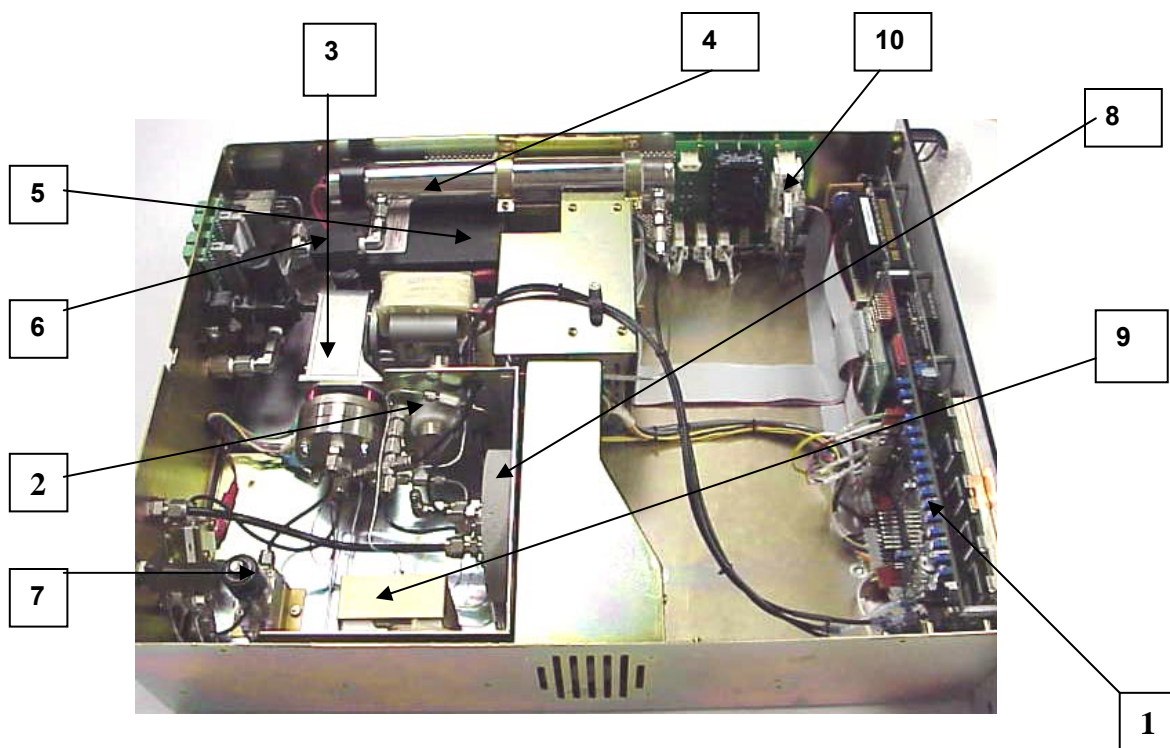

Major Internal Components

1. **Electronics:** Includes instrument electronics. (See Main Electronic Board)
2. **NO/NO_x Solenoid Valve:** Switches flow between the NO and NO_x mode.
3. **Optional Internal Sample Pump:** Provides sample to analyzer.
4. **Ozonator:** Contains UV Lamp.
5. **Ozonator High Voltage Supply:** Produces High Voltage to UV lamp.
6. **Proportional Flow Pressure Regulator:** Regulates flow of ozone.
7. **Proportional Flow Pressure Regulator:** Regulates flow of sample.
8. **Reaction Chamber & Detector Assembly:** See Figure 8.
9. **NO/NO_x Converter:** Converts NO₂ to NO for total NO_x
10. **Relay Control Board:** Provides AC Voltage to Heaters, Pump and UV Transformer.


8.3. Main Electronics Board (Potentiometers)

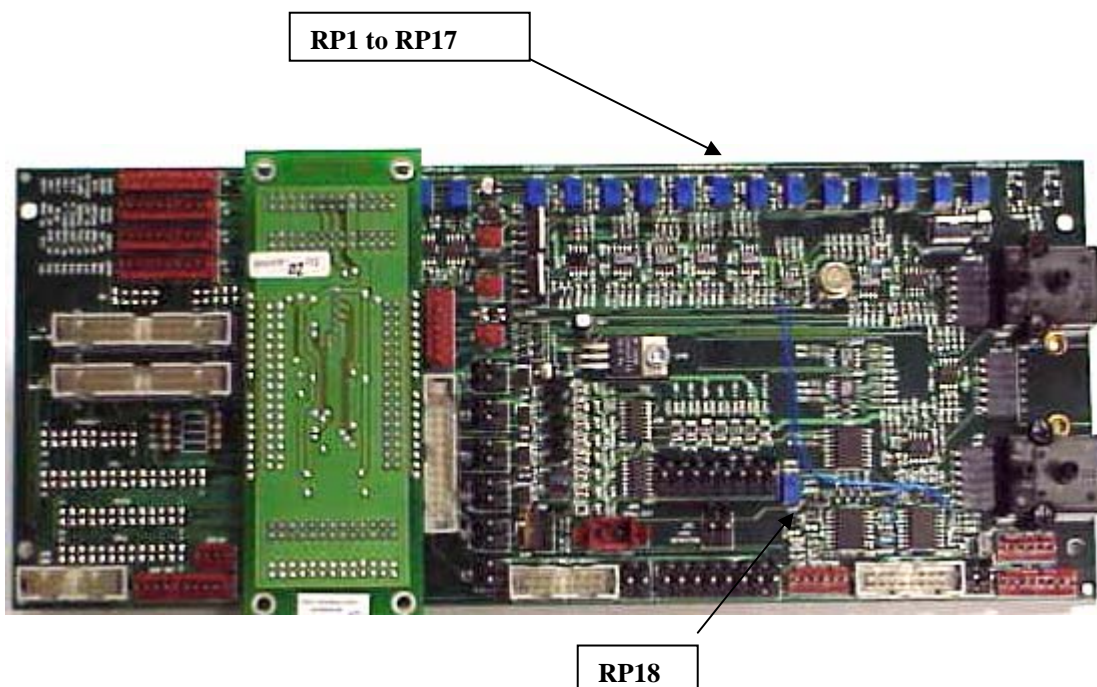

Main Electronic Board Potentiometers

| | | | |
|------------|-----------------------|-------------|-------------------------|
| RP1 | : EPC 9.5V Sample Set | RP10 | : Chiller Zero Temp Set |
| RP2 | : EPC 9.5V Air Set | RP11 | : Chiller Span Temp Set |
| RP3 | : O3 Cutoff | RP12 | : Chiller Temp Set |
| RP4 | : Cell Temp Set | RP13 | : 12VDC Adjust |
| RP5 | : Oven Temp Set | RP13 | : Sample Pressure Set |
| RP6 | : Pump Temp Set | RP15 | : Air Pressure Set |
| RP7 | : Converter Temp Set | RP16 | : Not Used |
| RP8 | : O2 Temp Set | RP17 | : Not Used |
| RP9 | : NH3 Temp Set | RP18 | : Coarse Zero Adjust |

NOTE: Potentiometers are clearly labeled on both sides of the PCB.

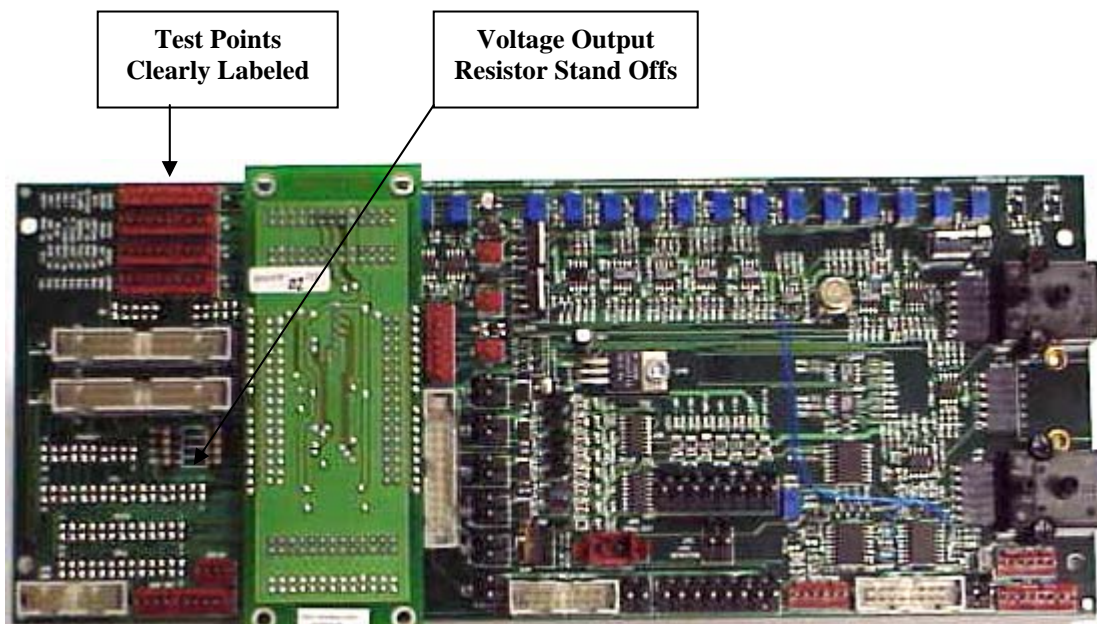8.4. Main Electronics Board (Connectors)


Main Electronic Board Connectors

| | |
|------------------------------------|---|
| J1 : Test Points | J2 : Test Points |
| J3 : Test Points | J4 : EPC Air Valve |
| J5 : Test Points | J6 : Digital Output 2 (DIDO Board) |
| J7 : EPC Sample | J8 : Sample Transducer |
| J9 : Aux Back Panel | J10 : Spare Digital Output |
| J11 : Diluter | J12 : Main Back Panel |
| J13 : Digital Input 2 | J13 : NO/NOx Valve |
| J15 : Diluter Transducer | J16 : Span Valve |
| J17 : Digital Output 1 | J18 : Zero |
| J19 : Aux Power | J20 : Air Transducer |
| J21 : Sample Overflow Valve | J22 : Daisy Chain Input 1 (DIDO Board) |
| J23 : Wet/Dry Valve | J24 : Chiller Out |
| J25 : + 5 Volt Detector | J26 : Spare Analog Input |
| J27 : Chiller Temp Sense | J28 : Spare Back Panel |
| J29 : Spare Digital Input | J30 : Daisy Chain Output (DIDO Board) |
| J31 : Fan Power | J32 : Relay Board |
| J33 : Chiller Power | J34 : Power |
| J35 : Detector | J36 : O2 Detector |
| J37 : Thermocouple | J38 : RTD |
| JP1 : PGA Zero | |


NOTE: Connections are clearly labeled on the PCB

--

8.5. Reaction Chamber

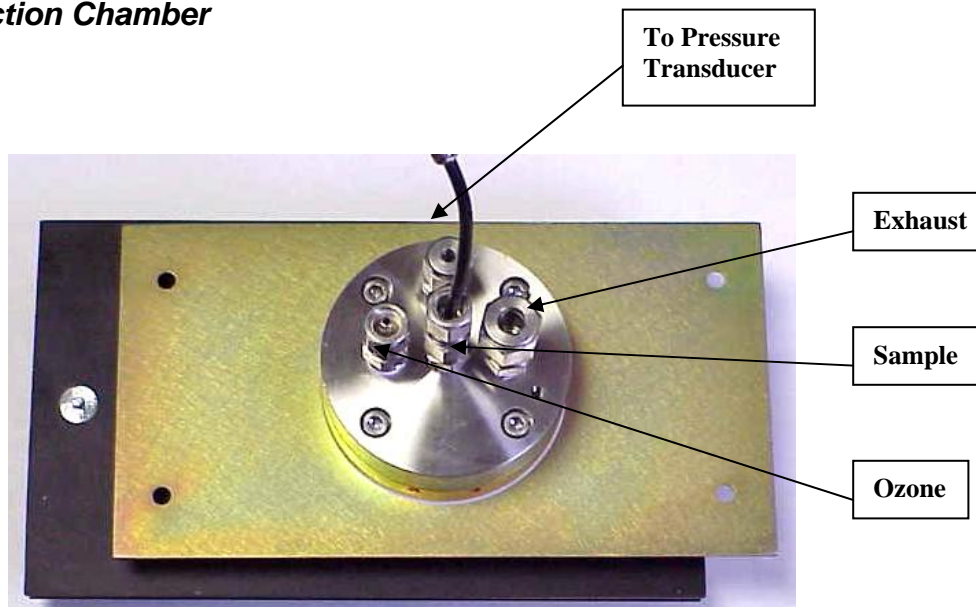


Reaction Chamber Assembly (Oven Side)

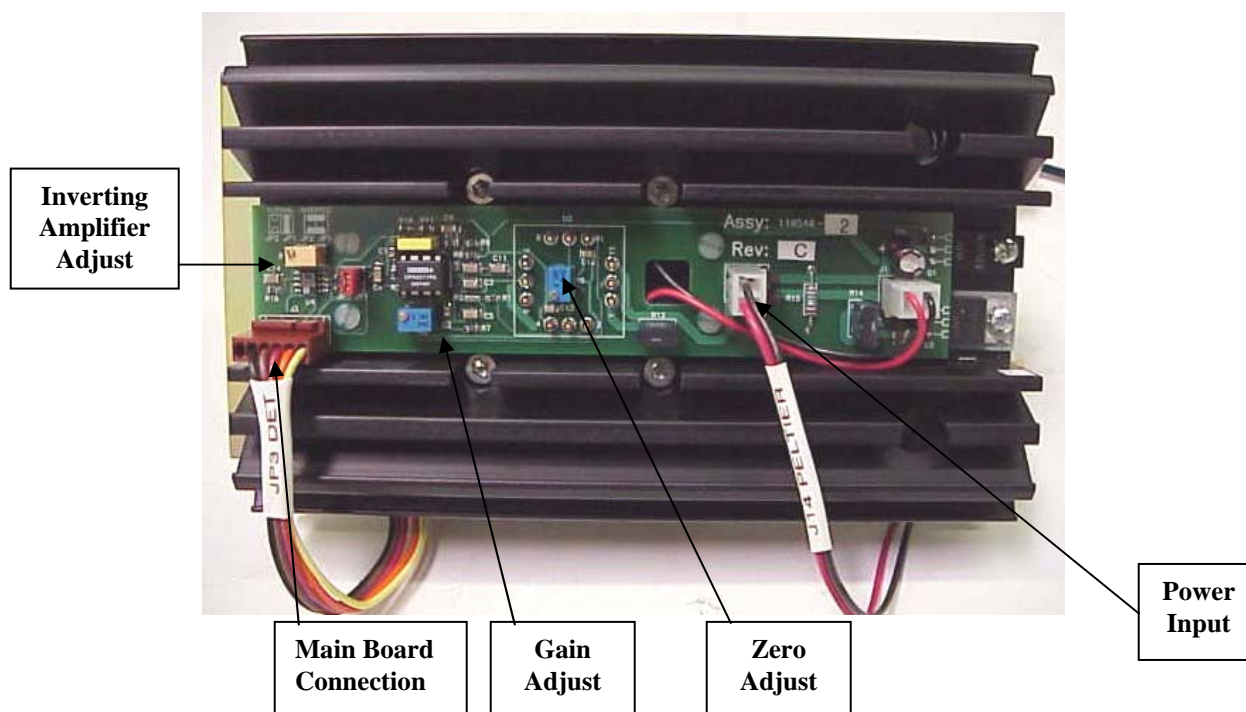

Reaction Chamber Pre-Amplifier

8.6. Relay Board Connections

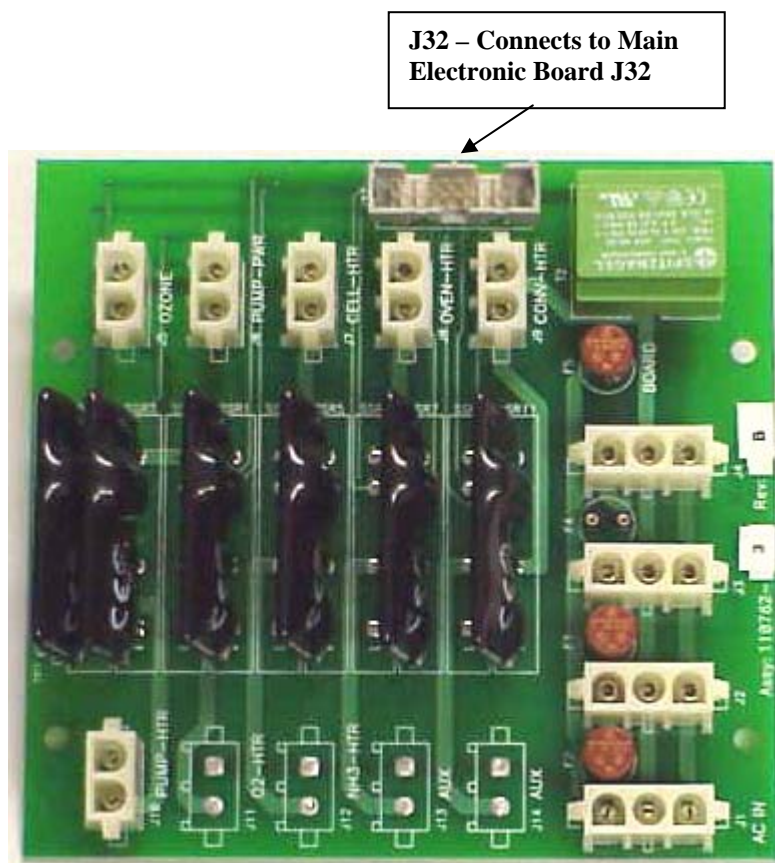

Relay Board Connections

| | |
|-----------------------------------|----------------------------------|
| J1 : AC Input | J2 : Power Supply 1 |
| J3 : Power Supply 2 | J4 : Power Supply 3 |
| J5 : Ozone Lamp | J6 : Pump Power |
| J7 : Cell Heater | J8 : Oven Heater |
| J9 : Converter Heater | J10 : Pump Heater |
| J11 : Optional O2 Heater | J12 : Optional NH3 Heater |
| J13 : Aux | J13 : Aux |

9. Operation

9.1. Preparation for Operation

Check that the external plumbing and wiring have been connected correctly, as described in this manual.

NOTE: The internal ozone generator requires approximately 1 hour of continuous operation for the analyzer to achieve full zero and span calibration stability.

9.2. Operation

1. **Power On:** Turn ON the power switch on the rear panel. The digital display should illuminate.
2. **Introduce Ozone Supply (Air or O₂):** Adjust the cylinder output pressure to 25 PSIG. The internal air pressure is factory set to deliver the air pressure required for optimum analyzer performance as indicated in the factory settings screen.
3. **Air or O₂ Pressure Settings:** Check the air pressure setting by referring to the diagnostic screen to check air pressure. The pressure should read as indicated in the factory settings screen.
4. **Zero Adjustment:** Flow zero gas through the instrument by selecting the calibration screen and select either manual or auto calibrate. **NOTE:** The instrument may also be operated by an external computer or by remote contact closures.
5. **Span Adjustment:** Flow span gas through the instrument by selecting the calibration screen and select either manual or auto calibrate.. **NOTE:** The instrument may also be operated by an external computer or by remote contact closures.

NOTE: The correct calibration gas values must be entered. The instrument is available from the factory with four ranges.

6. **NO/NO_x Function:** The analyzer switches the NO_x converter in and out of the sample stream and is controlled from the measurement screen. In the NO mode, the sample bypasses the converter and the resultant analysis produces the value of NO (Only) in the sample. In the NO_x mode, the sample passes through converter and the resultant analysis produces the value of NO_x (NO + NO₂) in the sample. The analyzer will also display the values of NO_x, NO and NO₂. The NO mode may be switched in and out remotely by a contact closure or computer. Remote control wiring is terminated in the rear panel connector. (See Appendix).

7. **Sample Pressure Check:** With sample gas flowing through the instrument, check the sample pressure setting by referring to the diagnostic screen. The sample pressure should read as indicated in the factory pressure settings screen.
8. **Sample Pump:** If the analyzer is supplied with the optional internal sample pump, it is always on in the measure mode. It is turned off during calibration and may be manually turned off by putting the analyzer in standby.
9. **Sample Line:** Make certain the sample line is flushed before connecting to the analyzer sample inlet.
10. **Instrument Power:** Turn instrument power on and allow the reaction chamber and NOx converter to stabilize before turning on the sample pump and/or connecting the sample line.
11. **Sampling System:** Prepare and check the sample system. Check the sample pressure as indicated in the factory settings screen.
12. **Air or O₂ Pressure:** Check the Air/O₂ pressure for proper setting as indicated in the factory setting screen. Readjust internal pressure as required. Note: Cylinder pressure should be set at 25 PSIG.
13. **Zero & Span Calibration:** Zero and span adjustment should be checked every 24 hours by either manual or automatic calibrations.
14. **Reaction Chamber Assembly:** Dust, water droplets, or mist entering the reaction chamber assembly may cause drift due to contamination. If the calibration procedures fails to bring the instrument to zero, check the chamber for contamination.

9.3. Shut Down Procedure

1. Turn off the zero, span and air/O₂ cylinders.
2. If the analyzer contains the optional internal sample pump, disconnect the sample line from the rear inlet port. Do **NOT** turn off the sample pump or analyzer power at this point.
3. Allow the analyzer to draw in room air for approximately 5 minutes, or flush out any remaining sample which may cause condensation as the analyzer cools.
4. Turn off the optional internal sample pump by setting the analyzer to standby.
5. Turn off the analyzer power.
6. Back-flush the sample line (and filter) of any sample before disconnecting and powering down.

10. Functional Description

10.1. Operating Principle

The California Analytical Model 600 CLD Analyzer utilizes the chemiluminescent method of determination of oxides of nitrogen (NO or NO_x) in a sample gas. In the NO mode, the NO in the sample is quantitatively converted to NO₂ by gas phase oxidation with molecular ozone produced by the UV reaction of cylinder air. Generally, 10 to 15 percent of these NO₂ molecules are elevated to an electronically-excited state. This reaction is immediately followed by reversion to a non-excited state and emission of photons. The photons impinge on a photodiode detector (PHOTODIODE) which generates a low DC current directly proportional to the NO contained in the sample gas. This current is amplified by a precision electrometer and presented to digital panel meter and recorder output. In the NO_x mode, the sample is first routed to the NO_x converter where the NO₂ component is reduced to NO. The complete sample is analyzed by the PHOTODIODE as above.

10.2. Reaction Chamber

The sample and ozone are delivered to the reaction chamber via the unique regulated flow system described below. The sample and ozone are mixed together at the center of the chamber where the reaction takes place. The sample is vented from the chamber through a 1/8 inch stainless steel tube. The chamber contains a red filter which is sealed with an integral O Ring. The chamber assembly is O Ring mounted to the PHOTODIODE. The complete chamber and PHOTODIODE assembly is housed in an RFI shielded enclosure.

10.3. Flow System

The basic function is to deliver highly regulated flows of sample and air or O₂ to the ozonator and reaction chamber assemblies. The EPC valve delivers approximately 15 PSIG to a pre-set capillary and consequently accurately predetermines the ozone flow rate. The air supply cylinder should be set to 25 psig. The sample is presented to the reaction chamber via a precision, factory set electronically controlled proportional pressure valve through a capillary. This pressure is factory set at approximately 3.85 PSIG. A close coupled bypass capillary minimizes "dead volume" and improves response time. Sample inlet pressure and regulated air pressures are monitored by internal pressure transducers and presented in PSIG via the diagnostics screen. NOTE: The correct pressures are determined by the factory for optimum analyzer performance and measured by N.I.S.T. traceable standards. They are recorded on the Factory Settings Screen.

10.4. Main Electronics Board

The main electronics board contains the instrument power supplies and required instrument electronics. A single transformer provides power to the main circuit board and includes provisions for 110/220 VAC at 50/60 Hz input.

10.5. Relay Board

The relay circuit board contains the logic circuitry required to control and switch the AC power to the required heaters and sample pump.

11. Reaction Chamber**11.1. Disassembly Procedure**

- a. Shut off ALL gas flow.
- b. Remove power from the instrument.
- c. Remove the top cover retaining screws.
- d. Remove all 4 tubes from the 4 way cross.
- e. Remove the 4 screws securing the photodiode and reaction chamber from the oven.
- f. Remove the photodiode electrical connector from the main circuit board.
- g. Remove the chiller connection from the photodiode/reaction chamber.
- h. Separate the photodiode and heat sink assembly from the reaction chamber by removing the 4 Allen screws from the front of the heat sink. Save the 2 black rubber "O" rings.
- i. Separate the mounting plate and the glass filter from the reaction chamber. Save the 2 Teflon spacers and "O" ring.
- j. Separate the manifold from the gold reaction chamber. NOTE the position of the holes in the Teflon gasket relative to the assembly screw holes. The large hole is ozone.

11.2. Assembly Procedure

- a. Wash the reaction chamber glass filter and manifold separately in detergent using a test tube brush. Be careful of the sample tube in the manifold. Do not use abrasives.
- b. Dry by blowing clean with dry nitrogen.
- c. Reassemble the chamber assembly in reverse order per the above. Make certain the sample tube is centered when assembling the manifold to the reaction chamber.

12. Troubleshooting

12.1. Ozone Air/O₂ Supply

The Air/O₂ flow is controlled by an EPC valve. It requires 25 psig cylinder supply pressure and is factory set to deliver approximately 10 to 20 psig to the ozone capillary. This pressure may be monitored by the diagnostics display. The flow rate from the capillary is very low and will require a bubble flow meter to accurately determine proper flow.

12.2. Sample Supply

The sample flow is controlled by an adjustable electronic proportional pressure valve. This valve requires a 10 to 25 PSIG sample supply pressure to deliver the proper pressure to the sample capillary. This pressure may be monitored by the diagnostics meter at any time after inlet sample has been applied. The sample flow rate from the capillary is very low and will require a bubble flow meter to determine proper flow rate. If the pressure is properly set, and a clogged capillary is suspected, replace the sample capillary.

NOTE: If the analyzer contains an optional internal sample pump, the introduction of a pressurized sample gas in excess of 1.5 PSIG will damage the pump.

12.3. NO/NO_x Converter

Several published test procedures require periodic NO_x efficiency tests to be performed on the converter to determine NO₂ to NO conversion efficiency utilizing a NO_x generator. The CAI Model NO_xGen may be used for this procedure. A short test using NO₂ calibration gas is also defined in the U.S Federal Register, Title 40, Part 86.332.79 (e).

13. Drawings

13.1. AK Protocol

13.2. Rear Panel Connections

13.3. Flow Diagram

13.4. Block Electrical Drawing

13.1 Serial Interface and AK-Commands

The serial interface enables remote control of the Model 600 analyzer by a master computer. It is implemented as an RS232 V24 interface and meets all requirements of the AK protocol.

A 9-pin male connector at the back of the unit is used to connect a master computer with the following pin assignment:

Pin 3 = Txd (transmit)

Pin 2 = Rxd (receive)

Pin 5 = Gnd (ground)

Interface Parameters

Baud rate: 9600, 4800, 2400, 1200, 600, 300 baud

Data bits: 7 or 8

Stop bit: 1 or 2

Don't care: 1 byte, adjustable (e.g. 32)

Parity: Even, odd, none

XON/XOFF: Active or not active

General AK Requirements

- 1) If the command message contains no error, the acknowledge message contains the echo of the function code and the error status number (1 to 9).
- 2) If the transfer was faulty or the function code unknown, the answer contains four question marks (example. "???? 0").
- 3) If the displayed value is not valid, a "#" is placed in front of the measured value (example: "AIKG 0 #9999").
- 4) If a control or adjusting command is sent via the serial interface while the measuring device is in "Manual" mode, it sends an answer like "SLIN 0 K0 OF".
- 5) If a channel does not exist, the answer for control and adjusting commands is e.g. "ATEM 0 3 NA" in which 3 is the number of the sub-channel.
- 6) If the device is busy with a running function (SLIN, for example), every arriving control command is ignored (except SRES and STBY); and the response message is e.g. "SMAN 0 BS. If in the mode "SINT" an additional "SINT KO" command is received, the integrator is reset to 0 and the integration is restarted.
- 7) If the command message contains data that the measuring device cannot process ("ESYZ K0 ABC", for example), the response message is "ESYZ 0 SE". A syntax error is recognized if the data does not match the expected format or if the parameters do not fit the expected size.
- 8) Numbers are in floating-point format with decimal point. The decimal point can be dropped for integers.
- 9) If you switch from "Manual" to "Remote" at the device, it remains in "Manual" mode until a "SREM K0" is received by the control computer. On the display, this mode is indicated by REME" (Remote enable) on the status line. In manual mode, query commands via the serial interface are possible at any time.

AK Protocol Format

The master computer and the Model 600 analyzer communicates via the RS232 serial link. The Model 600 analyzer acts as a "slave" and only responds to commands.

Serial Interface Parameters:

- 1) Baud from 300 to 9600 bps, can be selected via the display.
- 2) 7 or 8 data bits, 1 or 2 stop bits, and the parity (yes/no).
- 3) The data transmission is full duplex (no echo) with XON/XOFF protocol.
- 4) The "don't-care" byte" (byte 2) is adjustable (factory setting 20H).

Command Format:

<STX> 02H Example: ASTZ K0
 don't care any byte (default 20H)
 function code code 4 byte long (e.g., ASTZ)
 space 20H 20H
 channel N° always "K0" for the analyzer
 space 20H (only if followed by data, otherwise <ETX>)
 data data bytes (depending on the command)
 <ETX> 03H

Answer Format:

<STX> 02H Example: STZ 0 SREM STBY
 don't care adjustable, factory setting 20H
 function code same code as command package (e.g., ASTZ)
 space 20H
 status 0 without error or 1 to 9 when error (see also ASTF command)
 space 20H (only if followed by data, otherwise <ETX>)
 data parameter (depending on the command)
 <ETX> 03H

Scans

AKON: Measured concentration value

| Command | Response | Description |
|----------|-------------------------|--|
| _AKON_K0 | _AKON_s_z.z_y.y_x.x_w.w | Measured concentration value is responded z.z:current Measured Value y.y:NO x.x: NO2 w.w:Nox y.y,x.x,w.w are only used in dual measure mode. Otherwise "O.O" will be returned |

AEMB: Set measuring range

| Command | Response | Description |
|----------|------------|--------------------------------------|
| _AEMB_K0 | _AEMB_s_Mn | Current measuring range is responded |

AMBE: Measuring range limit

| Command | Response | Description |
|-------------|-------------------------------------|---|
| _AMBE_K0 | _AMBE_s_M1_w.w_M2_x.x_M3_y.y_M4_z.z | All existing measuring range limits are responded |
| _AMBE_K0_Mn | _AMBE_s_Mn_z.z | Range limit of Range Mn is responded |

AKAK: Calibration gas concentrations

| Command | Response | Description |
|------------|-------------------------------------|---|
| AKAK_K0 | _AKAK_s_M1_w.w_M2_x.x_M3_y.y_M4_z.z | All existing calibration gas values are responded |
| AKAK_K0_Mn | AKAK_s_Mn_z.z | Calibration gas value of Range Mn is responded |

AMBU: Upper and lower range switchover values for autorange

| Command | Response | Description |
|----------|---|---|
| _AMBU_K0 | _AMBU_s_M1_w.w_W.W_M2_x.x_X.X_M3_y.y_Y.Y_M4_z.z_Z.Z | Lower and upper range switchover value of autorange are responded |

ASTZ: Normal device status

| Command | Response | Description |
|----------|-----------------------------------|----------------------------|
| _ASTZ_K0 | _ASTZ_s_SREM_STBY__SENO_SARE_SDRY | Device status is responded |

Possible states:

| | | | | |
|-----------------|---------------------------------------|-------------------|------------------------|----------------------|
| SREM: remote | STBY: standby | SENO: NO mode | SARE: Autorange on | SDRY: Chiller on |
| SMAN: manual | SPAU: pause | SMAN: NOx mode | SARA: Autorange off | SWET: Chiller off |
| | SMGA: measuring gas | | | |
| | SNGA: zero gas | | | |
| | SEGA: end gas | | | |
| | SATK SNGA: zero gas during autocal | | | |
| | SATK SEGA: end gas during autocal | | | |
| | SLIN: For compatibility only | | | |
| | SSPL: purging | | | |
| | SKOP: measure | | | |

ASTF: Error status

| Command | Response | Description |
|----------|--------------------------|-----------------------------------|
| _ASTF_K0 | _ASTF_s_f1_f2_f3_..._f15 | Current error number is responded |

Errors:

| | |
|----|-------------------------------|
| 1 | Sample Pressure Failure |
| 2 | Air Pressure Failure |
| 3 | Oven Temp Failure |
| 4 | Converter Temp Failure |
| 5 | Pump Temp Failure |
| 6 | Diode Temp Failure |
| 7 | Cell Temp Failure |
| 8 | Peltier Gas Temp Failure |
| 9 | Reaction Chamber Temp Failure |
| 10 | EPC Coil Sample Failure |
| 11 | EPC Coil Air Failure |
| 12 | Range Overflow |
| 13 | ADC Range Overflow |
| 13 | ADC Range Underflow |
| 15 | Range 1 is not calibrated |
| 16 | Range 2 is not calibrated |
| 17 | Range 3 is not calibrated |
| 18 | Range 4 is not calibrated |

AKEN: Device identification

| Command | Response | Description |
|----------|------------------------|------------------------------------|
| _AKEN_K0 | _AKEN_s_devicename | Device identification is responded |
| _AKEN_K1 | _AKEN_s_model | Device model |
| _AKEN_K2 | _AKEN_s_serialno | Device serial number |
| _AKEN_K3 | _AKEN_s_airpressure | Suggested input air pressure |
| _AKEN_K4 | _AKEN_s_samplepressure | Suggested input sample pressure |

ARMU: Rawvalue

| Command | Response | Description |
|----------|-------------|--|
| _ARMU_K0 | _ARMU_s_z.z | Raw value before linearization and offset-span-correction is responded |

ATEM: Temperatures

| Command | Response | Description |
|------------|---------------------|---|
| _ATEM_K0 | _ATEM_s_z.z_y.y_... | All Temperatures in degrees celsius are responded |
| _ATEM_K0_x | _ATEM_s_z.z | Temperature of x in degrees celsius is responded |

Description of x:

| X | CLD |
|---|-----------------------|
| 1 | Oven Temp |
| 2 | Converter Temp |
| 3 | Pump Temp |
| 4 | Diode Temp |
| 5 | Cell Temp |
| 6 | Peltier Temp |
| 7 | Reaction Chamber Temp |

ADRU: Pressures

| Command | Response | Description |
|------------|---------------------|-----------------------------|
| _ADRU_K0 | _ADRU_s_z.z_y.y_... | All pressures are responded |
| _ADRU_K0_x | _ADRU_s_z.z | Pressure of x is responded |

Description of x:

| | |
|---|----------------------------|
| 1 | Sample Pressure |
| 2 | Air Pressure |
| 3 | Sample EPC Coil Voltage |
| 4 | Air/Ozone EPC Coil Voltage |

ADUF: Flows

| Command | Response | Description |
|------------|---------------------|-------------------------|
| _ADUF_K0 | _ADRU_s_z.z_y.y_... | All flows are responded |
| _ADUF_K0_x | _ADRU_s_z.z | Flow of x is responded |

Description of x:

| | |
|---|-------------|
| 1 | Sample Flow |
| 2 | Air Flow |

AGRD: Polynom coefficients

| Command | Response | Description |
|-------------|---------------------------|--|
| _AGRD_K0_Mn | _AGRD_s_Mn_a0_a1_a2_a3_a4 | Polynom coefficients of range Mn are responded |

AANG: Deviation from zero point after autocalibration

| Command | Response | Description |
|----------|---|---|
| _AANG_K0 | _AANG_s_M1_z.z_da_dr_ M2_z.z_da_dr_ M3_z.z_da_dr_ M4_z.z_da_dr_ | Deviation from zero point after autocalibration |

AAEG: Deviation from end point after autocalibration

| Command | Response | Description |
|----------|---|--|
| _AAEG_K0 | _AANG_s_M1_z.z_da_dr_ M2_z.z_da_dr_ M3_z.z_da_dr_ M4_z.z_da_dr_ | Deviation from end point after autocalibration |

AFDA: Purge and Autocalibration times

| Command | Response | Description |
|---------------|---------------------|---|
| _AFDA_K0_SATK | _AFDA_s_z_y_x_w_Z.Z | Autocalibration times: z: Purge time y: Calibration time x: Total Calibration time w: Verify time (z,y,x,w in seconds) |
| _AFDAKO_SSPL | AFDA_s_z.z | Purge time will be responded |

APAR: Request Autocalibration tolerance values

| Command | Response | Description |
|---------------|-------------------------|---|
| _APAR_K0_SATK | _APAR_s_z.z_y.y_x.x_w.w | Autocalibration tolerance value(%): z.z: Range 1 y.y: Range 2 x.x: Range 3 w.w: Range 4 |

AKAL: Deviations from calibration

| Command | Response | Description |
|-----------|--|--|
| _AKAL_K0_ | _AKAL_s_M1_z.z_y.y_x.x_w.w _AKAL_s_M2_z.z_y.y_x.x_w.w _AKAL_s_M3_z.z_y.y_x.x_w.w _AKAL_s_M4_z.z_y.y_x.x_w.w | Deviation (ppm): z.z: Zero gas relative to last calibration y.y: Zero gas factory calibration x.x: Span gas relative to last calibration w.w: Span gas factory calibration |

ASYZ: Respond System Time

| Command | Response | Description |
|-----------|-----------------------|---|
| _ASYZ_K0_ | _ASYZ_s_yymmdd_hhmmss | Respond system time yymmdd:year, month,day (each 2 characters wide, no spaces) hhmmss:hour,minutes,seconds) |

AT90: Respond Lowpass filter time

| Command | Response | Description |
|-----------|-----------|---|
| _AT90_K0_ | _AT90_s_t | Respond lowpass filter time t=filter time in seconds |

ADAL:Diagnostic alarm limits

| Command | Response | Description |
|------------|----------------------------------|--------------------------|
| _ADAL_K0 | _ADAL_s_a1.min_a1.max_.._f12.max | All alarms are responded |
| _ADAL_K0_x | _ADAL_s_x.min_x.max | Alarm limits of x |

Alarm Limits:

| | |
|----|----------------------------|
| 1 | Sample Pressure |
| 2 | Air Pressure |
| 3 | Oven Temp |
| 4 | Converter Temp |
| 5 | Pump Temp |
| 6 | Diode Temp |
| 7 | Cell Temp |
| 8 | Peltier Gas Temp |
| 9 | EPC Coil Sample Voltage |
| 10 | EPC Coil Air/Ozone Voltage |
| 11 | Reserved |
| 12 | Sample Content |

ATCP: Query TCP/IP settings

| Command | Response | Description |
|----------|--|---|
| _ATCP_K0 | _ATCP_s_ddd.ddd.ddd.ddd _ATCP_s_yyy.yyy.yyy.yyy _ATCP_s_xxxx | ddd: TCP/IP Address yyy: TCP/IP subnet mask xxxx: TCP/IP port |

AENT: Query calibration gas flow setting

| Command | Response | Description |
|----------|-----------|---|
| _AENT_K0 | _AENT_s_x | x=10: Calibration through sample gas inlet (pump) y=11: Calibration through zero/span valves |

Control commands**SRES:** Reset

| Command | Response | Description |
|----------|----------|-------------|
| _SRES_K0 | _SRES_s | Reset |

SPAU: Pause

| Command | Response | Description |
|----------|----------|-------------|
| _SPAU_K0 | _SPAU_s | Pause mode |

STBY: Standby

| Command | Response | Description |
|----------|----------|--------------|
| _STBY_K0 | _STBY_s | Standby mode |

SNGA: Open valve for zero gas calibration

| Command | Response | Description |
|-------------|----------|---|
| _SNGA_K0 | _SNGA_s | Open valve for zero gas calibration of actual measuring range |
| _SNGA_K0_Mn | _SNGA_s | Open valve for zero gas calibration of range Mn |

SEGA: Open valve for end gas calibration

| Command | Response | Description |
|-------------|----------|--|
| _SEGA_K0 | _SEGA_s | Open valve for end gas calibration of actual measuring range |
| _SEGA_K0_Mn | _SEGA_s | Open valve for end gas calibration of range Mn |

SSPL: Purge Analyzer with zero gas

| Command | Response | Description |
|----------|----------|--|
| _SSPL_K0 | _SSPL_s | Open valve for zero gas and purge the analyzer |

SLIN: Linearization mode

| Command | Response | Description |
|----------|----------|--|
| _SLIN_K0 | _SLIN_s | Change status to SLIN (only for compatibility) |

SKOP: Converter Check

| Command | Response | Description |
|----------|----------|---|
| _SKOP_K0 | _SKOP_s | Change status to SKOP and activate sample pump (only for compatibility) |

SWET: Chiller off – Wet mode measuring

| Command | Response | Description |
|----------|----------|--------------------|
| _SWET_K0 | _SWET_s | Switch chiller off |

SDRY: Chiller on – Dry mode measuring

| Command | Response | Description |
|----------|----------|-------------------|
| _SDRY_K0 | _SDRY_s | Switch chiller on |

SATK: Start automatic calibration

| Command | Response | Description |
|-------------|----------|--|
| _SATK_K0 | _SATK_ | Start automatic calibration of all ranges |
| _SATK_K0_Mn | _SATK_s | Start automatic calibration using range Mn |

SEMB: Set measuring range

| Command | Response | Description |
|-------------|----------|--|
| _SEMB_K0_Mn | _SEMB_s | Set measuring range Autorange is disabled |

SARE: Auto range on

| Command | Response | Description |
|----------|----------|-------------------|
| _SARE_K0 | SARE_s | Set auto range on |

SARA: Auto range off

| Command | Response | Description |
|----------|----------|-------------------|
| _SARA_K0 | _SARA_s | Set autorange off |

SREM: Remote mode for AK-commands

| Command | Response | Description |
|----------|----------|---------------------------|
| _SREM_K0 | _SREM_s | Set device in remote mode |

SMAN: Manual control to control device manually

| Command | Response | Description |
|----------|----------|---------------------------|
| _SMAN_K0 | _SMAN_s | Set device in manual mode |

SMGA: Start measuring

| Command | Response | Description |
|----------|----------|--|
| _SMGA_K0 | _SMGA_s | Start measuring Turn on pump for sample gas |

SNKA: Saves measured value as new offset.

| Command | Response | Description |
|----------|----------|---|
| _SNKA_K0 | _SNKA_s | Saves measured value of actual range as new offset if zero valve is opened |

SEKA: Saves measured value as new span value

| Command | Response | Description |
|----------|----------|--|
| _SEKA_K0 | _SEKA_s | Saves measured value of actual range as new span value if span valve is opened |

SENO: Converter off

| Command | Response | Description |
|----------|----------|--|
| _SENO_K0 | _SENO_s | Set converter off Only NO is measured |

SNOX: Converter on

| Command | Response | Description |
|----------|----------|---|
| _SNOX_K0 | _SNOX_s | Set converter on All kinds of NOx are measured |

SNO2: Converter on

| Command | Response | Description |
|----------|----------|--|
| _SNO2_K0 | _SNO2_s | Activates dual measure mode. Analyzer switches periodically between NO and NOx mode and displays NO, NO2, NOx |

SFGR: Reset calibration settings to factory defaultConverter on

| Command | Response | Description |
|----------|----------|---|
| _SFGR_K0 | _SFGR_s | Reset all calibration settings to their factory settings |

SENT: Set calibration gas flow

| Command | Response | Description |
|------------|----------|---|
| _SENT_K0_x | _SENT_s | x=10: Calibration through sample gas inlet (pump) y=11: Calibration through zero/span valves |

Settings**EKAK:** The four span gas concentration values are set

| Command | Response | Description |
|--------------------------------------|----------|--------------------|
| _EKAK_K0_M1_w.w_M2_x.x_M3_y.y_M4_z.z | _EKAK_s | Set end gas values |

EMBE: The four measuring range end values are set

| Command | Response | Description |
|--------------------------------------|----------|------------------|
| _EMBE_K0_M1_w.w_M2_x.x_M3_y.y_M4_z.z | _EMBE_s | Set range limits |

EMBU: The upper and the lower range switchover for autorange are set

| Command | Response | Description |
|--|----------|--|
| _EMBU_K0_M1_w.w_W.W_M2_x.x_X.X_M3_y.y_Y.Y_M4_z.z_Z.Z | _EMBU_s | Set lower and upper range switchover limits |

EKEN: Set new device identification

| Command | Response | Description |
|--------------------------|----------|--|
| _EKEN_K0_new device-name | _EKEN_s | Set new device identification Maximum length of device name are 40 characters |

NOTE: To change device identification, you must first rename the device to "RESET".
Now a name up to 40 letters can be given.

NOTE: The device name must not have any blanks between, f.e. "CAI CLD" is not allowed. You can use underline, i.e. "CAI_CLD".

EGRD: Set polynom coefficients

| Command | Response | Description |
|----------------------------|----------|--------------------------------------|
| _EGRD_K0_Mn_a0_a1_a2_a3_a4 | _EGRD_s | Set polynom coefficients of range Mn |

EFDA: Set autocalibration and purge times

| Command | Response | Description |
|-----------------------|----------|--|
| _EFDA_K0_SATK_z_y_x_w | _EFDA_s | Set autocalibration times: z= Purge time y=Calibration time x=Total calibration time w=Verify time (z,y,x,w in seconds) |
| _EFDA_K0_SSPL_z | _EFDA_s | Set analyzer purge time to z seconds |

EPAR: Set autocalibration tolerance values

| Command | Response | Description |
|-----------------------------------|----------|--|
| _EPAR_K0_SATK_z.z_y.y_x.x_w. w | _EPAR_s | Autocalibration Tolerance value (%): z.z= Range 1 y.y= Range 2 x.x= Range 3 w.w= Range 4 |

ESYZ: Set System Time

| Command | Response | Description |
|------------------------|----------|--|
| _ESYZ_K0_yymmdd_hhmmss | _ESYA_s | Respond system time: yymmdd:year, month,day (each 2 characters wide, no spaces) hhmmss:hour,minutes,seconds) |

ET90: Set Lowpass Filter Time

| Command | Response | Description |
|------------|----------|---|
| _ET90_K0_t | _ET90_s | Set lowpass filter time: t= filter time in seconds |

EDAL:Diagnostic alarm limits

| Command | Response | Description |
|-----------------------------------|----------|-----------------------|
| _EDAL_K0_a1.min_a1.mas_..._a12max | _EDAL_s | Set all alarm limits |
| _EDAL_K0_x_x.min_xmax | _EDAL_s | Set alarm limits of x |

Alarm Limits:

| | |
|----|----------------------------|
| 1 | Sample Pressure |
| 2 | Air Pressure |
| 3 | Oven Temp |
| 4 | Converter Temp |
| 5 | Pump Temp |
| 6 | Diode Temp |
| 7 | Cell Temp |
| 8 | Peltier Gas Temp |
| 9 | EPC Coil Sample Voltage |
| 10 | EPC Coil Air/Ozone Voltage |
| 11 | Reserved |
| 12 | Sample Content |

ETCP: Set TCP/IP Parameters

| Command | Response | Description |
|--------------------------|----------|---|
| _ETCP_K0_zzz.zzz.zzz.zzz | _ETCP_s | zzz= TCP/IP address |
| _ETCP_K0_yyy.yyy.yyy.yyy | | yyy= TCP/IP subnet mask |
| _ETCP_K0_xxxx | | xxxx= TCP/IP port |
| | | All changes take effect after next power on cycle |

Abbreviations used

| | |
|-------------|---|
| Mn | : Measuring range number |
| M1 .. M4 | : Measuring Range 1 .. 4 |
| w.w .. Z.Z. | : Numerical value |
| x | : Number |
| t | : Numeric integer value |
| a0 .. a4 | : Polynom coefficients |
| s | : Status |
| yyymmdd | : Date of format year, month and day with 2 characters each and no spaces |
| hhmmss | : Time of format hour, minute and second with 2 characters each and no spaces |

13.2 Rear Panel Connectors

28 Pin Main Connector Assignments:


| Signal Type | Analog |
|-------------|---------------------------------------|
| | pin # |
| A Output | 1 GND (analog) Realtime 0-10 VDC |
| A Output | 2 Maximum |
| A Output | 3 NO 0-10 VDC Maximum NOx 0-10 VDC |
| A Output | 4 Maximum NO2 0-10 VDC |
| A Output | 5 Maximum |
| | Digital |
| D Output | 6 GND (Digital) |
| D Output | 7 Sense AutoRange |
| D Output | 8 Sense Range 1 |
| D Output | 9 Sense Range 2 |
| D Output | 10 Sense Range 3 |
| D Output | 11 Sense Range 4 |
| D Input | 12 Set AutoRange |
| D Input | 13 Control Range 1 |
| D Input | 13 Control Range 2 |
| D Input | 15 Control Range 3 |
| D Input | 16 Control Range 4 |
| D Input | 17 Auto Cal |
| D Input | 18 Calibrate |
| D Input | 19 Zero |
| D Input | 20 Span |
| D Input | 21 Pump |
| D Output | 23 Span Gas Flow |
| D Output | 24 Sample Gas Flow |
| D Output | 25 Local/Remote |
| D Output | 26 Read Cal Mode |
| D Output | 27 Reserved |
| D Output | 28 Reserved |

28 Pin Auxiliary Connector Assignments:


| Signal Type | Analog |
|-------------|------------------------|
| | pin # |
| Spare | |
| A Input | 1 GND (analog) |
| A Input | 2 External Analog 1 |
| A Input | 3 External Analog 2 |
| A Input | 4 Spare analog |
| A Input | 5 Spare analog |
| | Alarms |
| | Digital |
| D Output | 6 GND (Alarm) |
| D Output | 7 General Alarm |
| D Output | 8 Ch. 1 Conc. 1 Alarm |
| D Output | 9 Ch. 1 Conc. 2 Alarm |
| D Output | 10 Reserved |
| D Output | 11 Reserved |
| D Output | 12 Reserved |
| D Output | 13 Reserved |
| D Output | 13 Reserved |
| D Output | 15 GND (Alarm) |
| D Output | 16 Calibration Alarm 1 |
| D Output | 17 Reserved |
| D Output | 18 Reserved |
| D Output | 19 Reserved |
| D Output | 20 Read Wet Mode |
| D Output | 21 Read OverFlow |
| D Input | 23 Set Wet Mode |
| D Input | 24 Set OverFlow |
| D Input | 25 Set NO Mode |
| DI/DO | 26 Spare |
| DI/DO | 27 Spare |
| DI/DO | 28 Spare |




Note: All digital outputs and inputs are 0-5 VDC ONLY. All analog inputs are 0-10 VDC ONLY. Connecting analog outputs to existing current loop systems or voltage loop systems ***WILL DAMAGE*** the instrument.

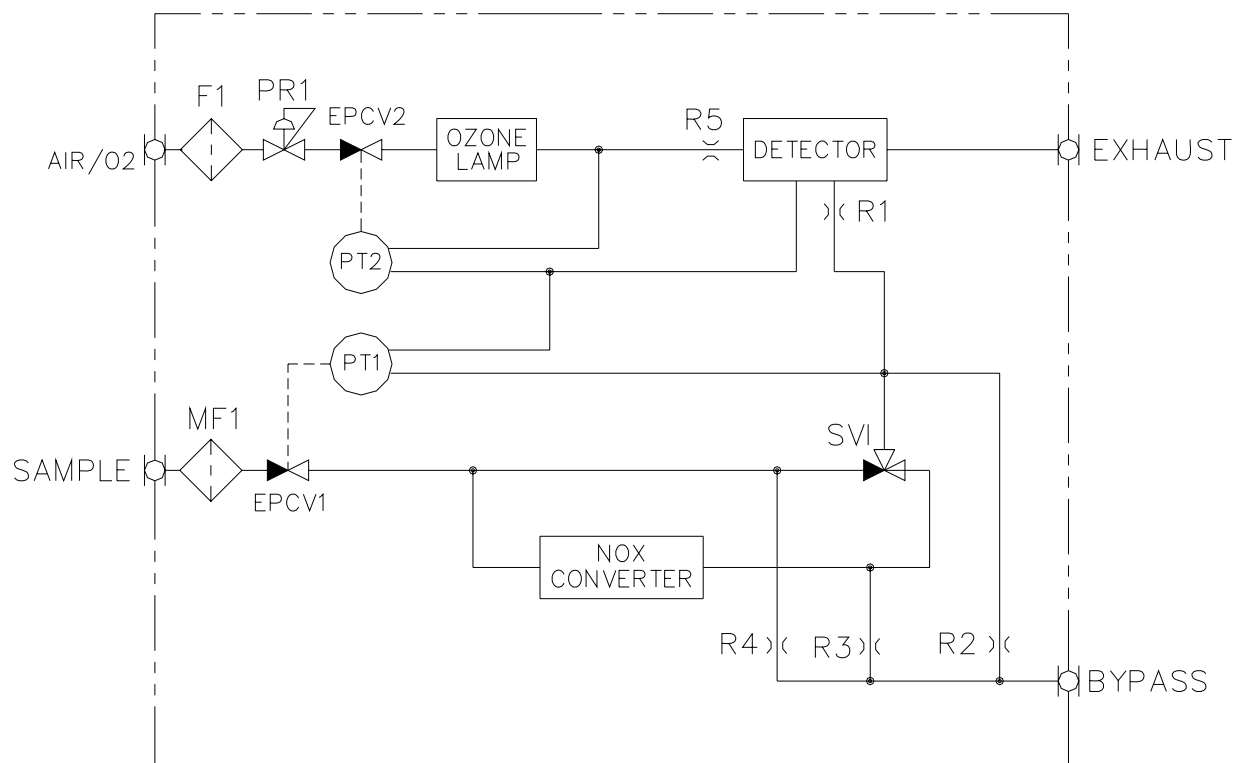13.3 MODEL 600 CLD Flow Diagrams


NOTE: Standard Analyzer

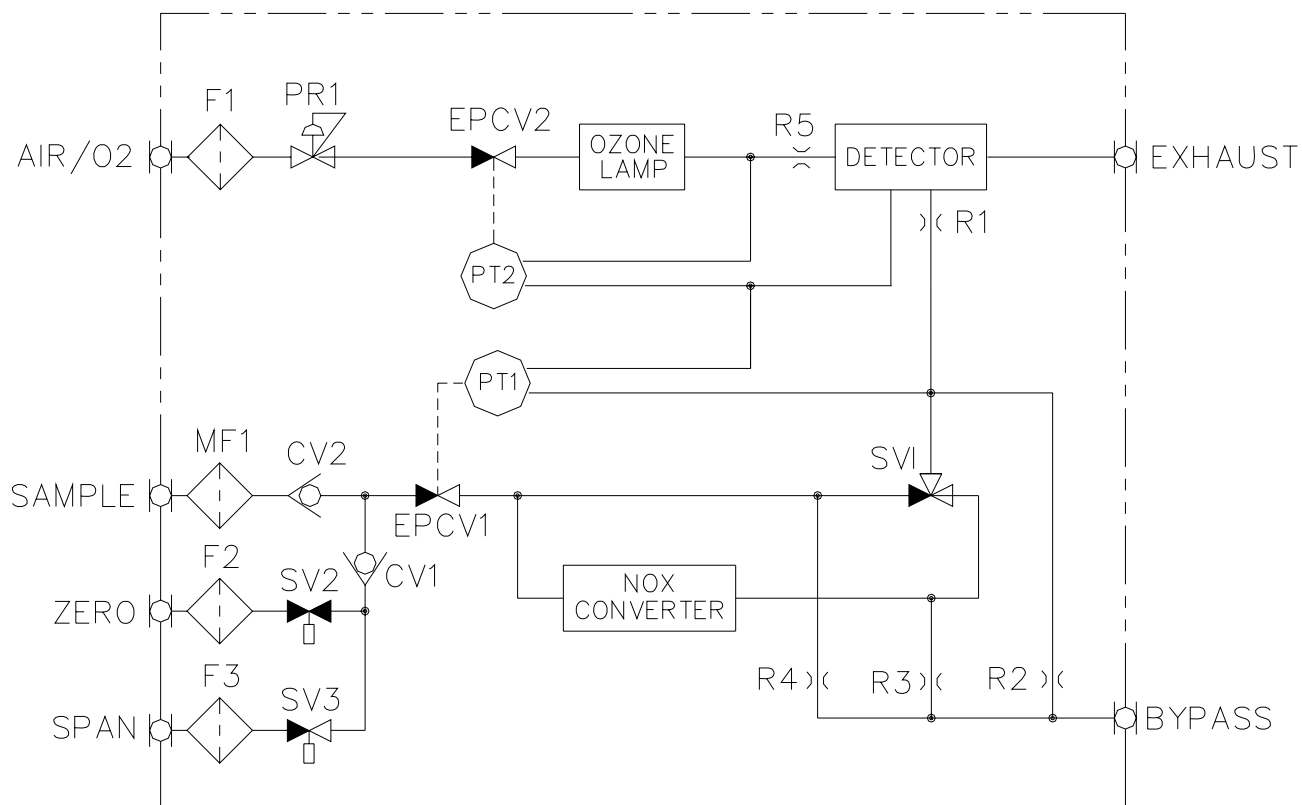

NOTE: Standard Analyzer with Optional Zero/Span Solenoids

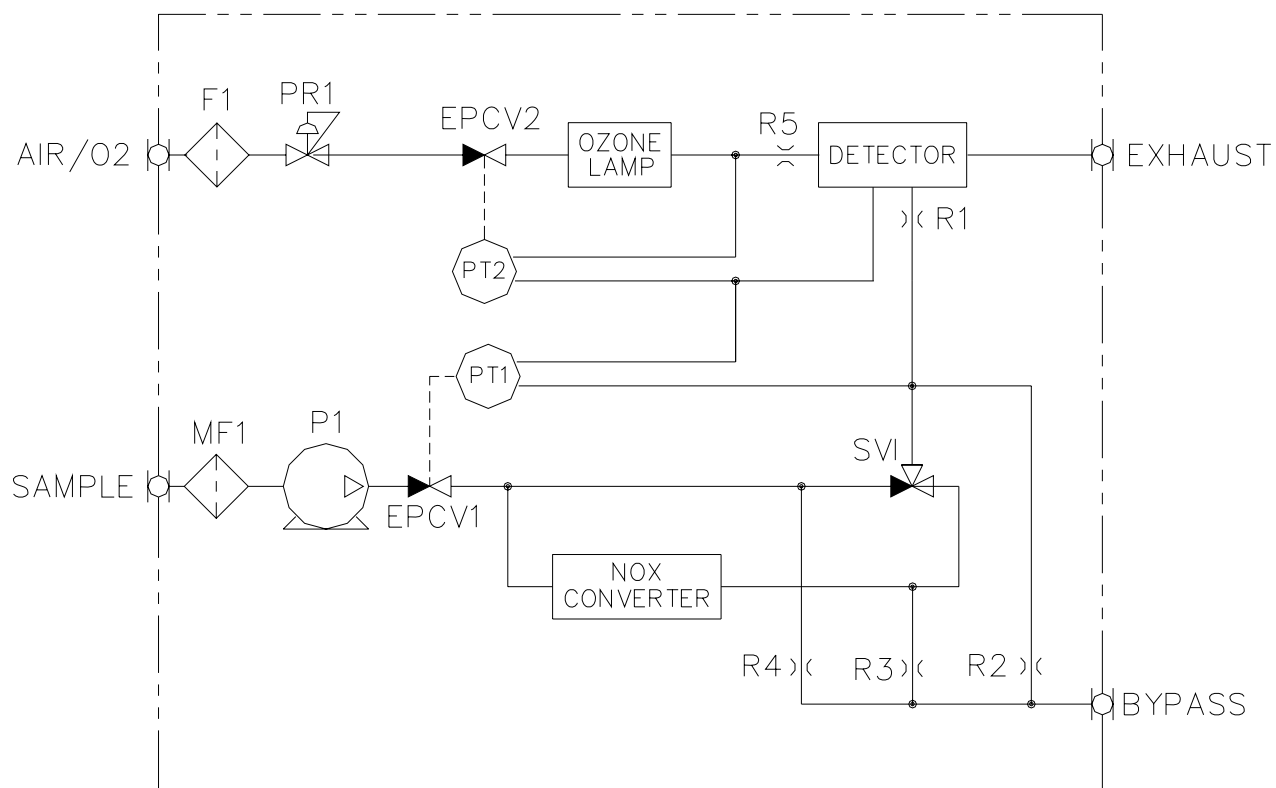

NOTE: Standard Analyzer with Optional Sample Pump

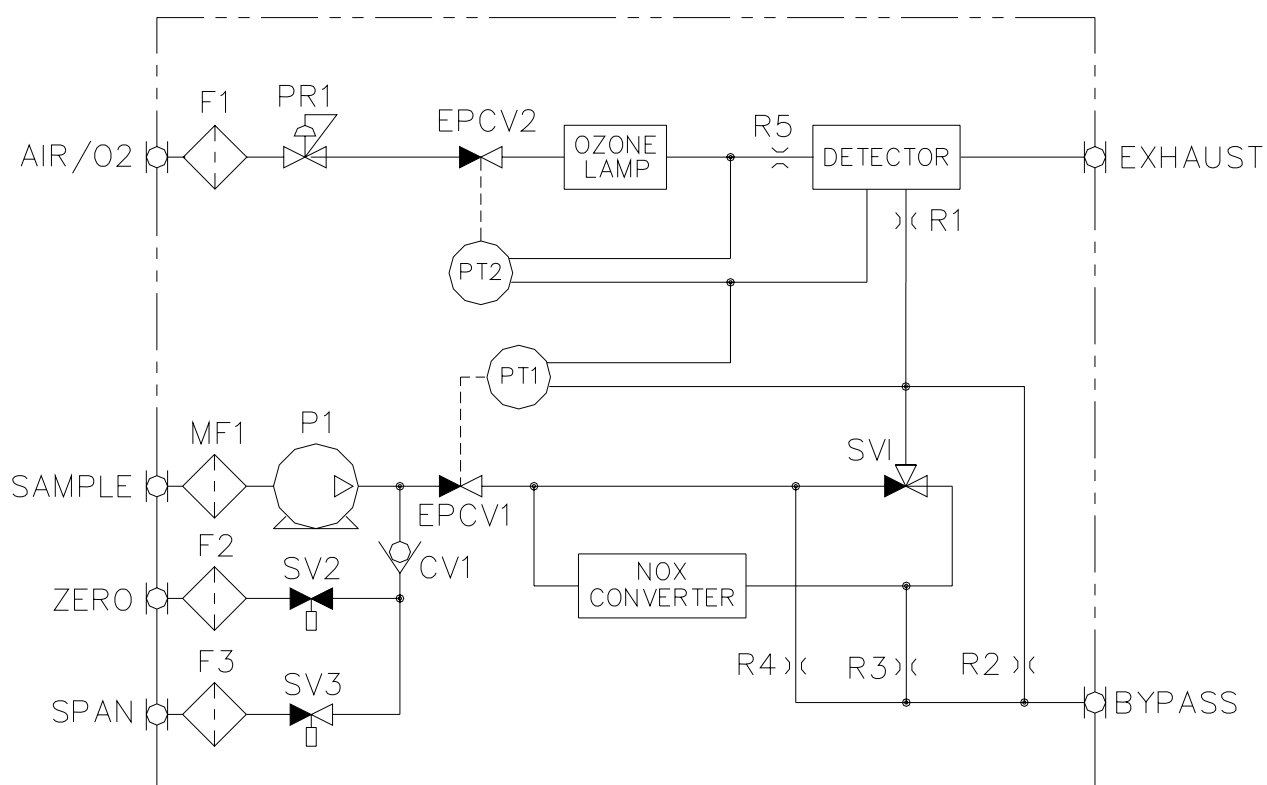

NOTE: Standard Analyzer with Optional Sample Pump and Zero/Span Solenoids

[illegible]

13.5 Starting With SERIAL NUNBER UO6081

1.0 INTRODUCTION

The Model 600 CLD Series Instruments starting with Serial Number U06061 have several new Hardware and Software features.

The Hardware includes the use of a new memory system, isolation of the analog output signals and 15 relays that are used to buffer the many new digital output signals that are now available.

The available digital signals consist of a SERVICE group, to externally monitor a number of parameters that provide insight for preventative maintenance and diagnostics decisions.

A second STATUS group, is provided to define the operation of the instrument such as Spanning, Zeroing, Calibrating and the Current Range (1, 2, 3, 4, AUTO).

The Software includes modifications to existing functions, changes to the Measurement screen, additional Short-Cut Keys and several New Functions that are listed as follows:

- MEASUREMENT**

| | |
|---------------------|----------------|
| Over Range | 888888 |
| C1/C3 | Removed |
| Zero | F5 |
| Span | F6 |
| Range Limits | F8 |
| Span Values | F9 |
| Outputs | F10 |

Note: The operator can use these Short-Cut Keys or continue to use existing procedures.

- NEW FUNCTIONS**

| | |
|----------------------------------|--|
| Auto Startup | F5, F7, F7 |
| ALARMS | F5, F7, (Use F6 to toggle ON/OFF) |
| Offsets& Gains | F4, F3, F5 |
| D/A Calibration | F5, F7, F8 |
| Cal Analog Outputs | F5, F8, (Use F8 to toggle ON/OFF) |
| Ignite ON Power Up | F5, F8, (Use F9 to toggle ON/OFF) |
| Save Data Archiving Time | F5, F7, F1, F5 (Use ENTER to change recording time) |
| User Digital Outputs | F5, F9 |
| Modifications | |
| Saved/Not good | F4, F2, F1 or F2 (To flow Zero or Span Gas) |
| Re-Set Calibration Values | F4, F5 |

2.0 OPERATION OF MEASUREMENT KEYS

Note: The ← & → Keys continue to be used to view a complete list of menu items.


2.1 Over range 888888

In the Sampling mode only, any value that exceeds the “software range” by more than 10% will be displayed as 888888.

Note: If the ranges have not been modified then the original factory physical ranges and the software ranges are the same.

2.2 Zero: Select the required range then press **F5**.

Note: For instruments with an internal Zero Solenoid select Calibration by Valves. (**F5**, **F2**, **F4**)

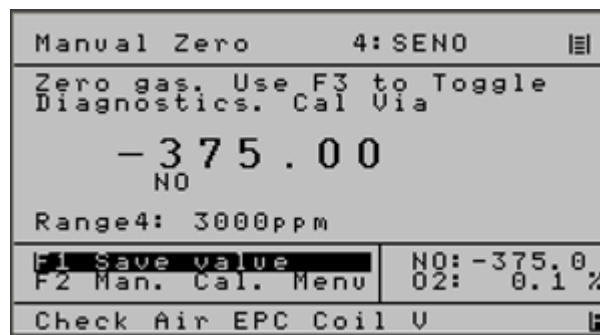

Zero Gas will be enabled and the observed results can be used to evaluate instrument performance.

Press **F3** to toggle to the Diagnostic screen for additional information

Press **F1** to save the value and complete a ZERO calibration for this channel and range

Press **MAIN**, **F1** to return to the MEASUREMENT screen to select other channels and ranges and repeat the process or press **F2** to return to the Manual Measuring screen

2.3 Span: Select the required range then press **F6**.

Note: For instruments with an Internal Span Solenoid select Calibration by Valves. (**F5**, **F2**, **F4**)

| | | | |
|--|--|-------------------|--|
| Manual Span | | 5: SWET | |
| Span gas. Use F3 to Toggle Diagnostics. Cal Via | | | |
| - 375.00 | | | |
| NO | | Conc.: 1860.00000 | |
| Range4: 3000ppm | | | |
| F1 Save value | | NO: -375.0 | |
| F2 Man. Cal. Menu | | O2: 0.1 % | |

Span Gas will be enabled and the observed results can be use to evaluate instrument performance.

Press **F3** to toggle to the Diagnostic screen for additional information

Press **F1** to save the new value and complete the SPAN calibration for this Range.

Note that the span gas value used for this range is highlighted and can be changed if necessary. Use the Enter key and the numeric keys as required

2.4 Range Limits: (**F8**);

| | | | |
|---------------------------|---------|------------|--|
| Setup | | 5: SWET | |
| Upper range limits [ppm] | | | |
| *** MUST be ASCENDING *** | | | |
| Range 1: | 3.00 | | |
| Range 2: | 30.00 | | |
| Range 3: | 300.00 | | |
| Range 4: | 3000.00 | | |
| *Do NOT Exceed Max Range* | | | |
| | | NO: -375.0 | |
| | | O2: 0.1 % | |
| Wed Jan 03 23:24:23 2001 | | | |

The analyzer is factory configured with 4 Physical Ranges of 3, 30, 300, 3,000 PPM.

The operator can change the number of ranges and select a more convenient full scale concentration if required.

Note: Do not exceed the maximum range of 3,000 PPM and always use ascending order as shown.

- Example a) For a single range instrument, set Range 1 to the desired value and all others to zero
- b) For a two range instruments, set Range 1 to the lowest value, Range 2 to the highest value and the others to zero.

2.5 Span Values: F9

| | | | |
|-----------------------------|----------|------------|--|
| Setup | | 1: SARA | |
| Span gas conc. range limits | | | |
| Range1: | 2.250 | 3.00 | |
| Range2: | 18.500 | 30.00 | |
| Range3: | 186.500 | 300.00 | |
| Range4: | 1860.000 | 3000.00 | |
| 02: | 0.000 | | |
| | | NO: -375.0 | |
| | | O2: 0.1 % | |
| Check Sample EPC Coil V | | | |

Use to define the concentration of the span gas that will be used to calibrate each range.

Note: The span gas value used for this range is highlighted and can be changed if necessary. Use the Enter key and the numeric keys as required.

2.6 Outputs: F10

| | | | |
|-------------------------|----------|-----------|--|
| assignment | | 4: SNOX | |
| Output Signal | | | |
| 1 | NO | | |
| 2 | NOx | | |
| 3 | NO2 | | |
| 4 | RealTime | | |
| | | NOx: 11.3 | |
| CAI Tel. (800) 959 0959 | | | |

Use the ↑ to select the desired Output

Use this screen to define the signals and their location that will be monitored by a remote reordering device.

3.0 NEW FUNCTIONS

3.1 Auto Start Up: (F5, F7, F7)

| | | |
|--------------------|-----------|---|
| Auto Startup | 5: SWET | ⏏ |
| •Auto Startup | : Off | |
| •Wait for [min] | : 1 | |
| •Calibrations | : 2 | |
| •Startrange | : 4 | |
| •Access Level | : 0 | |
| •Remote/Manual | : Manual | |
| •NO/NOx-Mode | : NO | |
| MAIN SAVE | NOx: 22.5 | |
| BACK SAVE | | |
| Check Air Pressure | | ⏏ |

Wait: The time delay in minutes before proceeding. If **Zero** is used the instrument will wait until all warnings are cleared to continue.

Calibrations: The number of attempts to complete a successful calibration as required in the operator defined Deviation Tables. If calibration is not successful the instrument will continue reporting results using the last completed calibration.

The analyzer can be configured use the previous calibration by selecting zero Calibrations.

Starting Range: When all defined actions are completed the analyzer will return to the Measurement Screen and at the range specified.

Access Level: The final access level

Remote/Manual: The final operating MODE


NO/NOx Mode: The final measurement MODE.

3.2 Alarms On/Off: (F5, F7)

All key analyzer parameters are stored in a secure memory location and retained when power is removed. In the event of an unexpected power failure it may be desirable to change some parameters until an operator can resume control.

This screen may be used to establish several desirable special instrument start-up parameters that define how the analyzer recovers from loss of AC power.

When enabled this screen will define the following:

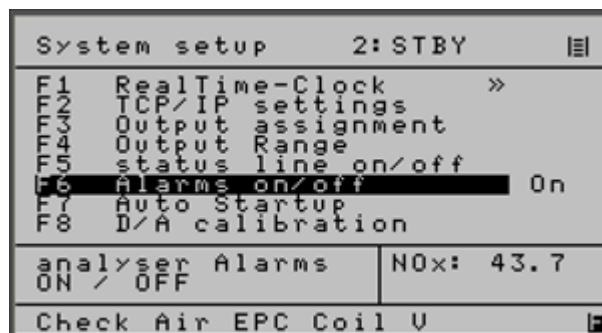

Use **F6** to toggle ON/OFF

The instrument has an extensive library of warning messages that will aid in the identification of various anomalous events and are displayed at the bottom of the screen. These messages will assist in Diagnostics and indicating the need for preventative maintenance

This screen provides an **option** to disable these messages during initial start-up or as may be desired for a particular application.

3.3 Offset & Gain: (F4, F3, F5)

| Offset and Gain | | |
|------------------------------|--------|----------|
| Zero/Span applied correction | | |
| | Offset | Gain |
| Range1: | 0.00 | 1.00 |
| Range2: | 0.00 | 1.00 |
| Range3: | 0.00 | 1.00 |
| Range4: | 0.00 | 1.00 |
| | | N0: -0.8 |
| Thu May 03 12:40:03 2007 | | |

This screen can be used to provide an additional means to display calibration deviations.

The OFFSET is the value stored during zero calibration.

The GAIN is the value stored during span gas calibration using the operator defined calibration gas.

An increasing or decreasing change to the OFFSET or GAIN when used in conjunction with "Deviations" will provide insight to both short and long term changes to system performance

3.4 Calibrate the Digital Outputs: F5, F7, F8 D/A Calibration F5, F7, F3 Output Assignment

| output scaling | | |
|-----------------------------|--------|---------|
| Calibrate D/A Outputs | | |
| Output | Offset | Gain |
| 1 | 0.0000 | 1.2210 |
| 2 | 1.8000 | 0.7000 |
| 3 | 1.8000 | 0.7000 |
| 4 | 1.8000 | 0.7000 |
| **Set Output-F3-Calibrate** | | |
| MAIN SAVE | F1: | 0% FS |
| BACK SAVE | F5: | 100% FS |
| CAI Tel. (800) 959 0959 | | |

Overview

The 600 CLD Series is designed to provide three analog outputs that can be configured as 0-10v, 4-20 ma, or 0-20 ma. With this version the outputs can also be configured to include an additional 1.0 volt and 5.0 volt output and a **calibration** capability.

The instrument can be configured to provide either voltage or current signals.

This screen is used to select the scaling of the current (ma) or voltage(1,5,10) range that is required by a remote recording device. The outputs can be calibrated to exactly match the results obtained on a PLC, Recorder, Data Logger or other remote recording device that may be connected to the analyzer.

The operator will first select the **OUTPUT ASSIGNMENT** screen and choose the output that is set to be calibrated. All outputs of interest may be selected. When calibration is completed, the operator will return the outputs to their original assignment.

The **D-A CALIBRATION** screen will be then be used to complete the calibration procedure. This screen provides a section that is used to record the zero signal corrections (zero offset) and a second area to record the 100% signal corrections (Gain) for each of the four output signals that may be defined to develop a voltage or current signal. Since this is a Digital to Analog conversion, the calibration will require the completion of a simple “trial and error” procedure. The operator will observe the results of a “zero or full scale (Gain) signal generated by the analyzer to the remote recording device and select a correction factor. The operator will save this value and then observe the results on the attached remote recording device.

The process of selection and saving for “zero” and “span” will be repeated until a satisfactory calibration is achieved. For 0-1V, 0-5V, 0-10V and a 0-20 ma outputs the Offset and Gain values are independent and do not interact. With the 4-20 ma output, the “Offset (zero)” and “Gain (span)” values interact and may require a few more trials.

The following is a table of typical values:

| OUTPUT | OFFSET | GAIN |
|---------------|---------------|-------------|
| 0-20 ma | 0.000 | 0.927 |
| 4-20 ma | 1.820 | 0.740 |
| 0-1 V | 1.300 | 0.820 |
| 0-5 V | 1.100 | 0.820 |
| 0-10 V | 1.050 | 0.820 |

- **Procedure**

3.4.1 From the Main Menu press **F5,F7,F3**, to obtain following screen:

| assignment | | 1: SARA | |
|------------|------------|------------|--|
| Output | Signal | | |
| 1 | NO | | |
| 2 | NOx | | |
| 3 | NO2 | | |
| 4 | SamplePres | | |
| | | NO: -375.0 | |
| | | O2: 0.1 % | |

3.4.2. Use the ↑ to highlight the outputs that require calibration.

Note: In the above example only Output 1 will be calibrated. Record the name of these signals, they will be restored.

3.4.3. Press enter to provide access to all the menu of signals that are available. (Real Time, THC, CH₄, Calibration, Sample Pressure, etc.)

3.4.4. Select Calibration and press **ENTER** to complete the selection.

Note: Any or all of the four outputs can be selected for calibration. This screen will not be used again until calibration has been completed.

3.4.5. Press **BACK** to return to the SYSTEM SETUP screen (**F5**, **F7** from main menu)

3.4.6. Press **F8** to obtain the following screen

| output scaling | | | |
|-----------------------------|--------|----------|----|
| Calibrate D/A Outputs | | | |
| Output | Offset | Gain | mA |
| 1 | 1.8000 | 0.7000 | |
| 2 | 1.8000 | 0.7000 | |
| 3 | 1.8000 | 0.7000 | |
| 4 | 1.8000 | 0.7000 | |
| **Set Output-F3-Calibrate** | | | |
| MAIN / BACK SAVE | | NOx: 0.1 | |
| FS F1: 0% F5: 100% | | | |
| CAI Tel. (800) 959 0959 | | | |

F5, F7, F8

3.4.7. Use the ↑ to select the desired output press **ENTER**.

3.4.8. Press **F1** to select a ZERO signal and observe the results on the remote device

- 3.4.9. Change the offset value press **BACK** to save the new value.
- 3.4.10. Press **F8** to return to the D-A Calibration screen and note the results on the remote device.
- 3.4.11. Repeat steps 8.0 thru 10.0 until a satisfactory ZERO calibration is achieved.
- 3.4.12. Complete steps 8.0 thru 10.0 for each of the remaining outputs that require calibration.
- 3.4.13. Press **F5** to produce a full scale (100%) signal.
- 3.4.14. Use the arrow keys to position the cursor at the require GAIN value.
- 3.4.15. Observe the results on the remote device and make a correction to the GAIN value for the output of interest. Press **BACK** to save this new value
- 3.4.16. Press **F8** to return to the D-A calibration screen
- 3.4.17. Observe the results on the remote device and repeat the steps to change the GAIN value by repeating steps 14.0 thru 16.0 as needed for each output.
- 3.4.18. Return to the OUTPUT Assignment screen **F5, F7, F3** from the main menu and change the output signals from CALIBRATE to their original values as defined in step 3.4.2.

3.5 Cal Analog Output: (F5, F8,)

| Measure setup 4: SEND | |
|------------------------------------|---------------------|
| F1 | Lowpass filter Tc |
| F2 | Purge time |
| F3 | T + P compensation |
| F4 | Set dual-mode times |
| F8 | Cal Analog Out On |
| F9 | Ignite on power up |
| set analog out during cal NO: 42.5 | |
| Tue May 15 08:47:49 2007 | |

Use F8 to toggle on/off

This will provide improved versatility and control of the NO and NOx output signals. When NO or NOx are assigned to specific outputs. The CAL ANALOG output can be enabled by the operator and the MODE selected at the Measuring Screen will be impressed at the selected output.

3.6 Save Data Archiving Time (F5, F7, F1, F5)

| Setup RTC 3: SMAN | |
|-----------------------------|--------------------------|
| F1 | Set time |
| F2 | Set autocalibration time |
| F3 | Select calibration range |
| F4 | Autocalibration on/off |
| F5 | Save Data Archiving Time |
| F10 | Show time |
| Archive Time(secs) 10 | |
| Archive Interval NO: -286.8 | |
| 0 = Off F5: SAVE 02: 0.1 % | |

Use ENTER to change recording time

3.7 User Digital Outputs

- Overview

The 600 HCD Series of instruments have 15 solid state, optically coupled, isolated relays that can be programmed by the operator to indicate the status of numerous digital conditions

The available digital signals consist of a SERVICE group, that can be used to externally monitor a number of conditions to aid in preventative maintenance and diagnostics.

A second STATUS group, is provided and is used to define the operation of the instrument such as Spanning, Zeroing, Calibrating and the current Range (1, 2, 3, 4, AUTO) etc.

The individual output signals can be operator selected and set to a **HOLD** or **CLEAR** mode.

In the **HOLD** mode an activated signal is retained until the operator returns to the **Digital Output Screen** and selects the appropriate output signal and performs a manual CLEAR.

In the **CLEAR** mode the signal will automatically change state when the microprocessor detects that the noted condition no longer exists.

The operator can select from the following the desired **SERVICE** or **STATUS** items that are to be digitally monitored.

- **SERVICE**

| Text | Signal | Min | Max | Definition | |
|---------|----------------------------|-----|-----|------------|-----------------|
| SampleP | | | 3.8 | 3.9 | Sample pressure |
| AirP | Sample Pressure ,Check | | 14 | 16 | Air pressure |
| OvenT | Air Pressure ,Check | | 84 | 86 | Oven temp |
| ConvT | Oven Temp ,Check | | 204 | 208 | Converter temp |
| PumpT | Converter Temp ,Check | | 84 | 102 | Pump temp |
| DiodeT | Pump Temp ,Check | | 1.1 | 2.1 | Diode temp |
| CellT | Diode Temp ,Check | | 65 | 67 | Cell temp |
| DryT | Cell Temp ,Check | | 2 | 7 | Peltier temp |
| O2T | Peltier Gas Temp ,Check | | 1 | 7 | Not used |
| SEV | O2 Temp ,Check | | 1 | 5 | Not used |
| AEV | EPC Coil Sample ,Check | | 2 | 8 | Sample EPC Coil |
| OR | EPC Coil Air ,Check | | 2 | 8 | Ozone EPC Coil |
| AO | Range overflow | | | | |
| AU | ADC Range Overflow | | | | |
| R1NC | ADC Range Underflow | | | | |
| R2NC | Range 1 is not calibrated | | | | |
| R3NC | Range 2 is not calibrated | | | | |
| R4NC | Range 3 is not calibrated | | | | |
| RCP | Range 4 is not calibrated | | | | |
| LoC | Low concentration Warning | | | | |
| HiC | High concentration Warning | | | | |
| NH3T | | | | | NOT Used |
| OFF | All Alarms OFF | | | | |


-
-

• STATUS

GenAlarm General Alarm
 InRem In Remote
 CalAlarm Cal Alarm

| | |
|---------|--------------|
| AutoR | AutoRange |
| R1 | Range 1 |
| R2 | Range 2 |
| R3 | Range 3 |
| R4 | Range 4 |
| InCal | In Calibrate |
| Zero | In Zero |
| Span | In Span |
| Sample | In Sample |
| InNO | In NO Mode |
| InWet | In Wet Mode |
| InOflow | In Overflow |
| InNH3 | In NH3 Mode |

TYPICAL RELAY

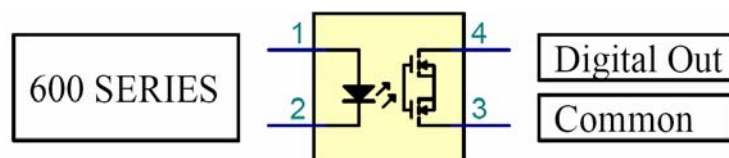

These contacts (3, 4) will drive continuously up to 500 MA using a customer voltage supply that does not to exceed 60 VDC.

- **OPERATION**

Use **F5**, **F9** to select the first eight outputs.

Use the \updownarrow to select the desired output.

Press ENTER and use \updownarrow to select desired item.

| User DO I | |
|-----------------------------|---------|
| 1 | SampleP |
| 2 | Pump1 |
| 3 | Off |
| 4 | Off |
| 5 | Off |
| 6 | Off |
| 7 | Off |
| 8 | Off |
| ***** | |
| *Set Unused * | |
| *Channels Off* | |
| ***** | |
| F1 9 to 15 DO's & NOxa 0.1d | |
| MAIN/BACK to SAVE | |
| Mon Jan 01 03:06:07 2001 | |

Press **F1** to observe the remaining seven outputs

Program as desired per the above

| User DO II | |
|--------------------------|---------|
| 9 | FilterT |
| 10 | Off |
| 11 | Off |
| 12 | Off |
| 13 | Off |
| 14 | Off |
| 15 | Off |
| ***** | |
| *Set Unused * | |
| *Channels Off* | |
| ***** | |
| F2 Hold / Clear | |
| F1 1 to 8 DO's NOx: 0.1 | |
| MAIN/BACK to SAVE | |
| Mon Jan 01 03:05:49 2001 | |

4.0 CHANGES TO EXISTING FUNCTIONS

4.1 Saved or NOT GOOD

During Manual Calibration the following screens will be displayed to indicate the instruments response to the value of the zero or span gas using the amount that the operator defined in the deviation table.

| Manual Zero | |
|--|-----------|
| Zero gas. Use F3 to Toggle Diagnostics. Cal Via Valves ***Saved*** | |
| - 0 . 0 1 | |
| NOx | |
| Range2: 30.00ppm | |
| F1 Save value | NOx: -0.0 |
| F2 Man. Cal. Menu | |
| CAI Tel. (800) 959 0959 | |

| Manual Span | |
|---|----------|
| Span gas. Use F3 to Toggle Diagnostics. Cal Via Valves *NOT GOOD* | |
| 0 . 0 0 | |
| NOx Conc.: 28.46000 | |
| Range2: 30.00ppm | |
| F1 Save value | NOx: 0.0 |
| F2 Man. Cal. Menu | |
| Mon Jan 01 00:05:06 2001 | |

The above is shown using Zero Gas

From Measurement use: **F5 or F6**

From Main Menu use: **F4, F2, F1 or F2**

| Setup | | |
|--------------------------|----------|----------|
| Deviations[%] | | |
| | absolute | relative |
| Range1: | 10.00 | 10.00 |
| Range2: | 10.00 | 10.00 |
| Range3: | 10.00 | 10.00 |
| Range4: | 10.00 | 10.00 |
| | | NOx: 0.0 |
| Mon Jan 01 00:04:32 2001 | | |

F5, F2, F3

Note :This screen is used by the operator to define the maximum acceptable limits of the Zero and Span gas for both Manual and Automatic Calibrating.

4.2 Reset Calibration Values

When the re-set calibrations value function is used all recorded deviations will be set to zero

| Def cal value | |
|---|-----|
| Are you sure to reset calibration values to default values? | |
| F1 | yes |
| F2 | no |
| NOx: -0.0 | |
| Mon Jan 01 00:01:36 2001 | |

F4,F5

| Span gas devs | | | |
|--------------------------|-----------------|------------|------|
| Span gas devs. verifying | | | |
| | measured values | deviations | |
| | | Uar | %FS |
| Range1: | 0.0 | 0.00 | 0.00 |
| Range2: | 0.0 | 0.00 | 0.00 |
| Range3: | 0.0 | 0.00 | 0.00 |
| Range4: | 0.0 | 0.00 | 0.00 |
| | | NOx: -0.0 | |
| Mon Jan 01 00:02:08 2001 | | | |

F4 F3, F4

(Used to observe Auto Cal Results)

| Span gas devs | | |
|--------------------------|------|-----------|
| Span gas deviations [%]: | | |
| | abs | rel |
| Range1: | 0.00 | 0.00 |
| Range2: | 0.00 | 0.00 |
| Range3: | 0.00 | 0.00 |
| Range4: | 0.00 | 0.00 |
| | | NOx: -0.0 |
| CAI Tel. (800) 959 0959 | | |

F4, F3, F2

(Used to observe Manual Cal results)

The above are the new deviations after the operator elects to re-set the calibration values

TABLE OF CONTENTS

ADDENDUM

13.5 STARTING WITH SERIAL #U06081

| <u>Section</u> | <u>Title</u> | <u>Page</u> |
|-----------------------|--|--------------------|
| 1.0. | Introduction ----- | |
| 2.0. | Operation of Measurement Keys ----- | |
| | 2.1. Over Range----- | |
| | 2.2. Diagnostics----- | |
| | 2.3. Zero----- | |
| | 2.4. Span ----- | |
| | 2.5. Range Limits----- | |
| | 2.6. Span Values ----- | |
| | 2.7. Outputs----- | |
| 3.0 | New Functions ----- | |
| | 3.1. Auto Start Up ----- | |
| | 3.2. Alarms On/Off ----- | |
| | 3.3. Offset & Gain ----- | |
| | 3.4. Calibrate The Analog Outputs ----- | |
| | 3.5. Save Data Archiving Time ----- | |
| | TABLE A ----- | |
| | 3.6. User Digital Outputs ----- | |
| | TABLE B----- | |
| | TABLE C----- | |
| | 3.7. Cal Analog Output----- | |
| 4.0. | Changes To Existing Functions ----- | |
| | 4.1. Saved or Outside Limits ----- | |
| | 4.2. Calibration Deviations----- | |
| | 4.3. Flow Zero or Span Gas----- | |
| | 4.4. Reset Calibration Values ----- | |
| | TABLE D----- | |

13.5 Starting With SERIAL NUMBER U06081

1.0 INTRODUCTION

The Model 600 CLD Series Instruments starting with Serial Number U06081 have several new Hardware and Software features.

The Hardware includes the use of a new memory system, isolation of the analog output signals and 15 relays that are used to buffer the many new digital output signals that are now available. **SEE TABLE D**

The available digital signals consist of a SERVICE group, to externally monitor a number of parameters that provide insight for preventative maintenance and diagnostics decisions.

A second STATUS group, is provided to define the operation of the instrument such as Spanning, Zeroing, Calibrating and the Current Range (1, 2, 3, 4, AUTO).

The Software includes modifications to existing functions, changes to the Measurement screen, additional Short-Cut Keys and several New Functions that are listed as follows:

- **MEASUREMENT**

| | |
|---------------------|---------------|
| Over Range | 888888 |
| Zero | F5 |
| Span | F6 |
| Range Limits | F8 |
| Span Values | F9 |
| Outputs | F10 |

Note: The operator can use these Short-Cut Keys or continue to use existing procedures.

- **NEW FUNCTIONS**

| | |
|---------------------------------|---|
| Auto Startup | F5, F7, F7 |
| ALARMS | F5, F7, (Use F6 to toggle ON/OFF) |
| Offsets& Gains | F4, F3, F5 |
| D/A Calibration | F5, F7, F8 |
| Cal Analog Outputs | F5, F8, (Use F8 to toggle ON/OFF) |
| Save Data Archiving Time | F5, F7, F1, F5 (Use ENTER to change record time) |
| User Digital Outputs | F5, F9 |

- **MODIFICATIONS**

| | |
|-----------------------|--|
| Saved/Not good | F4, F2, F1 or F2 (To flow Zero or Span Gas) |
|-----------------------|--|

2.0 OPERATION OF MEASUREMENT KEYS

Note: The ← &→ Keys continue to be used to view a complete list of menu items.

2.1 Over Range 888888

In the MEASUREMENT mode only, any value that exceeds the “range” by more than 10% will be displayed as 888888.

2.2 Diagnostics: Use **F3** to toggle between MEASUREMENT and DIAGNOSTIC.**2.3 Zero:** From the MEASUREMENT Screen select the required range for calibration then press **F5**.

Note: For instruments with an internal Zero Solenoid select Calibration by Valves. (**Main, F5, F2, F4**)

| | |
|--|---|
| Manual Zero 4:SEN0 | Manual Zero 3:SMAN |
| Zero gas. Use F3 to Toggle Diagnostics. Cal Via | Zero gas. Use F3 To Toggle Diagnostics Cal Via *Using Current Temp & Press* |
| -375.00 NO | 0.08 NOx |
| Range4: 3000ppm | Range1: 100.0ppm |
| F1 Save value NO:-375.0 F2 Man. Cal. Menu O2: 0.1 % | F1 Save value NOx: 0.1 F2 Man. Cal. Menu O2: 0.0 % |
| Check Air EPC Coil V | Thu May 15 06:15:11 2008 |

2 versions

Zero Gas will be enabled and the observed results can be used to evaluate instrument performance.

Press **F3** to toggle to the Diagnostic screen for additional information

Press **F1** to save the value and complete a ZERO calibration for this channel and range

Press **MAIN, F1** to return to the MEASUREMENT screen to select other channels and ranges and repeat the process. Or press **F2** to return to the Manual Calibration Screen

2.4 Span: From the MEASUREMENT Screen select the required range then press **F6**.

Note: For instruments with an Internal Span Solenoid select Calibration by Valves. (**Main, F5, F2, F4**)

| | |
|--|-------------------|
| Manual Span 5: SWET | |
| Span gas. Use F3 to Toggle Diagnostics. Cal Via | |
| - 375.00 | |
| NO | Conc.: 1860.00000 |
| Range4: 3000ppm | |
| F1 Save value | NO: -375.0 |
| F2 Man. Cal. Menu | O2: 0.1 % |

| | |
|---|-----------------|
| Manual Span 5: SWET | |
| Span gas. Use F3 To Toggle Diagnostics Cal Via *Using Current Temp & Press* | |
| 0.00 | |
| NOx | Conc.: 90.00000 |
| Range1: 100.0ppm | |
| F1 Save value | NOx: 0.0 |
| F2 Man. Cal. Menu | O2: 0.0 % |
| CAI Tel. (800) 959 0959 | |

2 versions

Span Gas will be enabled and the observed results can be used to evaluate instrument performance.

Press **F3** to toggle to the Diagnostic screen for additional information

Press **F1** to save the new value and complete the SPAN calibration for this Range.

Note: The span gas value used for this range is highlighted and can be changed if necessary. Use the Enter key and the numeric keys as required

Press **MAIN, F1** to return to the MEASUREMENT screen to select other ranges and repeat the process or press **F2** to return to the Manual Calibration screen

2.5 Range Limits: F8 From the MEASUREMENT Screen;

| Setup 5: SWET | |
|---------------------------|---------|
| Upper range limits [ppm] | |
| *** MUST be ASCENDING *** | |
| Range 1: | 3.00 |
| Range 2: | 30.00 |
| Range 3: | 300.00 |
| Range 4: | 3000.00 |
| *Do NOT Exceed Max Range* | |
| NO: -375.0 | |
| O2: 0.1 % | |
| Wed Jan 03 23:24:23 2001 | |

| Setup 1: SARA | |
|-------------------------------|-----------|
| Upper range limits [ppm] | |
| *** MUST be ASCENDING *** | |
| Range1: | 100.00 |
| Range2: | 500.00 |
| Range3: | 1000.00 |
| Range4: | 3000.00 |
| Maximum Range Limit: 3000.000 | |
| F1 Save with new | NOx: 0.1 |
| autorange Up/Down | O2: 0.0 % |
| CAI Tel. (800) 959 0959 | |

2 versions

The standard analyzer is factory configured with 4 Physical Ranges of 3, 30, 300, 3,000 PPM.

The optional high level analyzer is factory configured with 4 Physical Ranges of 5, 50, 500, 500

The operator can change the number of ranges and select a more convenient full scale concentration if required.

Note: Do not exceed the maximum range set by the factory and always use ascending order as shown.

- Example a) For a single range instrument, set Range 1 to the desired value and all others to zero
- b) For a two range instruments, set Range 1 to the lowest value, Range 2 to the highest value and the others to zero. .

2.6 Span Values: F9 From the MEASUREMENT Screen

| Setup 1: SARA | |
|-----------------------------|------------------|
| Span gas conc. range limits | |
| Range1: | 2.250 3.00 |
| Range2: | 18.500 30.00 |
| Range3: | 186.500 300.00 |
| Range4: | 1860.000 3000.00 |
| O2: | 0.000 |
| NO: -375.0 | |
| O2: 0.1 % | |
| Check Sample EPC Coil U | |

Use to define the concentration of the span gas that will be used to calibrate each range.

Note: The span gas value used for this range is highlighted and can be changed if necessary. Use the Enter key and the numeric keys as required.

2.7 Outputs: F10 From the MEASUREMENT Screen

| | | | |
|-------------------------|----------|-----------|--|
| assignment | | 4: SNOX | |
| Output | Signal | | |
| 1 | NO | | |
| 2 | NOx | | |
| 3 | NO2 | | |
| 4 | RealTime | | |
| | | NOx: 11.3 | |
| CAI Tel. (800) 959 0959 | | | |

Use the ↑ to select the desired Output. Press Enter to select
Use the ↓ to select the desired Signal. Press Enter to select

Use this screen to define the signals and their location that will be monitored by a remote reordering device.

Note: NO₂ only active in the “Switching Mode”

SEE TABLE D

3.0 NEW FUNCTIONS

3.1 Auto Start Up: (Main, F5, F7, F7)

| | | | |
|--------------------|---|-----------|--|
| Auto Startup | | 5: SWET | |
| •Auto Startup | : | Off | |
| •Wait for [min] | : | 1 | |
| •Calibrations | : | 2 | |
| •Startrange | : | 4 | |
| •Access Level | : | 2 | |
| •Remote/Manual | : | Manual | |
| •NO/NOx-Mode | : | NO | |
| MAIN SAVE | | NOx: 22.5 | |
| BACK SAVE | | | |
| Check Air Pressure | | | |

All key analyzer parameters are stored in a secure memory location and retained when power is removed

In the event of an unexpected power failure it may be desirable to change some parameters until an operator can resume control.

This screen may be used to establish several desirable special instrument start-up parameters that define how the analyzer recovers from loss of AC power

When enabled this screen will define the following:

Wait: The time delay in minutes before proceeding. If “Zero” is selected the instrument will wait until all temperature warnings are cleared.

Calibrations: The number of attempts to complete a successful calibration as required in the operator defined Deviation Tables. If calibration is not successful the instrument will continue reporting results using the last completed calibration.

The analyzer will utilize the last completed calibration by selecting zero for Calibrations.

Starting Range: When all defined actions are completed the analyzer will return to the Measurement Screen and at the range specified.

Access Level: The user level at Start Up.

Remote/Manual: Put the analyzer in either “Remote” or “Manual” at Start Up


NO/NOx Mode: Put the analyzer in NO or NOx or NO2 mode at Start Up.

3.2 Alarms On/Off: (Main, F5, F7)

All key analyzer parameters are stored in a secure memory location and retained when power is removed. In the event of an unexpected power failure it may be desirable to change some parameters until an operator can resume control.

This screen may be used to establish several desirable special instrument start-up parameters that define how the analyzer recovers from loss of AC power.

When enabled this screen will define the following:

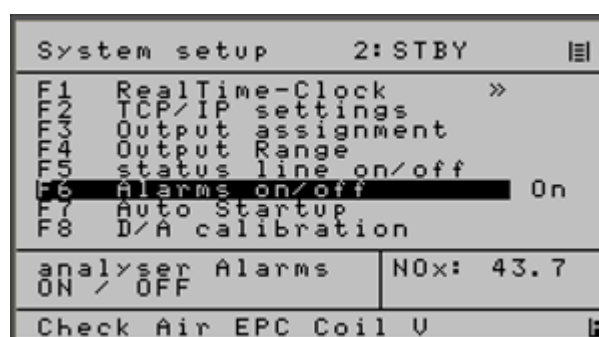

Use **F6** to toggle ON/OFF

The instrument has an extensive library of warning messages that will aid in the identification of various anomalous events and are displayed at the bottom of the screen. These messages will assist in diagnostics and indicating the need for preventative maintenance

This screen provides an **option** to disable these messages during initial start-up or as may be desired for a particular application.

3.3 Offset & Gain: (Main, F4, F3, F5)

| Offset and Gain | | | |
|------------------------------|--------|------|------|
| Zero/Span applied correction | | | |
| | Offset | Gain | |
| Range1: | 0.00 | 1.00 | |
| Range2: | 0.00 | 1.00 | |
| Range3: | 0.00 | 1.00 | |
| Range4: | 0.00 | 1.00 | |
| | | N0: | -0.8 |
| Thu May 03 12:40:03 2007 | | | |

This screen can be used to provide an additional means to display calibration deviations.

The OFFSET is the value stored during zero calibration.

The GAIN is the value stored during span gas calibration using the operator defined calibration gas.

An increasing or decreasing change to the OFFSET or GAIN when used in conjunction with "Deviations" will provide insight to both short and long term changes to system performance

Note: Reset calibration values will reset OFFSETS & GAINS to zero and 1 respectively

3.4 Calibrate the Analog Outputs: (Main, F5, F7, F8) D/A Calibration (Main, F5, F7, F3) Output Assignment

- **Overview**

The 600 CLD Series is designed to provide four analog outputs that can be configured as 0-1 VDC, 0-5 VDC, 0-10 VDC, 4-20 mA, or 0-20 mA.

This screen (**Main, F5, F7, F8**) is used to select the scaling of the current (mA) or voltage (1, 5, 10) range that is required by a remote recording device. The outputs can be calibrated to exactly match the results obtained on a PLC, Recorder, Data Logger or other remote recording device that may be connected to the analyzer.

The operator will first select the **OUTPUT ASSIGNMENT (Main, F5, F7, F3)** screen and choose the output that is set to be calibrated. By selecting "calibrate" as the output all outputs of interest may be selected. When calibration is completed, the operator will return the outputs to their original assignment.

The **D-A CALIBRATION** screen will then be used to complete the calibration procedure.

| | | | |
|-----------------------------|----------|--------|-----------|
| output scaling 5:SDRY | | | ≡ |
| Calibrate D/A Outputs | | | |
| Output | Offset | Gain | mA |
| 1 | 1.8000 | 0.7000 | |
| 2 | 1.8000 | 0.7000 | |
| 3 | 1.8000 | 0.7000 | |
| 4 | 1.8000 | 0.7000 | |
| **Set Output-F3-Calibrate** | | | |
| MAIN / | BACK | SAVE | NO: 0.1 |
| FS F1: 0% | F5: 100% | | O2: 0.0 % |
| Mon Jan 01 01:13:03 2001 | | | ☐ |

This screen provides a section that is used to record the zero signal corrections (zero offset) and a second area to record the 100% signal corrections (Gain) for each of the four output signals that may be defined to develop a voltage or current signal.

Since this is a Digital to Analog conversion, the calibration will require the completion of a simple “trial and error” procedure. The operator will observe the results of a “zero” or “span”(Gain) signal generated by the analyzer to the remote recording device and select a correction factor. The operator will save this value and then observe the results on the attached remote recording device.

The process of selection and saving for “zero” and “span” will be repeated until a satisfactory calibration is achieved. For 0-1VDC, 0-5VDC, 0-10VDC and a 0-20 mA outputs the Offset and Gain values are independent and do not interact. With the 4-20 mA output, the “Offset (zero)” and “Gain (span)” values interact and may require a few more trials.

The following is a table of typical values that can be used for start points for offsets and gains for different output types

| OUTPUT | OFFSET | GAIN |
|---------|--------|-------|
| 0-20 ma | 0.000 | 0.927 |
| 4-20 ma | 1.820 | 0.740 |
| 0-1 V | 1.300 | 0.820 |
| 0-5 V | 1.100 | 0.820 |
| 0-10 V | 1.050 | 0.820 |

- Procedure

3.4.1 From the Main Menu press **F5,F7,F3**, to obtain following screen:

| assignment | | 1: SARA | |
|------------|------------|------------|--|
| Output | Signal | | |
| 1 | NO | | |
| 2 | NOx | | |
| 3 | NO2 | | |
| 4 | SamplePres | | |
| | | NO: -375.0 | |
| | | O2: 0.1 % | |

3.4.2. Use the ↑ to highlight the outputs that require calibration.

3.4.3. Press enter so you provide access to all the menu of signals that are available.
(Real Time, NO, NOx, Calibration, Sample Pressure, etc.)

3.4.4. Select Calibration and press **ENTER** to complete the selection.

Note: In the example below only Output 1 will be calibrated

Record the name of these signals, they will be restored.

| assignment | | 1: SARA | |
|-------------------------|------------|-----------|--|
| Output | Signal | mA | |
| 1 | Calibrate | | |
| 2 | NOx | | |
| 3 | NO2 | | |
| 4 | SamplePres | | |
| | | NO: 0.1 | |
| | | O2: 0.0 % | |
| CAI Tel. (800) 959 0959 | | | |

Note: Any or all of the four outputs can be selected for calibration. This screen will not be used again until calibration has been completed.

3.4.5. Press **BACK** to return to the SYSTEM SETUP screen (**Main, F5, F7**)

3.4.6. Press **F8** to obtain the following screen

| | | | |
|-----------------------------|--------|----------|------|
| output scaling | | | ⏏ |
| Calibrate | D/A | Outputs | |
| Output | Offset | Gain | mA ↕ |
| 1 | 1.8000 | 0.7000 | |
| 2 | 1.8000 | 0.7000 | |
| 3 | 1.8000 | 0.7000 | |
| 4 | 1.8000 | 0.7000 | |
| **Set Output-F3-Calibrate** | | | |
| MAIN / BACK SAVE | | NOx: 0.1 | |
| FS F1: 0% F5: 100% | | | |
| CAI Tel. (800) 959 0959 | | | |

Main, F5, F7, F8

3.4.7. Use the \uparrow to select the desired output press **ENTER**.

3.4.8. Press **F1** to select a ZERO signal and observe the results on the remote device

3.4.9. Change the offset value press **BACK** to save the new value.

3.4.10. Press **F8** to return to the D-A Calibration screen and note the results on the remote device.

3.4.11. Repeat steps 3.4.8 thru 3.4.10 until a satisfactory ZERO calibration is achieved.

3.4.12. Complete steps 3.4.8 thru 3.4.10 for each of the remaining outputs that require calibration.

3.4.13. Press **F5** to produce a full scale (100%) signal.

3.4.14. Use the arrow keys to position the cursor at the require GAIN value.

3.4.15. Observe the results on the remote device and make a correction to the GAIN value for the output of interest. Press **BACK** to save this new value

3.4.16. Press **F8** to return to the D-A calibration screen

3.4.17. Observe the results on the remote device and repeat the steps to change the GAIN value by repeating steps 3.4.14 thru 3.4.16 as needed for each output.

3.4.18. Return to the OUTPUT Assignment screen **F5, F7, F3** from the main menu and change the output signals from CALIBRATE to their original values as defined in step 3.4.1.

3.5 Save Data Archiving Time (Main, F5, F7, F1, F5)

Archive Time is the Time in seconds between each set of data points. If “zero” no data is stored in the SEC data files. The SEC data files are in .CSV format for direct import into Excel. CAI can provide the tools necessary to download these files.

| Setup RTC | | 3: SMAN | |
|------------------|--------------------------|------------|----|
| F1 | Set time | | |
| F2 | Set autocalibration time | | |
| F3 | Select calibration range | | |
| F4 | Autocalibration on/off | | |
| F5 | Save Data Archiving Time | | |
| | Archive Time(secs) | | 10 |
| F10 | Show time | | |
| Archive Interval | | N0: -286.8 | |
| 0 = Off F5:SAVE | | 02: 0.1 % | |

Use ENTER to change recording time

SEE TABLE A

TABLE A**600 SERIES CLD DATA ARCHIVE FILES**

Time,
Date,
Month,
Year,
Error Index,
TimeStamp,
NO Conc,
NO2 Conc,
NOx Conc,
NH3 Conc,
Concentration,
Detector Volts,
Range,
Auto / Manual,
Span Gas,
Offset,
Gain,
Sample Pressure,
Sample Flow,
Sample EPV Volts,
Air Pressure,
Ozone Flow,
Ozone EPC Volts,
Diode Temp,
Cell Temp,
Oven Temp,
Pump Temp,
Converter Temp,
Dryer Temp,
O2 Detector Temp,
NH3 Conv Temp,
O2 Concentration,
O2 Detector Volts,
O2 Offset,
O2 Gain,
Wet / Dry,
Meas Mode,
Local / Remote,
Converter

3.6 User Digital Outputs

- **Overview**

The 600 CLD Series of instruments have 15 solid state, optically coupled, isolated relays that can be programmed by the operator to indicate the status of numerous digital conditions

The available digital signals consist of a SERVICE group, that can be used to externally monitor a number of conditions to aid in preventative maintenance and diagnostics. **SEE TABLE B & D**

A second STATUS group, is provided and is used to define the operation of the instrument such as Spanning, Zeroing, Calibrating and the current Range (1, 2, 3, 4, AUTO) etc. **SEE TABLE C & D**

The individual output signals can be operator selected and set to a **HOLD** or **CLEAR** mode.

In the **HOLD** mode an activated signal is retained until the operator returns to the **Digital Output Screen** and selects the appropriate output signal and performs a manual CLEAR. After performing a Clear Operation, the operator must press F2 again to put the outputs back into the Hold mode.

| User | DO I | 4: SENO | |
|--------------------------|--------------|----------------|--|
| 1 | Off | | |
| 2 | Off | | |
| 3 | Off | | |
| 4 | Off | ***** | |
| 5 | Off | *Set Unused * | |
| 6 | Off | *Channels Off* | |
| 7 | Off | ***** | |
| F2 | Hold / Clear | | |
| F1 | 8 to 15 DO's | NO: 0.1 | |
| MAIN/BACK | to SAVE | O2: 0.0 % | |
| Mon Jan 01 01:40:39 2001 | | | |

In the **CLEAR** mode the signal will automatically change state when the microprocessor detects that the noted condition no longer exists.

| User | DO II | 5: SDRY | |
|-------------------------|-------------|----------------|--|
| 8 | Off | | |
| 9 | Off | | |
| 10 | Off | | |
| 11 | Off | ***** | |
| 12 | Off | *Set Unused * | |
| 13 | Off | *Channels Off* | |
| 14 | Off | ***** | |
| 15 | Off | | |
| F1 | 1 to 7 DO's | NO: 0.1 | |
| MAIN/BACK | to SAVE | O2: 0.0 % | |
| CAI Tel. (800) 959 0959 | | | |

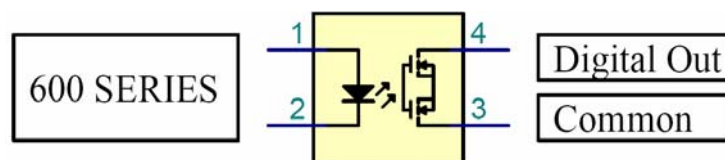
The operator can select from the following the desired **SERVICE** or **STATUS** groups that are to be digitally monitored.

TABLE B

| Index | Service Group | User DO Screen Display |
|--------------|--|-------------------------------|
| 0 | | |
| 1 | Sample Pressure Failure | SampleP |
| 2 | Air Pressure Failure | AirP |
| 3 | Oven Temp Failure (hclld only) | OvenT |
| 4 | Converter Temp Failure | ConvT |
| 5 | Pump Temp Failure | PumpT |
| 6 | Diode Temp Failure | DiodeT |
| 7 | Cell Temp Failure | CellT |
| 8 | Peltier Gas Temp Failure (hclld only) | DryT |
| 9 | O2 Temp Failure (650 only) | O2T |
| 10 | EPC Coil Sample Failure | SEV |
| 11 | EPC Coil Air Failure | AEV |
| 12 | Range overflow | OR |
| 13 | ADC Range Overflow | AO |
| 14 | ADC Range Underflow | AU |
| 15 | Range 1 is not calibrated | R1NC |
| 16 | Range 2 is not calibrated | R2NC |
| 17 | Range 3 is not calibrated | R3NC |
| 18 | Range 4 is not calibrated | R4NC |
| 19 | Reaction chamber pressure | RCP |
| 20 | Low concentration Warning | LoC |
| 21 | High concentration Warning | HiC |
| 22 | NH3 Converter Temp Failure (605 only) | NH3T |
| 23 | dummy text for RTC | Off |
| 24 | General Alarm | GenAlarm |
| 26 | Cal Alarm | CalAlarm |

TABLE C

| Index | Status Group | User DO Screen Display |
|-------|-------------------------|------------------------|
| 25 | In Remote | InRem |
| 27 | AutoRange | AutoR |
| 28 | Range 1 | R1 |
| 29 | Range 2 | R2 |
| 30 | Range 3 | R3 |
| 31 | Range 4 | R4 |
| 32 | In Calibrate | InCal |
| 33 | In Zero | Zero |
| 34 | In Span | Span |
| 35 | In Sample | Sample |
| 36 | In NO Mode | InNO |
| 37 | In NOx Mode (605 only) | InNOx |
| 38 | In Wet Mode (HCLD only) | InWet |
| 39 | In Overflow (not used) | InOflow |
| 40 | In NH3 Mode (605 only) | InNH3 |

TYPICAL RELAY

These contacts (3, 4) will drive continuously up to 500 MA using a customer voltage supply that does not to exceed 60 VDC.

- **OPERATION**

Use **(Main, F5, F9)** to select the first seven outputs.

Use the \uparrow to select the desired output.

Press ENTER and use \uparrow to select desired item.

Press ENTER to save selection

Note: The 600 CLD has 15 user selectable isolated digital outputs from the list of 40 in **TABLE B & C**

| User | DO | I | 4: SENO | |
|--------------------------|---------|------|----------------|--|
| 1 | Off | | | |
| 2 | Off | | | |
| 3 | Off | | | |
| 4 | Off | | ***** | |
| 5 | Off | | *Set Unused * | |
| 6 | Off | | *Channels Off* | |
| 7 | Off | | ***** | |
| F2 | Hold | / | Clear | |
| F1 | 8 to 15 | DO's | NO: 0.1 | |
| MAIN/BACK | to | SAVE | O2: 0.0 % | |
| Mon Jan 01 01:40:39 2001 | | | | |

Press **F1** to observe the remaining eight outputs
Program as desired per the above

| User | DO | II | 5:SDRY | ⏏ |
|-----------|--------|----------|----------------|---|
| 8 | Off | | | |
| 9 | Off | | | |
| 10 | Off | | | |
| 11 | Off | | ***** | |
| 12 | Off | | *Set Unused* | |
| 13 | Off | | *Channels Off* | |
| 14 | Off | | ***** | |
| 15 | Off | | | |
| F1 | 1 to 7 | DO's | NO: 0.1 | |
| MAIN/BACK | to | SAVE | O2: 0.0 % | |
| CAI Tel. | (800) | 959 0959 | | |

3.7 Cal Analog Output: (Main, F5, F8,)

| | |
|---------------------------|---------------------|
| Measure setup 4: SENO | |
| F1 | Lowpass filter Tc |
| F2 | Purge time |
| F3 | T + P compensation |
| F4 | Set dual-mode times |
| F8 | Cal Analog Out On |
| F9 | Ignite on power up |
| set analog out during cal | |
| NO: 42.5 | |
| Tue May 15 08:47:49 2007 | |

Use F8 to toggle on/off

Normally, during Auto Cal the Sample and Hold Outputs NO, NOx and NO2 are held at the last process value. If Cal Analog Output is set "On" then the values are not held, and the Real Time value is Output.

4.0 CHANGES TO EXISTING FUNCTIONS

4.1 Saved or Outside Limits

During Manual Calibration the following screens will be displayed to indicate the instruments response to the value of the zero or span gas using the amount that the operator defined in the deviation table.

| | |
|--|-----------|
| Manual Zero 3: SMAN | |
| Zero gas. Use F3 To Toggle Diagnostics Cal Via *Using Current Temp & Press* ***Saved Current*** | |
| NOx - 0.00 | |
| Range1: 100.0ppm | |
| F1 Save value | NOx: -0.0 |
| F2 Man. Cal. Menu | O2: 0.0 % |
| Thu May 15 06:16:17 2008 | |

| | |
|---|-----------|
| Manual Span 1: SARA | |
| Span gas. Use F3 To Toggle Diagnostics Cal Via *Using Current Temp & Press* Outside Deviation Limits | |
| NOx 0.08 | |
| Range1: 100.0ppm Conc.: 90.000000 | |
| F1 Save value | NOx: 0.1 |
| F2 Man. Cal. Menu | O2: 0.0 % |
| Thu May 15 06:18:22 2008 | |

The above is shown using Zero Gas

From Measurement use: **F5 “Zero” or F6 “Span”**

From Main Menu use: **F4, F2, F1 “Zero” or F2 “Span”**

4.2 Calibration Deviations.

MAIN, F5, F2, F2 Deviations, F3 Measuring Deviations.

| Setup 5: SWET | | | Setup 1: SARA | | |
|---------------|----------|----------|------------------------|------|--|
| Deviations[%] | | | Measuring deviation[%] | | |
| | absolute | relative | | | |
| Range1: | 10.00 | 10.00 | Range1: | 1.00 | |
| Range2: | 10.00 | 10.00 | Range2: | 1.00 | |
| Range3: | 10.00 | 10.00 | Range3: | 1.00 | |
| Range4: | 10.00 | 10.00 | Range4: | 1.00 | |
| NOx: -0.0 | | | NOx: 0.0 | | |

Note: These screens are used by the operator to define the maximum acceptable limits of the Zero and Span gas for both Manual and Automatic Calibrating.

4.3 Flow Zero or Span

Some analyzers have the above and the ability to flow Zero and Span Gas.

| Manual calib 3: SMAN | | | Manual calib 1: SARA | | |
|-------------------------|---------------|--|--------------------------|---------------|--|
| F1 Cal zero | 0.08 | | F1 Cal zero | 0.08 | |
| F2 Cal span | ZERO | | F2 Cal span | SPAN | |
| F3 Toggle zero | | | F3 Toggle zero | | |
| F4 Toggle span | | | F4 Toggle span | | |
| F5 Measurements | | | F5 Measurements | | |
| F6 Cal zero 02 | | | F6 Cal zero 02 | | |
| F7 Cal span 02 | | | F7 Cal span 02 | | |
| F8 Range Up ** | F9 Range Down | | F8 Range Up ** | F9 Range Down | |
| Range1: 100.0ppm | NOx: 0.1 | | Range1: 100.0ppm | NOx: 0.1 | |
| | O2: 0.0 % | | | O2: 0.0 % | |
| CAI Tel. (800) 959 0959 | | | Thu May 15 06:22:23 2008 | | |

The above is shown using Zero/SPAN Gas

From Cal Screen use: **F2** or **Main** or **Back**

From Main Menu use: **F4, F2**

4.4 Reset Calibration Values

When the re-set calibrations value function is used all recorded deviations will be set to zero

| Def cal value | |
|---|-----|
| Are you sure to reset calibration values to default values? | |
| F1 | yes |
| F2 | no |
| NOx: -0.0 | |
| Mon Jan 01 00:01:36 2001 | |

Main, F4,F5

| | | | | |
|--------------------------|----------|------------|------|---|
| Span gas devs | | | | ⌂ |
| Span gas devs. verifying | | | | |
| | measured | deviations | | |
| | values | Var | %FS | |
| Range1: | 0.0 | 0.00 | 0.00 | |
| Range2: | 0.0 | 0.00 | 0.00 | |
| Range3: | 0.0 | 0.00 | 0.00 | |
| Range4: | 0.0 | 0.00 | 0.00 | |
| | | NOx: -0.0 | | |
| Mon Jan 01 00:02:08 2001 | | | | |

Main,F4 F3, F4

(Used to observe Auto Cal Results)

| Span gas devs | | |
|--------------------------|------|-----------|
| Span gas deviations [%]: | | |
| | abs | rel |
| Range1: | 0.00 | 0.00 |
| Range2: | 0.00 | 0.00 |
| Range3: | 0.00 | 0.00 |
| Range4: | 0.00 | 0.00 |
| | | NOx: -0.0 |
| CAI Tel. (800) 959 0959 | | |

Main, F4, F3, F2

(Used to observe Manual Cal results)

The above are the new deviations after the operator elects to re-set the calibration values

TABLE D**600 SERIES CLD IO CHART****28 PIN MAIN CONNECTOR ASSIGNMENTS**

AO = Analog Output, OC= Open Collector, SV = Solenoid Valve TTL = Transistor Logic

| OPTO I/O | Signal Type | 600 CLD/HCLD Analog | | Levels |
|---------------------|------------------------|--------------------------------|------------------------------|----------------------------------|
| ALG 1 | | pin # | | |
| COM | A Output | 1 | GND (Isolated analog) | <u>Isolated AI</u> |
| 0 | A Output | 2 | User Defined AO-1 | 1v,5v,10v,mA |
| 1 | A Output | 3 | User Defined AO-2 | 1v,5v,10v,mA |
| 2 | A Output | 4 | User Defined AO-3 | 1v,5v,10v,mA |
| 3 | A Output | 5 | User Defined AO-4 | 1v,5v,10v,mA |
| DIG 1 | | Digital | | |
| COM | D Output | 6 | GND (Digital) | |
| 0 | D Output | 7 | Sense Auto Range | TTL-low true |
| 1 | D Output | 8 | Sense Range 1 | TTL-low true |
| 2 | D Output | 9 | Sense Range 2 | TTL-low true |
| 3 | D Output | 10 | Sense Range 3 | TTL-low true |
| | D Output | 11 | Sense Range 4 | TTL-low true |
| 5 | D Input | 12 | Set Auto Range | |
| 6 | D Input | 13 | Control Range 1 | |
| 7 | D Input | 14 | Control Range 2 | |
| 8 | D Input | 15 | Control Range 3 | |
| 9 | D Input | 16 | Control Range 4 | |
| 10 | D Input | 17 | Auto Cal | |
| 11 | D Input | 18 | Calibrate | |
| 12 | D Input | 19 | Zero | |
| 13 | D Input | 20 | Span | |
| 14 | D Input | 21 | Sample | |
| 15 | SPARE | | | |
| DIG 2 | | | | |
| 0 | D Output | 22 | Zero Gas Flow | OC (24vdc if internal SV) |
| 1 | D Output | 23 | Span Gas Flow | OC (24vdc if internal SV) |
| 2 | D Output | 24 | Sample Gas Flow | OC (24vdc if internal SV) |
| 3 | D Output | 25 | Local/Remote | TTL-low true |
| 4 | D Output | 26 | Read Cal Mode | TTL-low true |
| 5 | D Output | 27 | Reserved | |
| 6 | D Output | 28 | Reserved | |

TABLE D (CONT)**600 SERIES CLD IO CHART****28 PIN AUXILLIARY CONNECTOR ASSIGNMENTS****NO = Normally Open**


| OPTO I/O | Signal Type | CLD Analog | | LEVELS |
|------------------|----------------|----------------|---------------------------------|---|
| ALG 1 | Spare | pin # | | |
| COM | A Input | 1 | GND (analog) | |
| 4 | A Input | 2 | External Analog 1 | 0-10V |
| 5 | A Input | 3 | External Analog 2 | 0-10V |
| 6 | A Output | 4 | GND (Isolated analog) | |
| 7 | D Output | 5 | Relay RTN 1 | 9,10,11,12 use RTN 1 |
| DIG 3 | Alarms | Digital | | <u>Status go CLOSED when active</u> <u>Alarms go OPEN when present</u> |
| COM | D Output | 6 | Relay RTN 2 | 13,14,15,16 use RTN 2 17,18,19,20 use RTN 3 21,27,28 use RTN 4 |
| 0 | D Output | 7 | Relay RTN 3 | |
| 1 | D Output | 8 | Relay RTN 4 | |
| 2 | D Output | 9 | User Defined NO Relay 1 | |
| 3 | D Output | 10 | User Defined NO Relay 2 | |
| 4 | D Output | 11 | User Defined NO Relay 3 | |
| 5 | D Output | 12 | User Defined NO Relay 4 | |
| 6 | D Output | 13 | User Defined NO Relay 5 | |
| 7 | D Output | 14 | User Defined NO Relay 6 | |
| 8 | D Output | 15 | User Defined NO Relay 7 | |
| 9 | D Output | 16 | User Defined NO Relay 8 | |
| 10 | D Output | 17 | User Defined NO Relay 9 | |
| 11 | D Output | 18 | User Defined NO Relay 10 | |
| 12 | D Output | 19 | User Defined NO Relay 11 | |
| 13 | D Output | 20 | User Defined NO Relay 12 | |
| 14 | D Output | 21 | User Defined NO Relay 13 | |
| 15 | D Output | 22 | Reserved Do Not Connect | |
| DIG 2 | | | | |
| 7 | D Input | 23 | Set Wet Mode | |
| 8 | D Input | 24 | Set Overflow | |
| 9 | D Input | 25 | Set NO Mode | |
| 10 | D Input | 26 | Set Remote | |
| 11 | D Output | 27 | User Defined NO Relay 14 | |
| 12 | D Output | 28 | User Defined NO Relay 15 | |

Communication is via a Serial or Ethernet (TCP/IP) interface. The Ethernet interface can carry streamed UDP data via the TCP/IP port to a separate configured UDP port.

1.1. *Serial Interface and AK-Commands*

The serial interface enables remote control of the analyzer by a master computer. It is implemented as an RS232 V24 interface and meets all requirements of the AK protocol.

A 9-pin male connector at the back of the unit is used to connect a master computer with the following pin assignment:

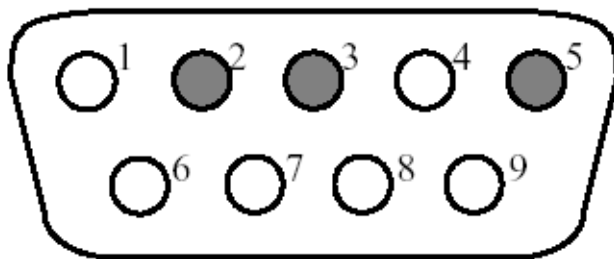

Pin 2 = Rxd (receive)

Pin 3 = Txd (transmit)

Pin 5 = Gnd (ground)

Figure Error! No text of specified style in document.-1 **Serial Interface**

1.2. *Interface Specifications*

| | |
|------------------|-------------|
| Speed: | 9600 bps |
| Character Length | 1 start bit |
| | 8 data bits |
| | 1 stop bit |
| Parity: | none |
| Handshake | no |

Ethernet. RJ45.

If connecting directly to a computer (without a Hub or Switch) a CROSSOVER Cable is required.

Protocol Description

Instruction command

| | Character | Explanation |
|-----------------------|-----------------|--|
| 1 st Byte | STX | ASCII code 02 |
| 2 nd Byte | Don't Care | Any ASCII code |
| 3 rd Byte | Function Code 1 | AK instruction e.g.: ASTF |
| 4 th Byte | Function Code 2 | |
| 5 th Byte | Function Code 3 | |
| 6 th Byte | Function Code 4 | |
| 7 th Byte | Blank | |
| 8 th Byte | K | |
| 9 th Byte | 0 | |
| 10 th Byte | Blank | |
| | D | AK instruction parameters, length is variable |
| | A | |
| | T | |
| | A | |
| nth Byte | ETX | ASCII code 03 |

Table Error! No text of specified style in document.-

1 Structure of an instruction command

Acknowledgement command

| | Character | Explanation |
|-----------------------|-----------------|---|
| 1 st Byte | STX | ASCII code 02 |
| 2 nd Byte | Don't Care | Any ASCII code |
| 3 rd Byte | Function Code 1 | Echo of the AK instruction command |
| 4 th Byte | Function Code 2 | |
| 5 th Byte | Function Code 3 | |
| 6 th Byte | Function Code 4 | |
| 7 th Byte | Blank | |
| 8 th Byte | K | |
| 9 th Byte | 0 | |
| 10 th Byte | Blank | |
| | D | AK acknowledgement parameters, length is variable |
| | A | |
| | T | |
| | A | |
| nth Byte | ETX | ASCII code 03 |

Table Error! No text of specified style in document.

-2 Structure of an acknowledgement command

Data Description

Each command begins with STX (Start of Text) in the first byte. The "don't care" byte can be any ASCII character. Generally, a blank or an underscore () is used for readability reasons. The four function bytes represent the AK command. A blank comes next, followed by K and the channel number. The analyzer is a single-channel device, and because of that the channel number is usually 0. For delimiting the command parameters from the channel number, another blank follows. This may be followed by command parameters with a variable length. Every command ends with the ETX (End of Text) character. The error status byte in the acknowledgment command signals if internal errors in the analyzer occurred. It is zero when no error appeared, and it is unequal zero when one or more errors occurred. Every time a change in the errors happens the error status byte is incremented by one, no matter if one or several errors disappear or are added. If it had the value 10, it would be reset to 1. The error status byte does not indicate the real number of errors. If the analyzer does not have errors, the error status byte contains the value 0.

In general, AK commands are subdivided into three classes:

- Control commands (Sxxx)
- Inquiry commands (Axxx)
- Configuration commands (Exxx)

Error Handling

It might happen that an unknown instruction is sent, that the analyzer is busy with a function which is not the desired one, or that an error occurred in the command parameters. Table 1-4 summarizes all errors that can appear upon any master instruction.

| Analyzers Acknowledgement | Explanation |
|----------------------------------|---|
| ???? f | Analyzer does not know the instruction sent. |
| xxxx f BS | Analyzer is busy with another function. |
| xxxx f SE | Syntax error within command parameters or incomplete command. |
| xxxx f NA | Requested function or data are not available. |
| xxxx f DF | Data error: The kind or number of given parameters are not valid. |
| xxxx f OF | Offline: Analyzer is offline, i.e. analyzer is in local Mode. Only inquiry commands and SREM (set analyzer in remote mode) are allowed. |

f is the error status byte.

xxxx is the function code of the command sent.

1.3. Scan Commands

1.3.1. AKON: Measured concentration value

| Command | Response | Description |
|----------------------|---|---|
| _AKON_K0 | _AKON_s_z.z_y.y_x.x_t _AKON_s_z.z_y.y_x.x_w.w_v.v_t _AKON_s_z.z_y.y_x.x_w.w_t | <u>NDIR</u> Measured concentration value of all 3 possible channels. <u>CLD / HCLD</u> z.z : current Measure Value. y.y : NO concentration. x.x : NO2 concentration. w.w : NOx concentration v.v : NH3 concentration Note. y.y, x.x, w.w are normally 0.0 when "Dual Measure" or "Analog Out during Cal" are not selected. v.v applies only to the NH3 option. <u>FID / HFID</u> z.z : current Measure Value. y.y : CH4 concentration. x.x : NMHC concentration. w.w : THC concentration Note. y.y, x.x, w.w are normally 0.0 when "Dual Measure" or "Analog Out during Cal" are not selected. t = Timestamp (1/10 sec) z.z : current Measure Value for channel m. <u>CLD / HCLD</u> channel 1 is O2 option. <u>CLD / FID / HCLD / HFID</u> . Current Measure Value for all 4 ranges based on their last calibration Offset and Gain. t = Timestamp (1/10 sec) |
| _AKON_Km | _AKON_s_z.z_t | |
| _AKON_K4 _AKON_K5 | _AKON_s_z.z_y.y_x.x_w.w_t | |

1.3.2. ARMU: Raw Engineering value

| Command | Response | Description |
|----------|-----------------------|--|
| _ARMU_K0 | _ARMU_s_z.z_y.y_x.x_t | Raw Engineering value before linearization and offset-span-correction for all 3 possible channels. These are the values used to calculate the polynomial coefficients. t = Timestamp (1/10 sec) |
| _ARMU_K0 | _ARMU_s_z.z_t | Raw Engineering value before linearization and offset-span-correction for channel m. These are the values used to calculate the polynomial coefficients. t = Timestamp (1/10 sec) |

1.3.3. ARAW: Raw Detector Volts

| Command | Response | Description |
|-----------------|-----------------------|---|
| _ARAW_K0 | _ARAW_s_z.z_y.y_x.x_t | Raw Detector Volts for all 3 possible channels. t = Timestamp (1/10 sec) |
| _ARAW_K0 | _ARAW_s_z.z_t | Raw Detector Volts for channel m t= Timestamp (1/10 sec) |

1.3.4. AEMB: Get measuring range

| Command | Response | Description |
|----------|-------------------|---|
| _AEMB_K0 | _AEMB_s_Mn_Mn_Mn | NDIR Current measuring range n of all 3 Channels. |
| | _AEMB_s_Mn | CLD / FID / HCLD / HFID Current range n. |
| _AEMB_Km | _AEMB_s_Mn | Current measuring range of channel m is responded |

1.3.5 AMBE: Measuring range limit

| Command | Response | Description |
|-------------|---|--|
| _AMBE_Km | _AMBE_s_M1_w.w _M2_x.x _M3_y.y _M4_z.z | All existing measuring range limits of channel m. NDIR m = 1, 2 or 3. CLD / FID / HCLD / HFID m = 0. |
| _AMBE_Km_Mm | _AMBE_s_Mn_z.z | Range limit n of channel m. n = 1 to 4. |

1.3.6. AKAK: Calibration gas concentrations

| Command | Response | Description |
|-------------|---|---|
| _AKAK_Km | _AKAK_s_M1_w.w _M2_x.x _M3_y.y _M4_z.z | All existing calibration gas values are responded for selected channel m NDIR m = 1, 2 or 3. CLD / FID / HCLD / HFID m = 0. |
| _AKAK_Km_Mn | AKAK_s_Mn_z.z | channel m calibration gas value of Range n. |

1.3.7. AMBU: Upper and lower range switchover values for auto range

| Command | Response | Description |
|-------------|---|---|
| _AMBU_Km | _AMBU_s_M1_w.w_W.W _M2_x.x_X.X _M3_y.y_Y.Y _M4_z.z_Z.Z | Lower and upper range switchover value of auto range for channel m NDIR m = 1, 2 or 3. CLD / FID / HCLD / HFID m = 0. |
| _AMBU_Km_Mn | _AMBU_s_Mn_w.w_W.W | Lower and upper range switchover value of auto range for channel m range n |

1.3.8 ASTZ: Normal device status

| Command | Response | Description |
|----------|---|--|
| _ASTZ_K0 | _ASTZ_s_K1_State1_State2_State3 _K2_State1_State2_State3 _K3_State1_State2_State3 | Respond device status for all channels |
| _ASTZ_Km | _ASTZ_s_State1_State2_State3 | Respond device status only for channel m |

Possible states:

NDIR.

| State 1 | State 2 | State 3 |
|------------------------------|---|---|
| SREM: remote SMAN: manual | STBY: standby SPAU: pause SMGA: measuring gas SNGA: zero gas SEGA: end gas SATK SNGA: zero gas during auto cal SATK SEGA: end gas during auto cal | SARE: auto range on SARA: auto range off |

CLD / HCLD.

| State 1 | State 2 | State 3 |
|------------------------------|---|---|
| SREM: remote SMAN: manual | STBY: standby SPAU: pause SMGA: measuring gas SNGA: zero gas SEGA: span gas SATK SNGA: zero gas during auto cal SATK SEGA: end gas during auto cal SSPL: Purging Overflow | SENO: NO mode SNOX: NOx mode SNO2: dual mode?? SNH3: triple mode. NH3 option ?? |

State 4
State 5
SARE: auto range on
SDRY: Chiller in
SARA: auto range off
SWET: Chiller out

FID / HFID.

| State 1 | State 2 | State 3 |
|------------------------------|---|---|
| SREM: remote SMAN: manual | STBY: standby SPAU: pause SMGA: measuring gas SNGA: zero gas SEGA: span gas SATK SNGA: zero gas during auto cal SATK SEGA: end gas during auto cal SSPL: Purging Overflow | SHCG: THC mode SCH4: CH4 mode SMNH: dual mode?? |

State 4
State 5
SARE: auto range on
?????
SARA: auto range off
?????

1.3.9. ASTF: Error status

| Command | Response | Description |
|----------|-------------------------|--|
| _ASTF_K0 | _ASTF_s_f1_f2_f3..._f10 | Current error numbers of all are responded |

Errors:

NDIR

| | | | |
|---|---------------------------|-------|-------------------------------------|
| 1 | Channel 1 Flow Failure | 8 | Channel 1 not calibrated |
| 2 | Channel 2 Flow Failure | 9 | Channel 2 not calibrated |
| 3 | Channel 3 Flow Failure | 10 | Channel 3 not calibrated |
| 4 | External Analog 1 Failure | 11-13 | Ch1...3: Low concentration warning |
| 5 | External Analog 2 Failure | 14-16 | Ch1...3: High concentration warning |
| 6 | Pressure Failure | 17-19 | Ch1...3: Temperature failure |
| 7 | Temperature Failure | 20-22 | Ch1...3: EPC Voltage failure |

CLD / HCLD

| | |
|----|-------------------------------|
| 1 | Sample Pressure Failure |
| 2 | Air Pressure Failure |
| 3 | Oven Temp. Failure |
| 4 | Converter Temp. Failure |
| 5 | Pump Temp. Failure |
| 6 | Diode Temp. Failure |
| 7 | Cell Temp Failure |
| 8 | Peltier Gas Temp Failure |
| 9 | Reaction Chamber Temp Failure |
| 10 | EPC Coil Sample Failure |
| 11 | EPC Coil Air Failure |
| 12 | Range Overflow |
| 13 | ADC Range Overflow |
| 14 | ADC Range Underflow |
| 15 | Range 1 is not calibrated |
| 16 | Range 2 is not calibrated |
| 17 | Range 3 is not calibrated |
| 18 | Range 4 is not calibrated |

FID / HFID

| | | | |
|---|-------------------------|----|----------------------------|
| 1 | No Flame | 9 | EPC Coil Sample Failure |
| 2 | Sample Pressure Failure | 10 | EPC Coil Air Failure |
| 3 | Air Pressure Failure | 11 | EPC Coil Fuel Failure |
| 4 | Fuel Pressure | 12 | Range Overflow |
| 5 | Burner Temp Failure | 13 | ADC Range Overflow |
| 6 | Oven Temp Failure | 14 | ADC Range Underflow |
| 7 | Cutter Temp Failure | 15 | Analyzer is not calibrated |
| 8 | Pump Temp Failure | | |

1.3.10. AKEN: Device identification

| Command | Response | Description |
|----------|------------------------|--|
| _AKEN_K0 | _AKEN_s_devicename | Device identification is responded |
| _AKEN_K1 | _AKEN_s_model | Device model |
| _AKEN_K2 | _AKEN_s_serial no | Device serial number |
| _AKEN_K3 | _AKEN_s_samplepressure | Sample pressure |
| _AKEN_K4 | _AKEN_s_Air pressure | Air Pressure (Not NDIR) |
| _AKEN_K5 | _AKEN_s_Fuel pressure | Fuel Pressure (Not NDIR or CLD / HCLD) |

1.3.11. ATEM: Temperatures

| Command | Response | Description |
|----------|--------------------------|---|
| _ATEM_K0 | _ATEM_s_z.z_y.y_x.x..... | <u>NDIR</u> Detector temperature of all 3 possible channels <u>CLD / HCLD</u> 1. Oven Temp 2. Converter Temp 3. Pump Temp 4. Diode Temp 5. Cell Temp 6. Chiller Temp 7. O2 Detector Temp 8. NH3 Converter Temp <u>FID / HFID</u> 1. Filter Temp 2. Burner Temp 3. Oven Temp 4. Cutter Temp 5. Pump Temp |
| _ATEM_Km | _ATEM_s_z.z | <u>NDIR</u> Detector temperature of channel m <u>CLD / HCLD / FID / HFID</u> Temperature of channel m as above. |

1.3.12. ADRU: Pressures/ Electronic Pressure Control Valve voltage

| Command | Response | Description |
|----------|--------------------------|--|
| _ADRU_K0 | _ADRU_s_z.z_y.y_x.x..... | All Pressures / Voltages. <u>NDIR</u> 1. Sample Pressure Ch1 2. Sample Pressure Ch2 3. Sample Pressure Ch3 4. EPC Volts Ch1 5. EPC Volts Ch2 6. EPC Volts Ch2 <u>CLD / HCLD</u> 1. Sample Pressure 2. Air Pressure 3. Sample EPC Volts 4. Ozone EPC Volts <u>FID / HFID</u> 1. Sample Pressure 2. Burner Air Pressure 3. Burner Fuel Pressure 4. Air Inject Pressure 5. Fuel Inject Pressure 6. Sample EPC Volts 7. Burner Air EPC Volts 8. Burner Fuel EPC Volts 9. Air Inject EPC Volts 10. Fuel Inject EPC Volts |
| _ADRU_Km | _ADRU_s_z.z | EPC voltage of channel m is returned in z.z |

1.3.13. ADUF: Flows

| Command | Response | Description |
|----------|---------------------|---|
| _ADUF_K0 | _ADUF_s_z.z_y.y_x.x | <u>NDIR</u> Sample gas flow of all 3 channels <u>CLD / HCLD</u> 1. Sample Flow 2. Air Flow <u>FID / HFID</u> 1. Sample Flow 2. Air Flow 3. Fuel Flow 4. Air Inject Flow???? 5. Fuel Inject Flow???? |
| _ADUF_Km | _ADUF_s_z.z | Sample gas flow of channel m |

1.3.14. AGRD: Polynomial coefficients

| Command | Response | Description |
|-------------|---------------------------|--|
| _AGRD_Km_Mn | _AGRD_s_Mn_a0_a1_a2_a3_a4 | Polynomial coefficients of channel m range n <u>NDIR</u> m = 1, 2 or 3 <u>CLD / HCLD</u> m = 0 m = 1 optional O2 coefficients <u>FID / HFID</u> m = 0 |

1.3.15. AKAL: Percent Deviations of last accepted calibration

| Command | Response | Description |
|-----------|---|--|
| _AKAL_Km_ | _AKAL_s_M1_z.z_y.y_x.x_w.w _M2_z.z_y.y_x.x_w.w _M3_z.z_y.y_x.x_w.w _M4_z.z_y.y_x.x_w.w | Percent Deviation of Ranges M1 to M4: z.z: Zero gas relative to last calibration y.y: Zero gas absolute to factory calibration. x.x: Span gas relative to last calibration. w.w: Span gas absolute to factory calibration <u>NDIR</u> m = 1, 2 or 3 <u>CLD / HCLD / FID / HFID</u> m = 0 |

1.3.16. AAOG: Applied Offset and Gains

| Command | Response | Description |
|-----------|---|--|
| _AAOG_Km_ | _AAOG_s_M1_z.z_y.y _M2_z.z_y.y _M3_z.z_y.y _M4_z.z_y.y _O2_z.z_y.y (CLD / HCLD) | Offset and Gain of Ranges M1 to M4: z.z: Offset y.y: Gain <u>NDIR</u> m = 1, 2 or 3 <u>CLD / HCLD / FID / HFID</u> m = 0 |

1.3.17. AANG: Verifying zero point deviation during auto calibration

| Command | Response | Description |
|----------|---|---|
| _AANG_Km | _AANG_s_M1_z.z_da_dr _M2_z.z_da_dr _M3_z.z_da_dr _M4_z.z_da_dr | Verifying deviations of Ranges M1 to M4 from zero point stored after auto calibration. Values: measured value (z.z), absolute dev (da), relative dev (dr) <u>NDIR</u> <u>m = 1, 2 or 3</u> <u>CLD / HCLD / FID / HFID</u> <u>m = 0</u> |

1.3.18. AAEG: Verifying span point deviation during auto calibration

| Command | Response | Description |
|----------|---|--|
| _AAEG_Km | _AANG_s_M1_z.z_da_dr _M2_z.z_da_dr _M3_z.z_da_dr _M4_z.z_da_dr | Verifying deviations of Ranges M1 to M4 from span point stored after auto calibration Values: measured value (z.z), absolute dev (da), relative dev (dr) <u>NDIR</u> <u>m = 1, 2 or 3</u> <u>CLD / HCLD / FID / HFID</u> <u>m = 0</u> |

1.3.19. AFDA: Auto calibration times and Purge time

| Command | Response | Description |
|---------------|-----------------|---|
| _AFDA_Km_SATK | _AFDA_s_z_y_x_w | Auto calibration times of channel m: z: Purge time y: Calibration time x: Total calibration time w: Verify time (z, y, x, w in seconds) <u>NDIR</u> <u>m = 1, 2 or 3</u> <u>CLD / HCLD / FID / HFID</u> <u>m = 0</u> |
| _AFDA_K0_SSPL | _AFDA_s_z | Purge time will be responded |

1.3.20. APAR: Auto calibration tolerance values

| Command | Response | Description |
|---------------|-------------------------|--|
| _APAR_Km_SATK | _APAR_s_z.z_y.y_x.x_w.w | During Autocal if all the values used in the average are within tolerance then Autocal moves to Verify. Auto calibration tolerance value (%): z.z: Range 1 y.y: Range 2 x.x: Range 3 w.w: Range 4 NDIR m = 1, 2 or 3 <u>CLD / HCLD / FID / HFID</u> m = 0 |

1.3.21. ASYZ: Respond System Time

| Command | Response | Description |
|-----------|-----------------------|--|
| _ASYZ_K0_ | _ASYZ_s_yymmdd_hhmmss | System time: yymmdd: year, month, day (each 2 characters wide, no spaces) hhmmss: hour, minute, second (each 2 characters wide, no spaces) |

1.3.22. AT90: Respond Lowpass filter time

| | Response | Description |
|-----------|---------------|--|
| _AT90_Km_ | _AT90_s_t.... | Low pass filter time in seconds of channel m T = filter time in seconds of channel m NDIR m = 0, 1, 2 or 3 if m = 0 then all channels. <u>CLD / HCLD / FID / HFID</u> m = 0 |

1.3.23. ADAL: Diagnostic alarm limits

| Command | Response | Description |
|------------|---|-------------------|
| _ADAL_K0 | _ADAL_s_a1.min_a1.max_... _a16.min_a16.max | All alarms limits |
| _ADAL_K0_x | _ADAL_s_x.min_x.max | Alarm limits of x |

Alarm Limits:**NDIR**

| | | | |
|---|---------------------------|-------|--------------------------------|
| 1 | Sample gas flow channel 1 | 7 | Temperature |
| 2 | Sample gas flow channel 2 | 8 | Sample concentration channel 1 |
| 3 | Sample gas flow channel 3 | 9 | Sample concentration channel 2 |
| 4 | External input 1 | 10 | Sample concentration channel 3 |
| 5 | External input 2 | 11-13 | Temperature channel 1...3 |
| 6 | Barometric -Pressure | 14-16 | EPC voltage channel 1...3 |

CLD / HCLD

| | |
|----|----------------------------|
| 1 | Sample Pressure |
| 2 | Air Pressure |
| 3 | Oven Temp. |
| 4 | Converter Temp. |
| 5 | Pump Temp. |
| 6 | Diode Temp. |
| 7 | Cell Temp |
| 8 | Peltier Gas Temp |
| 9 | EPC Coil Sample Voltage |
| 10 | EPC Coil Air/Ozone Voltage |
| 11 | Reserved |
| 12 | Sample Content |

FID / HFID

| | | | |
|---|--------------------|----|-------------------------|
| 2 | Sample Pressure | 8 | Pump Temperature |
| 3 | Air Pressure | 9 | Sample EPC Coil Voltage |
| 4 | Fuel Pressure | 10 | Air EPC Coil Voltage |
| 5 | Burner Temperature | 11 | Fuel EPC Coil Voltage |
| 6 | Oven Temperature | 12 | Sample Content |

1.3.24. AVER: Query Software version

| Command | Response | Description |
|----------|-------------------------------------|--|
| _AVER_K0 | _AVER_s_MAIN_z_3USER_y_OS MSR_x. | z:Main versionx.xxxx_ buildno_builddate_dd.mm.yyyy y: User version x.xxx_ buildno_builddate_dd.mm.yyyy x: OSMSR version x.xxx_ builddate_dd.mm.yyyy |

1.3.25. ATCP: Query TCP/IP settings

| Command | Response | Description |
|----------|--|---|
| _ATCP_K0 | _ADAL_s_zzz.zzz.zzz.zzz _yyy.yyy.yyy.yyy _xxxx | zzz: TCP/IP Address yyy: TCP/IP subnet mask xxxx: TCP/IP port |

1.3.26. AUDP: Query UDP data streaming parameter

| Command | Response | Description |
|----------|---|--|
| _AUDP_K0 | _AUDP_s_<UDPPort>_ <DataFrequency>_[<Mode>] _[<UDP_IP>] | UDP port: opened for connection DataFrequency: Transmission Frequency of the data in Hz Mode: A: ASCII Mode UDP_IP: Alternative IP address open for the UDP connection when it should not use the IP connected to the TCP/IP client |

1.3.27. AH2O: QueryH₂O correction parameter **NDIR**

| Command | Response | Description |
|----------|------------------------|--|
| _AH2O_Km | _AH2O_s_z.z.y.y_x.x_ww | m = Channel 1 to 3 z.z: Ext. analog 2 value y.y: Dry – voltage of A in with no water present x.x.: 1 st order coefficient w.w: 2nd order coefficient |

1.3.28. ACO2: Query CO₂ correction parameter **NDIR**

| Command | Response | Description |
|----------|-----------------------------|--|
| _ACO2_Km | _ACO2_s_z.z.y.y_x.x_w.w_v.v | m = Channel 1 to 3 z.z: Ext analog 1 value y.y: Offset – voltage of A in with no CO ₂ present x.x: Min A in – if A in is below this value no CO ₂ correction will be done. w.w: 1st order coefficient v.v.: 2nd order coefficient |

1.3.29. AAUX: 4 channels of 0 – 10 Volt Auxiliary Option **CLD / HCLD**

| Command | Response | Description |
|----------|-------------------------|---|
| _AAUX_K0 | _AAUX_s_z.z.y.y_x.x_w.w | z.z: Aux Channel 1 Volts y.y: Aux Channel 2 Volts x.x: Aux Channel 3 Volts w.w.: Aux Channel 4 Volts |

1.4. Control commands

1.4.1. SRES: Reset

| Command | Response | Description |
|----------|----------|-------------|
| _SRES_K0 | _SRES_s | Reset |

1.4.2. SPAU: Pause

| Command | Response | Description |
|----------|----------|-------------|
| _SPAU_K0 | _SPAU_s | Pause mode |

1.4.3. STBY: Standby

| Command | Response | Description |
|----------|----------|------------------------------------|
| _STBY_K0 | _STBY_s | Standby mode for CLD & FID |
| | | Standby mode for all channels NDIR |
| _STBY_Km | _STBY_s | Standby mode for channel m NDIR |

1.4.4. SNGA: Open valve for zero gas calibration

| Command | Response | Description |
|-------------|----------|--|
| _SNGA_K0 | _SNGA_s | Open Zero gas valve CLD & FID |
| | | Open all 3 Zero gas valves NDIR |
| _SNGA_Km | _SNGA_s | Open Zero gas valve of channel m NDIR |
| _SNGA_Km_Mn | _SNGA_s | Open Zero gas valve and set range to ' n ' CLD & FID |
| | | Open Zero gas valve of channel m and set range to ' n ' NDIR |

1.4.5. SEGA: Open valve for end gas calibration

| Command | Response | Description |
|-------------|----------|--|
| _SEGA_K0 | _SEGA_s | Open Span gas valve CLD & FID |
| | | Open all 3 Span gas valves NDIR |
| _SEGA_Km | _SEGA_s | Open Span gas valve of channel m NDIR |
| _SEGA_Km_Mn | _SEGA_s | Open Span gas valve and set range to ' n ' CLD & FID |
| | | Open Span gas valve of channel m and set range to ' n ' NDIR |

1.4.6. SSPL: Purge Analyzer with zero gas

| Command | Response | Description |
|----------|----------|---|
| _SSPL_K0 | _SSPL_s | Open Purge gas valve CLD & FID |
| | | Open Purge gas valve(s) and purge all Channels NDIR |

1.4.7. SATK: Start automatic calibration

| Command | Response | Description |
|-------------|----------|--|
| _SATK_K0 | _SATK_s | Start automatic calibration of all available ranges CLD & FID |
| | | Start automatic calibration of all available ranges and channels NDIR |
| _SATK_Km | _SATK_s | Not applicable CLD & FID |
| | | Start automatic calibration for channel m of all available ranges NDIR |
| _SATK_Km_Mn | _SATK_s | Start automatic calibration of range ' n ' CLD & FID |
| | | Start automatic calibration for channel m, Range ' n ' NDIR |

1.4.8. SEMB: Set measuring range

| Command | Response | Description |
|---------|----------|------------------------------------|
| | _SEMB_s | m = 0 CLD & FID |
| | | m = channel NDIR |
| | | Set measuring range to range ' n ' |
| | | Auto range will be disabled |

1.4.9. SARE: Auto range on

| Command | Response | Description |
|----------|----------|---|
| _SARE_K0 | _SARE_s | Set auto range CLD & FID |
| | | Set auto range on for all channels NDIR |
| _SARE_Km | _SARE_s | Set auto range on for channel m NDIR |

1.4.10. SARA: Auto range off

| Command | Response | Description |
|----------|----------|--|
| _SARA_K0 | _SARA_s | Set auto range off CLD & FID |
| | | Set auto range off for all channels NDIR |
| _SARA_Km | _SARE_s | Set auto range off for channel m NDIR |

1.4.11. SREM: Remote mode for AK-commands

| Command | Response | Description |
|----------|----------|---------------------------|
| _SREM_K0 | _SREM_s | Set device in remote mode |

1.4.12. SMAN: Manual control to control device manually

| Command | Response | Description |
|----------|----------|---------------------------|
| _SMAN_K0 | _SMAN_s | Set device in manual mode |

1.4.13. SMGA: Start measuring Turn on pumps if fitted

| Command | Response | Description |
|----------|----------|---|
| _SMGA_K0 | _SMGA_s | Open Sample gas valve CLD & FID |
| | | Open all 3 Sample gas valves NDIR |
| _SMGA_Km | _SMGA_s | Open sample gas valve of channel m NDIR |

1.4.14. SNKA: Saves measured value as new offset.

| Command | Response | Description |
|----------|----------|--|
| _SNKA_K0 | _SNKA_s | Saves measured value of actual range as offset if zero valve open CLD & FID K1 on CLD 650 saves measured value as offset of O2 channel if O2 zero valve open |
| _SNKA_Km | _SNKA_s | Saves measured value of actual range for each channel as new offset if zero valve of the channel is open NDIR Saves measured value of actual range of channel m as new offset if zero valve is open |

1.4.15. SEKA: Saves measured value as new span value.

| Command | Response | Description |
|----------|----------|--|
| _SEKA_K0 | _SEKA_s | Saves measured value of actual range as gain if span valve open CLD & FID K1 on CLD 650 saves measured value as gain of O2 channel if O2 span valve open |
| _SEKA_Km | _SEKA_s | Saves measured value of actual range for each channel as new gain if span valve of the channel is open NDIR Saves measured value of actual range of channel m as new gain if span valve is open |

1.4.16. SUDP: Start /Stop UDP data streaming

| Command | Response | Description |
|--------------|----------|--|
| _SUDP_K0_ON | _SUDP_s | Start Data streaming via the UDP channel. You need to configure the channel before with EUDP command |
| _SUDP_K0_OFF | _SUDP_s | Stop streaming via the UDP channel |

1.4.17. SVZS: Reset Offset to 0 and Gain to 1

| Command | Response | Description |
|----------|----------|--|
| _SVZS_K0 | _SVZS_s | Sets all range offsets to 0 and Gains to 1 for CLD & FID Not to be used on NDIR |
| _SVZS_Km | _SVZS_s | NDIR Sets Channel m range Offsets to 0 and Gains to 1 |

1.5. Settings

1.5.1. EKAK: Set the four span gas concentration values

| Command | Response | Description |
|--------------------------------------|----------|--|
| _EKAK_Km_M1_w.w_M2_x.x_M3_y.y_M4_z.z | _EKAK_s | m = 0 Set the span gas values CLD & FID Set the span gas values for channel m |

1.5.2. EMBE: Set the four measuring range full scale limits

| Command | Response | Description |
|--------------------------------------|----------|--|
| _EMBE_Km_M1_w.w_M2_x.x_M3_y.y_M4_z.z | _EMBE_s | m = 0 Set the range full scale limits CLD & FID Set the full scale limits for channel m |

1.5.3. EMBU: Set the upper and the lower range switchover for auto range

| Command | Response | Description |
|--|----------|---|
| _EMBU_Km_M1_w.w_W.W_M2_x.x_X.X_M3_y.y_Y.Y_M4_z.z_Z.Z | _EMBU_s | M = 0 set the lower and upper range switchover limits CLD & FID Set the lower and upper range switchover limits for channel m NDIR |

1.5.4. EKEN: Set new device identification and information

| Command | Response | Description |
|--------------------------|----------|---|
| _EKEN_K0_new device-name | _EKEN_s | Set new device identification Maximum length of device name is 40 characters |

NOTE: To change device identification, you must first rename the device to “RESET”. Now a name up to 40 characters can be given.

NOTE: The device name must not have any blanks between characters, e.g. “CAI CLD” is not allowed. You can use underscores, e.g.. “CAI_CLD”.

1.5.5. EGRD: Set the range polynomial coefficients

| Command | Response | Description |
|----------------------------|----------|---|
| _EGRD_Km_Mn_a0_a1_a2_a3_a4 | _EGRD_s | M = 0 Set the polynomial coefficients of range ‘ n ’ CLD & FID Set the polynomial coefficients for range ‘ n ’ of channel m NDIR |

1.5.6. EFDA: Set auto calibration and purge times

| Command | Response | Description |
|-----------------------|----------|--|
| _EFDA_Km_SATK_z_y_x_w | _EFDA_s | m = 0 for CLD & FID m = channel for NDIR Set auto cal. times : z: Purge time y: Calibration time x: Total calibration time w: Verify time (z, y, x, w in seconds) |
| _EFDA_K0_SSPL_z | _EFDA_s | Set analyzer purge time to z seconds |

1.5.7. EPAR: Set auto calibration tolerance values

| Command | Response | Description |
|-------------------------------|----------|--|
| _EPAR_Km_SATK_z.z_y.y_x.x_w.w | _EPAR_s | m = 0 CLD & FID m = channel NDIR Auto calibration tolerance value (%): z.z= Range 1 y.y= Range 2 x.x= Range 3 w.w= Range 4 |

1.5.8. ESYZ: Set System Time

| Command | Response | Description |
|------------------------|----------|---|
| _ESYZ_K0_yymmdd_hhmmss | _ESYA_s | Set system time: yymmdd: year, month, day (each 2 characters wide, no spaces) hhmmss: hour, minutes, seconds (each 2 characters, no spaces) |

1.5.9. ET90 Set Lowpass Filter Time

| Command | Response | Description |
|------------|----------|---|
| _ET90_Kn_t | _ET90_s | CLD / FID n = 0 Set lowpass filter time: t= filter time in seconds NDIR n = 0 Set lowpass filter time for all Channels: n = Channel Set lowpass filter time for Channel n t= filter time in seconds |

1.5.10. EDAL: Diagnostic alarm limits

| Command | Response | Description |
|-----------------------------------|----------|-----------------------|
| _EDAL_K0_a1.min_a1.mas_..._a12max | _EDAL_s | Set all alarm limits |
| _EDAL_K0_x_x.min_xmax | _EDAL_s | Set alarm limits of x |

Alarm Limits:**NDIR**

| | | | |
|---|-------------------|-------|--|
| 1 | Channel 1 Flow | 7 | Temperature Failure |
| 2 | Channel 2 Flow | 8 | Channel 1 sample concentration |
| 3 | Channel 3 Flow | 9 | Channel 2 sample concentration |
| 4 | External Analog 1 | 10 | Channel 3 sample concentration |
| 5 | External Analog 2 | 11-13 | Temperature alarm limits channel 1...3 |
| 6 | Pressure Failure | 14-16 | EPC voltage alarm limits channel 1...3 |

CLD / HCLD

| | |
|----|-------------------------------|
| 1 | Simple Pressure Failure |
| 2 | Air Pressure Failure |
| 3 | Oven Temp. Failure |
| 4 | Converter Temp. Failure |
| 5 | Pump Temp. Failure |
| 6 | Diode Temp. Failure |
| 7 | Cell Temp Failure |
| 8 | Peltier Gas Temp Failure |
| 9 | Reaction Chamber Temp Failure |
| 10 | EPC Coil Sample Failure |
| 11 | EPC Coil Air Failure |
| 12 | Range Overflow |
| 13 | ADC Range Overflow |
| 14 | ADC Range Underflow |
| 15 | Range 1 is not calibrated |
| 16 | Range 2 is not calibrated |
| 17 | Range 3 is not calibrated |
| 18 | Range 4 is not calibrated |

FID / HFID

| | | | |
|---|-------------------------|----|----------------------------|
| 1 | No Flame | 9 | EPC Coil Sample Failure |
| 2 | Sample Pressure Failure | 10 | EPC Coil Air Failure |
| 3 | Air Pressure Failure | 11 | EPC Coil Fuel Failure |
| 4 | Fuel Pressure | 12 | Range Overflow |
| 5 | Burner Temp Failure | 13 | ADC Range Overflow |
| 6 | Oven Temp Failure | 14 | ADC Range Underflow |
| 7 | Cutter Temp Failure | 15 | Analyzer is not calibrated |
| 8 | Pump Temp Failure | | |

1.5.11. ETCP: Set TCP/IP Parameters

| Command | Response | Description |
|---|----------|--|
| _ETCP_K0_zzz.zzz.zzz.zzz _yyy.yyy.yyy.yyy _xxxx | _EDAL_s | zzz= TCP/IP address yyy= TCP/IP subnet mask xxxx= TCP/IP port All changes take effect after next power on cycle |

1.5.12. EH2O Set H₂O correction parameters

| Command | Response | Description |
|--------------------|----------|--|
| _EH2O_Km_z.z_y.x.x | _EH2O_s | m = channel 1 to 3 z.z: dry – voltage of A in with no water present y.y: 1st order coefficient x.x: 2nd order coefficient |

1.5.13. ECO2 Set CO₂ correction parameters

| Command | Response | Description |
|------------------------|----------|--|
| _EH2O_Km_z.z_y.x.x_w.w | _EH2O_s | m = channel 1 to 3 z.z: offset – voltage of A in with no CO ₂ present y.y: Min A in – if A in is below this value no CO ₂ correction will be done x.x: 1st order coefficient w.w.: 2nd order coefficient |

1.5.14. EDUP Set UDP Data streaming parameters

| Command | Response | Description |
|--|----------|---|
| _EUDP_K0_<UDPPort>_<DataFrequency>_[<Mode>]_[<UDP_IP>] | _EUDP_s | Configure an UDP channel for data streaming of the measuring values via Ethernet UDP. Port: port for open the UDP connection DataFrequency: Frequency for transmit the data in Hz Mode: A: ASCII Mode (optional) UDP_IP: Alternative IP address for open the UDP connection when it should not use the IP of connected TCP/IP client (optional) |

1.5.15. Format of the streaming Data via UDP:

ASCII Mode:

The measuring values will be sent with ASCII signs. The format is:

<Sequence number>_x.x_y.y_z.z

The sequence number will be incremented with every data packet, which is sent.

The measuring values x.x, y.y and z.z will be sent like in AKON K0 telegram

1.6. *Abbreviations used*

| | |
|-----------|----------------------------------|
| Km | : K' + channel Number (→K1....K3 |
| Mn | : Measuring range number |
| M1... M4 | : Measuring Range 1... 4 |
| w.w...zz | : Numerical value |
| W.W...Z.Z | : Numerical value |
| T | : Numeric integer value |
| x | : Number |
| a0... a4 | : Polynomial coefficients |
| s | : Status |